



[www.invacare.co.uk]

[www.invacare.ie]

Customer Service Tel.: 01656 776222

Ireland Tel.: +353 1 8107084



INVACARE® Limited
Pencoed Technology Park - Pencoed -
Bridgend CF35 5HZ
Tel: +44 1 656 776 222
Fax: +44 1656 776 220
E-mail: uk@invacare.com
www.invacare.co.uk
Sales Order E-mail: ordersuk@invacare.com

INVACARE® Ireland Ltd
Unit 5 - Seatown Business Campus
Seatown Road - SWORDS - County Dublin
Tel. +353 1 8107084
Fax +353 1 8107085
Email: Ireland@invacare.com
www.invacare.ie

© 2011 Invacare International Sarl. All rights reserved.
All information quoted is believed to be correct at time of print.
Invacare reserves the right to alter product specifications without prior consultation.
The reproduction or use of this work, or of any substantial part of this work, in any form or by any electronic, mechanical or other means, including photocopying, scanning and recording, and in any information storage and retrieval system, including electronic transmission, is forbidden without the prior written permission of Invacare International Sarl.

Independent Living Catalogue - UK - 04/2011 - 1504613 (2011 Version 1)

Independent Living Catalogue 2011-2012



Independent Living Catalogue - 2011-2012



Invacare® | *Product Catalogues*

Invacare®
Pressure Reducing Cushions and Accessories



Invacare®
**Pressure Reducing
Cushions and
Accessories**

Invacare® Independent Living Catalogue
2011-2013



Invacare®
**Independent Living
Catalogue
2011-2012**

WALKING AIDS CATALOGUE
2011



Invacare®
Walking aids catalogue



Customer Service UK

Tel: 01656 776222
Fax: 01656 776200
E-mail: ordersUK@invacare.com

Customer Service Eire

Tel: (353) 181 070 84
Fax: (353) 181 070 85
E-mail: ireland@invacare.com

For more information on all our products including user guides and spare parts lists, please refer to our web site:

www.invacare.co.uk

www.invacare.ie



Yes, you can.®

Invacare® / Aquatec® | **Independent Living Range**

	Introduction Bathing & Hygiene	6-7
Bath lifts	Aquatec Orca, Orca F, Orca XL	8-11
Transfer aids	Aquatec Trans	12
	Aquatec Disk	12
	Aquatec Step	13
Shower commode chairs	Aquatec Ocean, Ocean XL	14-15
	Aquatec Ocean Vip, Ocean Vip XL	16
	Aquatec Ocean Dual VIP	17
	Aquatec Ocean E-VIP	18
	Accessories Ocean	19-21
Commode chairs	Invacare Glide about commode	22
	Invacare Omega Commode range	23-24
	Invacare Commodes and Toilet frames	25
Shower stools, shower chairs	Aquatec Dot	26
	Aquatec Disk on Dot	26
	Invacare I-Fit, I-Fit	27
	Invacare Alize range, Galaxy	28-29
	Aquatec Pico	30
	Aquatec Pico 3 in 1 Shower/Commode chair	31
Folding shower seats	Invacare Futura R8802	32
	Aquatec Sansibar	32
	Invacare Futura R8804-60	32
Bathboard	Invacare Marina	33
Bathmat	Invacare Bula bath mat	33
Toilet seat raisers	Invacare Finesse	34
	Invacare Maris	34
	Aquatec 90	34-35
	Aquatec 900	36-37
	Aquatec 90000	36-37
Walking Rollators	Introduction Rollators	38-39
	Invacare P429/2 Delta	40-41
	Invacare P452E/3 Banjo	42-43
	DOLOMITE JAZZ	44-45
	DOLOMITE LEGACY	46-47
	DOLOMITE SYMPHONY	48-49
	DOLOMITE MAXI+	50-51
	DOLOMITE ALPHA	52-53
	DOLOMITE STEP UP	54-57

Walking frames	Introduction Walking Aids	58-59
	<i>Invacare Asteria</i>	60-61
	<i>Invacare Brio</i>	60-61
	<i>Invacare Aventia</i>	60-61
	<i>Invacare P435B Escort</i>	60-61
	<i>Invacare Foria</i>	60-61
Crutches	<i>Invacare Crutches</i>	62-63
Canes	<i>Invacare Canes</i>	64-65
Paediatric Walking Aids range	<i>Invacare Pluto</i>	66
	<i>Invacare Fun Kids Crutch</i>	67
	<i>Invacare Fixi and Fox Rollator</i>	67
	<i>Invacare Walk-On walking frame</i>	67
	Lifter	Introduction Hoists and Slings
	<i>Invacare Birdie® Family</i>	70-73
	<i>Invacare Jasmine</i> 	74-75
	<i>Invacare Roze</i> 	76-77
	<i>Invacare Reliant 350</i> 	78-79
	<i>Invacare Robin</i>	80-83
	<i>Invacare EC-Track™</i>	84-87
	<i>Invacare Slings</i>	88-95
Beds	Introduction Beds	96-97
	<i>Invacare Etude® Basic</i>	98-101
	<i>Invacare Etude® Medley Family</i>	102-105
	<i>Invacare Etude® Classic</i>	106-109
	<i>Invacare SB 755</i>	110-115
	<i>Invacare Octave®</i>	116-117
Beds Accessories	<i>Invacare Beds Accessories</i>	118-127
Mattresses	Introduction Mattresses	128-129
	<i>Invacare Softform® Premier</i>	130-131
	<i>Invacare Softform® Excel</i>	132-133
	<i>Invacare Softform® Premier MaxiGlide</i>	134-137
	<i>Invacare Softform® Premier Active</i>	138-141
	<i>Invacare Softform® Premier Visco</i>	142-143
	<i>Invacare Softform® Odstock Wedge</i>	144-145
	<i>Invacare Softform® Heelpad</i>	146-147
	<i>Invacare Softform® Flexipad</i>	148-149
	<i>Invacare Propad® Overlay</i>	150-151
	<i>Invacare Propad® Premier Overlay</i>	152-153
	<i>Invacare Propad® Visco Overlay</i>	154-155
	<i>Invacare Propad® Leg Trough</i>	156-157

Cushions	Introduction Cushions	158-159
	<i>Invacare Flo-tech® Lite</i>	160-161
	<i>Invacare Flo-tech® Contour</i>	162-163
	<i>Invacare Flo-tech® Lite Visco</i>	164-165
	<i>Invacare Flo-tech® Contour Visco</i>	166-167
	<i>Invacare Flo-tech® Plus</i>	168-169
	<i>Invacare Flo-tech® Image</i>	170-171
	<i>Invacare Flo-tech® Solution</i>	172-173
	<i>Invacare Flo-tech® Back Support</i>	174-175
	<i>Invacare Flo-tech® Lumbar Pad</i>	176-177
	<i>Invacare Propad® Original</i>	178-179
	<i>Invacare Propad® Low Profile</i>	180-181
	<i>Invacare Propad® Premier</i>	182-183
	<i>Invacare Propad® Visco</i>	184-185
Additional Information	Technical Information - Pressure Area Care	186-187
	Pressure Ulcer Risk Assessment	188-189
	Pictograms	190-191

Invacare® / Aquatec® | **Bathing & Hygiene**



AQUATEC®



Yes, you can.®

Invacare® / Aquatec® | **Bathing & Hygiene**

Bath lifts	Aquatec Orca, Orca F, Orca XL	8-11
Transfer aids	Aquatec Trans	12
	Aquatec Disk	12
	Aquatec Step	13
Shower commode chairs	Aquatec Ocean, Ocean XL	14-15
	Aquatec Ocean Vip, Ocean Vip XL	16
	Aquatec Ocean Dual VIP	17
	Aquatec Ocean E-VIP	18
	Accessories Ocean	19-21
Commode chairs	Invacare Glide about commode	22
	Invacare Omega Commode range	23-24
	Invacare Commodes and Toilet frames	25
Shower stools, shower chairs	Aquatec Dot	26
	Aquatec Disk on Dot	26
	Invacare I-Fit, I-Fit	27
	Invacare Alize range, Galaxy	28-29
	Aquatec Pico	30
	Aquatec Pico 3 in 1 Shower/Commode chair	31
Folding shower seats	Invacare Futura R8802	32
	Aquatec Sansibar	32
	Invacare Futura R8804-60	32
Bathboard	Invacare Marina	33
Bathmat	Invacare Bula bath mat	33
Toilet seat raisers	Invacare Finesse	34
	Invacare Maris	34
	Aquatec 90	34-35
	Aquatec 900	36-37
	Aquatec 90000	36-37

Aquatec® | Orca Bath Lift Range



Hand Control
Battery contained in
the hand set



3 year warranty



Aquatec Bath Lifters - Tried and trusted by health-care professionals.

The *Aquatec Orca* family is comprised of superior quality and durable bath lifters that are designed to offer comfort, safety and reassurance to those most in need of support as they bathe. The *Orca* incorporates increased leg room, ease of disassembly and an ergonomic, floating hand control that contains the light weight battery. The *Orca* has an impressive weight limit of 140 kg (22 stone) as standard, while the *Orca XL* will support patients up to 170 kg (26.5 stone).

An emphasis on ergonomic design allows the *Orca* to sit lower and further towards the back of the bath, giving increased leg room for deeper and more relaxed bathing. The increased height of the backrest provides full support whilst upright or reclined.

AQUATEC®



Yes, you can.®



Aquatec Orca

Reclining backrest

- Quick release hand control
- Battery with Sony™ technology
- Auto unlock for quick removal from bath
- **1470754 - White**
- **1470764 - Blue**



Aquatec Orca F

An upright bath lift with a fixed backrest, that offers support for those who prefer to bathe upright. The integrated head support of the high backrest ensures a safe and comfortable bath.

- **Aquatec Orca F with fixed backrest**
- **1471074 - White**



Aquatec Orca XL

A powerful bath lifter with maximum client weight of 170 kg (26.5 stone).

The reinforced scissor mechanism combines safe lifting with a compact construction. Minimum seat height of just 60 mm; allows plenty of room for a deeper and more comfortable bath.

- **Aquatec Orca XL with reclining backrest**
- **1471068 - White**

Features and Options



Compact hand control with integral battery

Orca bath lifts incorporate a revolutionary lightweight battery contained within the ergonomic hand control. Raised profile, coloured buttons and an emergency stop feature allow simple, controlled operation to give reassurance and independence.



Safe hand control fitting

The hand control connects within seconds to the back of the Orca, providing a secure and completely watertight seal. Effortlessly released with just one hand.



Two-piece design

When pushed forward the **Orca's** backrest folds down and automatically unlocks. The backrest and base section can then be removed from the bath with ease. The 2 sections fit neatly into one another to give an extremely compact size for storage.



Textured seating surface

The **Orca** has a textured surface allowing use without covers if desired. The side flaps are hinged flush to the main seating surface and allow comfortable and safe side transfers, without entrapment of skin.



Supported transfers

For those who require additional support during transfer, the **Aquatec Trans** allows the client to slide onto the seat and swivel into the bath with the minimum of effort.

■ **Product Code: 4.03.002**



Refresh kit

The **Orca** refresh kit comprises a full set of covers, oyster suction cups, side flap diverters and user guide. Ideal for re-issue.

■ **Product code: 1471017**



Height adaptors

The **Orca** can be adapted to suit a full range of bath heights. Height adaptors can be quickly fitted to the bath lift base.



Large oyster cups

These extra large oyster cups allow the **Orca** to adhere to steeply curved and textured baths providing exceptional stability and safety. Smaller cups are available as an option.

Features and Options



Transport bag

A sturdy canvas bag with wheels for discreet transport and safe storage.

■ **Product code: 14586AQT**



Non-slip cover mats

Cover mats are available in white or blue and are easily removed for machine washing at 60°C.

■ **Product codes:**

■ **1525933 (12cm)**

■ **1525932 (7cm)**



Fast Fit Pommel

A cost effective and easy to fit option (available in 70 mm and 120 mm diameter in blue colour only).

■ **Product code:**

■ **1525933**

■ **1525932**



Easy to clean

The open, clean-cut design of the **Orca** combines with its smooth surfaces to offer effortless cleaning.


Technical data



Orca/Orca F	360 mm	140kg (22.5 stone)
Orca XL	360 mm	170 kg (26.5 stone)

* Further technical information is available in the **Orca** user guide, which can be found on www.invacare.co.uk

Colours

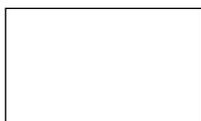


White



Blue

Aquatec® | **Bathing Transfer aids**

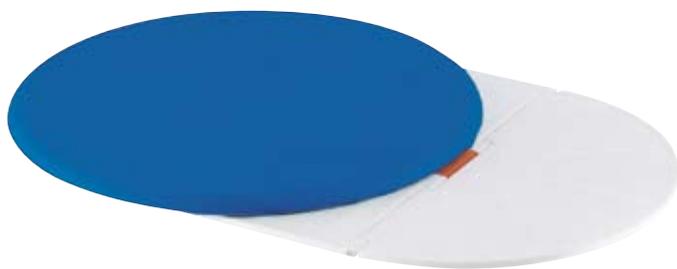

Aquatec **Trans**

- Compatible with the *Aquatec* range of bath lifters
- Large comfortable seating surface
- Safe: turntable is only unlocked by the weight of the user
- Seat surface turns and slides easily
- Soft, skin-friendly cover: easy-care, wipes clean

3 year warranty



Aquatec **Disk**

- Large comfortable seating surface
- Turns very easily
- Soft, skin-friendly cover: easy-care, wipes clean
- The *Aquatec Disk* is designed to be used in conjunction with the *Aquatec Dot* shower stool

3 year warranty

AQUATEC®



Yes, you can.®



Aquatec **Step**

- Stable frame
- Feet with non-slip rubber caps
- Non-slip, removable cover, machine-washable

3 year warranty

150kg/24 stone weight limit

Technical data

Order numbers	Aquatec Trans	Aquatec Disk	Aquatec Step
Turning and transfer aid, blue	4.03.002		
Turning and transfer aid, grey	4.03.008		
Rotating disk, blue		4.02.002	
Rotating disk, grey		4.02.008	
Step, blue			3.00.002
Dimensions			
Diameter	365 mm	365 mm	–
Length	530 mm	–	400 mm
Width	–	–	300 mm
Height	25 mm	25 mm	165 mm
More info			
Weight	2.2 kg	1.6 kg	1.5 kg
Load capacity	130 kg (21 stone)	130 kg (21 stone)	150 kg (24 stone)
Colours	Blue, grey	Blue, grey	Blue, white
Warranty	3 Years	3 Years	3 Years

Aquatec® | *Ocean family*



3 year warranty

150 kg weight limit
(23 stone)



Aquatec Ocean, Ocean XL, Ocean Vip, Ocean Vip XL, Ocean Dual Vip, Ocean E-Vip

The *Aquatec Ocean* series of shower commode chairs offers a flexible, modular solution for a wide range of clients.

AQUATEC®



Yes, you can.®



Aquatec **Ocean**

- Stable, stainless steel frame
 - Individually adjustable seat height: no tools required
 - Machine-washable adjustable backrest with quick-release buckles, can be adjusted to the user
 - Space-saving push bar to enable optimal positioning over the toilet
 - Ergonomic, swing-away armrests can be adjusted to two different height levels
 - Height-adjustable, swing-away and detachable footrests for simple transfer
 - Fold-away, textured footplates with heel loops provide a safe position
- The **Ocean** is available as follows:
- Standard version with castors, soft seat with insert, pan and lid.
Product code: 1517044
 - Self-propel model with 24" wheels, quick-release-axles and four castors, soft seat with insert, pan and lid.
Product code: 1517045
 - XL version (180 kg load capacity and more space between the armrests) with castors, soft seat with insert, pan and lid.
Product code: 1517058
 - All models can be retrofitted with 24" wheels and/or to the XL version

Options and features



Ocean 24" wheels
Self-propel model with 24" wheels, quick-release axles and four castors.



Ocean XL
More space between the armrests with a load capacity of up to 180 kg.



Individually adjustable seat height: no tools required.



Easy disassembly for mobile use, e.g. whilst travelling.



Stable, stainless steel frame.



Aquatec **Ocean Vip**

- Stable, stainless steel frame
- Seat tilt can be individually adjusted (tilt angle 0° - 35°)
- Stainless steel gas spring
- Extremely stable even at maximum tilt angle
- Seat height adjustable in 3 steps, no tools required
- Very comfortable, ergonomically contoured seat
- Machine-washable adjustable backrest with quick-release buckles – can be easily adjusted to the user

The **Ocean Vip** is available as follows:

- Standard version with castors, toilet pan holder and soft seat, insert, pan and lid.
Product code: 1517060
- **XL** version (more space between the armrests) with castors, soft seat with insert, pan and lid.
Product code: 1517061
- The **Ocean Vip** can be retrofitted with XL armrests

Options and features



Seat tilt

Tilt angle 0° - 35°.



Headrest

Height and depth-adjustable headrests included.



Directional lock

Easy to operate brakes on 3 castors, 1 castor with directional lock to help push the chairs down straight floors.



Adjustable seat height

Adjustable in 3 steps, no tools required.



Aquatec **Ocean Dual-Vip**

- Seat tilt can be individually adjusted (tilt angle 0° - 35°)
- Two additional stainless steel gas springs to lower the backrest (recline angle 10° - 50°)
- Second operating lever with safety lock to lower the backrest
- Particularly stable even at maximum tilt angle, anti-tippers included
- Machine-washable adjustable backrest with quick-release buckles – can be adjusted to the user

The **Ocean Dual Vip** is available as follows:

- Standard version with castors, soft seat with insert, pan and lid.
Product code: 1517062
- The **Ocean Vip** can be retrofitted with reclining backrests for the **Ocean Dual Vip** and be fitted with anti-tippers
- The **Ocean Dual Vip** can be retrofitted with XL armrests

Options and features



Seat surface 0°, backrest 25°



Full seat tilt 35°, back rest 50°



Anti-tippers included



Aquatec **Ocean E-Vip**

- Seat height can be steplessly adjusted between 500 and 1050 mm electrically
- Seat tilt can be electrically tilted forwards and backwards (tilt angle -5° to 35°)
- Particularly stable even at maximum tilt angle and maximum seat height
- Clearly laid out hand control for simple adjustment of the seat height and seat tilt (4 buttons only)
- Sealed electronic components: ideal for use in the shower
- Can be cleaned in automatic cleaning stations at temperatures up to 85 °C!
- Detachable battery pack can be recharged easily

The **Ocean E-Vip** is available as follows:

- Standard version with castors, soft seat with insert, pan and lid.
- **Product code: 1518424**
- The **Ocean E-Vip** can be retrofitted with XL armrests

Options and features



Maximum tilt angle and maximum seat height



Detachable battery pack



Clearly laid out hand control

Accessories Aquatec Ocean



Soft seat

With key-shaped cut out.
16341 for **Ocean / Ocean XL / Ocean E-Vip**
1470075 for **Ocean Vip / Ocean Vip XL / Ocean Dual Vip**



Soft seat insert

For soft seat with key-shaped cut out.
16343 for all models of the **Ocean** series



Soft seat

With small, oval opening.
16342 for all models of the **Ocean** series



Splash guard

Simply snap-on, for soft seat with key-shaped cut out.
16511 for all models of the **Ocean** series

Accessories Aquatec Ocean



Universal soft seat
Closed seat surface.
1526032 for all models of the *Ocean* series



Ergonomic special soft seat
Extra soft, textured surface
1531259 soft for *Ocean*, *Ocean XL*, *Ocean Vip*, *Ocean Vip XL*, *Ocean Dual Vip*



Safety bar
For additional safety. The surface with soft padding. Easy to assemble.
10167-10 for all models of the *Ocean* series. Length: 540 mm
1470078 XL safety bar for all XL versions of the *Ocean* series. Length: 620 mm



Safety belts
For safe positioning. Individually adjustable, high load capacity and easy to assemble. Machine-washable. For all models of the *Ocean* series.
1470081 Chest belt, grey
1470082 Pelvic belt, grey



Armrest locking device
To lock the armrest.
1470077 for all models of the *Ocean* series



XL back
Increases the width of the backrest by 80 mm. Retrofit kit for *Ocean*.
1534979
Not available until October 2011



Wider XL armrests
Increases the distance between the armrests by 40 mm on each side. Retrofit kit for *Ocean Vip*, *Ocean Dual Vip* and *Ocean E-Vip*
1470086 Wider XL armrests (pair)



24" wheels (set)
Incl. brakes and removable axles.
16379 for *Ocean / Ocean XL*



Amputation support (shown for Ocean)
Padded, both height and angle adjustable.
16616 for *Ocean / Ocean XL* (left or right)
1470083 left for *Ocean Vip*, *Ocean Vip XL*, *Ocean Dual Vip* and *Ocean E-Vip*
1470827 right for *Ocean Vip*, *Ocean Vip XL*, *Ocean Dual Vip* and *Ocean E-Vip*



Calf support (shown for Ocean)
Padded, both height and angle adjustable.
16615 for *Ocean / Ocean XL* (left or right)
1470084 left for *Ocean Vip*, *Ocean Vip XL*, *Ocean Dual Vip* and *Ocean E-Vip*
1470828 right for *Ocean Vip*, *Ocean Vip XL*, *Ocean Dual Vip* and *Ocean E-Vip*



Toilet pan
With handle and cover.
10230 for all models of *Ocean* series. Length: 460 mm



Head/Neckrest (neckrest shown)
1531465 Headrest for *Ocean Vip*, *Ocean Vip XL*, *Ocean Dual Vip* and *Ocean E-Vip*
1531608 Neckrest for *Ocean Vip*, *Ocean Vip XL*, *Ocean Dual Vip* and *Ocean E-Vip*


Toilet pan holder (not shown)
16335 for *Ocean / Ocean XL*
1470076 for *Ocean Vip / Ocean Vip XL / Ocean Dual Vip*
1471233 for *Ocean E-Vip*

Technical data


Dimensions	Aquatec <i>Ocean</i>	Aquatec <i>Ocean 24" wheels</i>	Aquatec <i>Ocean XL</i>
Seat width	480 mm (A)	480 mm (A)	480 mm (A)
Seat depth	450 mm (B)	450 mm (B)	450 mm (B)
	450 mm (B)	475 - 600 mm (C)	475 - 600 mm (C)
Total width	560 mm (D)	675 mm (D)	640 mm (D)
Total depth	900 mm (E)	1060 mm (E)	900 mm (E)
Total height	940 - 1070 mm (F)	940 - 1070 mm (F)	940 - 1070 mm (F)
Depth without footrests	700 mm (G)	855 mm (G)	700 mm (G)
Footplate depth	155 mm (H)	155 mm (H)	155 mm (H)
Height to armrests	680 - 850 mm (I)	680 - 850 mm (I)	680 - 850 mm (I)
Height, footrests to seat	255 - 495 mm (J)	255 - 495 mm (J)	255 - 495 mm (J)
Height, ground to lower seat edge:			
with holder for commode pan	390 - 520 mm (K)	390 - 520 mm (K)	390 - 520 mm (K)
without holder for commode pan	410 - 540 mm (K)	410 - 540 mm (K)	410 - 540 mm (K)
Height, seat to armrest	210 / 250 mm (L)	210 / 250 mm (L)	210 / 250 mm (L)
Width between the armrests	460 mm (M)	460 mm (M)	540 mm (M)
Width, cut out in seat plate	205 mm	205 mm	205 mm
Depth, cut out in seat plate	320 mm	320 mm	320 mm
Graduation, height adjustment	25 mm	25 mm	25 mm
Ø Castors	5" (127 mm)	5" (127 mm)	5" (127 mm)
Ø Drive wheels	-	24" (610 mm)	-
Width between the 24" wheels	-	540 mm (N)	-
More info			
Weight	14.5 kg	19.5 kg	14.5 kg
Load capacity	150 kg (23 stone)	150 kg (23 stone)	180 kg (28 stone)
Colour	White	White	White
Warranty	3 years	3 years	3 years

Heel loops are included in the scope of delivery. The toilet pan holder and toilet pan are not.


Ocean with castors

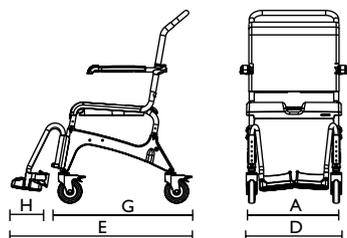

Ocean with 24" wheels

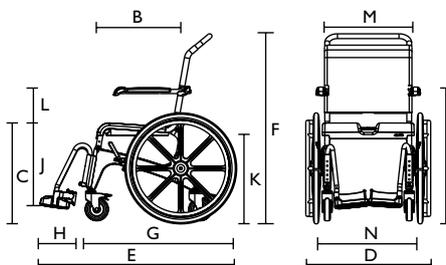


Ocean XL

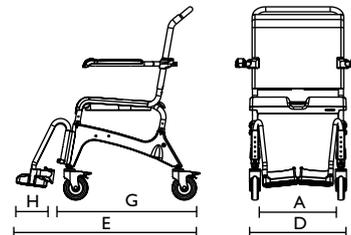


Technical data

Order numbers	Aquatec Ocean Vip with tilt	Aquatec Ocean Vip XL with tilt	Aquatec Ocean Dual Vip with tilt & recline	Aquatec Ocean E-Vip
Shower commode chair with castors	1470000	1470788	1470829	1471208
Shower commode chair with castors and toilet pan holder	1470712			1522373
Shower commode chair with castors, toilet pan holder and soft seat	1470713			
Dimensions				
Seat width	500 mm (A)	500 mm (A)	500 mm (A)	480 mm (A)
Seat depth	450 mm (B)	450 mm (B)	450 mm (B)	450 mm (B)
Seat height	510 - 610 mm (C)	510 - 610 mm (C)	510 - 610 mm (C)	500 - 1050 mm (C)
Total width	595 mm (D)	625 mm (D)	625 mm (D)	730 mm (D)
Total depth	1060 mm (E)	1060 mm (E)	1060 mm (E)	1010 mm (E)
Total height without adjustable headrest	1085 - 1185 mm (F)	1085 - 1185 mm (F)	1085 - 1185 mm (F)	1070 - 1620 mm (F)
Depth without footrests	845 mm (G)	845 mm (G)	845 mm (G)	1010 mm (G)
Footplate depth	160 mm (H)	160 mm (H)	160 mm (H)	155 mm (H)
Height to armrest	720 - 860 mm (I)	720 - 860 mm (I)	720 - 860 mm (I)	730 - 1320 mm (I)
Height, footrests to seat	345 - 465 mm (J)	345 - 465 mm (J)	345 - 465 mm (J)	345 - 465 mm (J)
Height, ground to lower seat edge: with holder for toilet pan	420 - 520 mm (K)	420 - 520 mm (K)	420 - 520 mm (K)	420 - 970 mm (K)
without holder for toilet pan	440 - 540 mm (K)	440 - 540 mm (K)	440 - 540 mm (K)	465 - 1015 mm (K)
Height, seat to armrest	215 / 245 mm (L)	215 / 245 mm (L)	215 / 245 mm (L)	230 - 270 mm (L)
Width between the armrests	475 mm (M)	555 mm (M)	475 mm (M)	470 mm (M)
Width, cut out in seat plate	205 mm	205 mm	205 mm	205 mm
Depth, cut out in seat plate	320 mm	320 mm	320 mm	320 mm
Tilt angle seat	0° - 35°	0° - 35°	0° - 35°	-
Tilt angle back	-	-	10° - 50°	-
Graduation, height adjustment	50 mm	50 mm	50 mm	stepless
Ø Castors	5" (127 mm)	5" (127 mm)	5" (127 mm)	5" (127 mm)
More info				
Weight	23.5 kg	23.5 kg	28.3 kg	39 kg
Load capacity	150 kg (23 stone)	150 kg (23 stone)	150 kg (23 stone)	150 kg (23 stone)
Colour	White	White	White	White
Warranty	3 years	3 years	3 years	3 years

Ocean Vip / Vip XL / Dual Vip

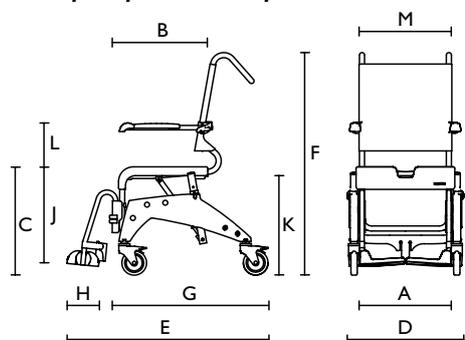

Ocean E-Vip

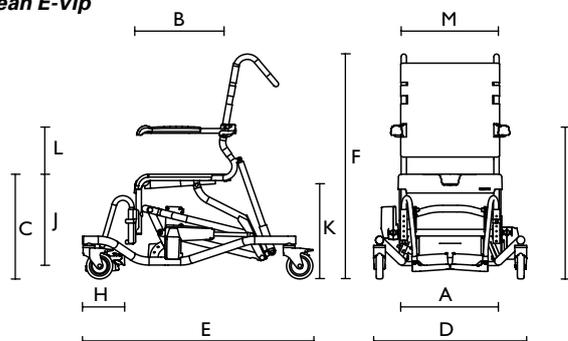

Invacare® | **Commode chairs**



Supplied with leg rests as standard



Invacare **Glideabout commode**

- Stable, robust frame made from powder-coated steel
- High corrosion protection
- Space-saving push bar for optimal positioning over the toilet
- Swing-away armrests and swing-away, detachable footrests
- Stable commode seat with integrated toilet pan holder
- Supplied with 4 braked castors as standard



Yes, you can.®



Invacare **Omega mobile**

- Compact modern design
- Castors with rear brakes for safety
- Padded back, seat and arms
- Integral foot support



Invacare **Omega classic**

- Static commode chair
- Space-saving: can be stacked
- Stable and corrosion-protected design with welded tubes
- Easy handling: commode pan can be removed from the side and from above



Invacare **Omega height adjustable**

- Height adjustable commode chair
- Space-saving: can be stacked
- Seat height adjustable in 7 steps from 510 - 660 mm
- Attractive design: fits well into a domestic setting
- Easy to clean



Invacare **Omega Piante**

- Foldable, Space saving commode
- Ideal for infrequent or travel use
- Contoured backrest
- Complete with sanitary pan


Technical data

Order numbers	Invacare Glideabout commode	Invacare Omega H750	Invacare H450 stationary	Invacare H450LA height adjustable	Invacare Omega H407
Commode chair, mobile	1525770				
Omega mobile commode chair with castors		1471044			
Static commode			1471164	1471166	
Omega height adjustable commode chair					1471010 (2 pack)


Dimensions					
Seat width	450 mm (A)	500 mm (A)	495 mm (A)	495 mm (A)	440 mm (A)
Seat depth	450 mm (B)	460 mm (B)	455 mm (B)	455 mm (B)	420 mm (B)
Seat height	480 mm (C)	480 mm (C)	500 mm (C)	510 - 660 mm (C)	450 mm (C)
Total width	575 mm (D)	635 mm (D)	670 mm (D)	670 mm (D)	580 mm (D)
Total depth	810 mm (E)	790 mm (E)	630 mm (E)	630 mm (E)	560 mm (E)
Total height	945 mm (F)	920 mm (F)	890 mm (F)	900 - 1050 mm (F)	910 mm (F)
Depth without footrests	545 mm (G)	670 mm (G)	-	-	-
Footplate depth	155 mm (H)	95 mm (H)	-	-	-
Height to armrests	755 mm (I)	700 mm (I)	690 mm (H)	700 - 850 mm (H)	660 mm (H)
Height, footrests to seat	345 - 460 mm (J)	370 mm (J)	-	-	-
Height, ground to lower seat edge: with holder for toilet pan	420 mm (K)	360 mm (K)	370 mm (I)	380 - 530 mm (I)	-
without holder for toilet pan	-	-	430 mm (J)	440 - 590 mm (J)	-
Height, floor to cover plate	525 mm (L)	505 mm (L)	525 mm (K)	535 - 685 mm (K)	475 mm (K)
Width between the armrests	470 mm (M)	510 mm (M)	510 mm (G)	510 mm (G)	460 mm (G)
Width, cut out in seat plate	220 mm	215 mm	215 mm	215 mm	215 mm
Depth, cut out in seat plate	260 mm	255 mm	250 mm	250 mm	255 mm
Ø Castors	5" (127 mm)	100 mm	-	-	-

More info					
Weight	14,0 kg	15,0 kg	8,3 kg	8,8 kg	8,5 kg
Load capacity	135 kg (22 stone)	130 kg (21 stone)	130 kg (21 stone)	130 kg (21 stone)	80 kg (13.5 stone)
Colour	Black anthracite	Lavender blue / Blue	Lavender blue / Blue	Lavender blue / Blue	Lavender blue
Warranty	2 years	2 years	2 years	2 years	2 years


Invacare **Glideabout commode**

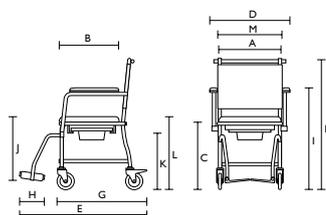


Invacare **Omega H750**

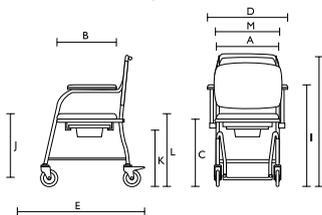

Invacare **H450 / H450LA**

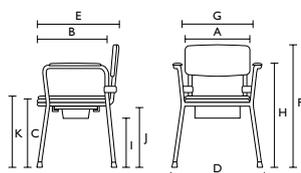

Invacare **Omega H407**

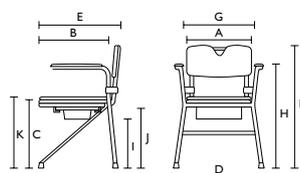



Invacare Styxo 9630E

- Height adjustable from 390 mm to 540 mm
- Easy to clean
- Complete with seat, sanitary pan and lid



Invacare Adeo C407A

- Height adjustable 640 - 740 mm
- Fits most standard toilets
- Complete with armrests




Invacare Izzo H340

- Height adjustable 450 - 550 mm
- Complete with swing away armrests
- Sanitary pan and lid included


Technical data

Order numbers	Invacare Styxo	Invacare Adeo	Invacare Izzo
Toilet frame	9630E	C407A	H340
Dimensions			
Seat width	380 mm	590 mm	420 mm
Seat depth	350 mm	440 mm	350 mm
Seat height	410-560 mm	640 - 740 mm	450 - 550 mm
Total width	620 mm	590 mm	530 mm
Total depth	490 mm	440 mm	490 mm
Total height	580 - 730 mm	740 mm	920 mm
More info			
Weight	5,5 kg	2kg	8.8kg
Load	135 kg (21 stone)	100kg (16 stone)	100kg (16 stone)
Colour	white / grey	Chrome	white / grey
Warranty	2 years	2 years	2 years


Invacare Styxo 9630E

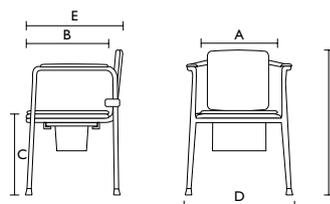

Invacare Adeo C407A

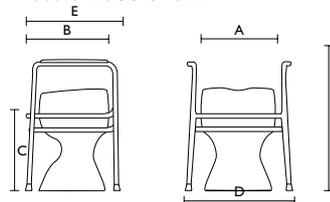

Invacare Izzo H340

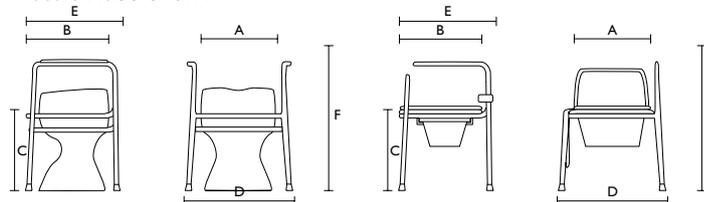

Aquatec® | Shower stools, shower chairs



3 year warranty

Aquatec **Dot**

- Round shape, fits into most showers
- Easily disassembled: ideal for travelling
- Textured seat surface for a safe seat
- Rubber feet provide a secure stand




Aquatec **Disk on Dot**

- Within seconds, the **Disk**, the shower stool **Dot** and an adapter are combined to create a shower stool with swivel seat
- High load capacity: 125 kg (20 stone)

Technical data

Order numbers	Aquatec Dot	Aquatec Disk on Dot Seat height 550 mm	Aquatec Disk on Dot Seat height 500 mm
Shower stool round, seat height 500 mm	10200-32		
Shower stool round, seat height 550 mm	10201-32		
Blue Disk + Dot + Adapter		9.01.002	9.00.002
Grey Disk + Dot + Adapter		9.01.008	9.00.008
Dimensions			
Seat width	400 mm (A)	400 mm (A)	400 mm (A)
Seat depth	400 mm (A)	400 mm (A)	400 mm (A)
Seat height	500 / 550 mm (B)	575 mm (B)	525 mm (B)
Total width	460 mm (C)	460 mm (C)	460 mm (C)
Total depth	460 mm (C)	460 mm (C)	460 mm (C)
Total height	500 / 550 mm (B)	575 mm (B)	525 mm (B)
More info			
Weight	2.5 kg	4.1 kg	4.1 kg
Load capacity	125 kg (19.5 stone)	125 kg (19.5 stone)	125 kg (19.5 stone)
Colour	White	Blue, grey	Blue, grey
Warranty	3 years	3 years	3 years
Accessories			
Disk , blue	4.02.002	-	-
Disk , grey	4.02.008	-	-
Adapter for Disk on Dot	15326AQT	-	-

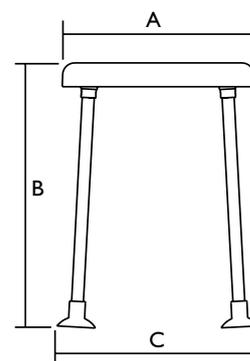

AQUATEC®



Yes, you can.®



Invacare® I-Fit 9780E

- Height Adjustable from 380 mm to 550 mm
- Large ergonomic shaped seat for comfort
- Easy to clean
- Light weight and easy to transport
- Stylish modular design
- Hygiene recess



**Substantial weight limit
180kg (28 stone)**


Invacare I-Fit 9781E

- Height Adjustable from 380 mm to 550 mm
- Removable backrest
- Ergonomic shaped seat with backrest for extra comfort
- Easy to clean
- Light weight and easy to transport


Technical data

Order numbers	Invacare I-Fit	Invacare I-Fit
Shower stool	9780E	
Shower chair		9781E
Dimensions		
Seat width	340 mm	340 mm
Seat depth	370 mm	370 mm
Seat height	380-550 mm	380-550 mm
Total width	500 mm	500 mm
Total depth	430 mm	430 mm
Total height	380-550 mm	760-930 mm
More info		
Weight	3,7 kg	4,5 kg
Load capacity	180 kg (28 stone)	180 kg (28 stone)
Colour	white	white
Warranty	2 years	2 years

Invacare I-Fit 9780E

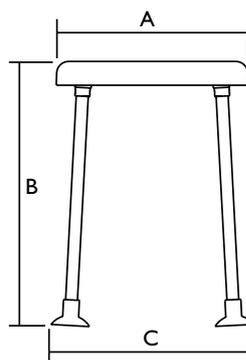

Invacare I-Fit 9781E

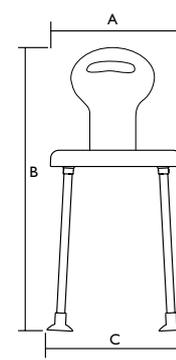



Invacare® Alize

- Stable shower stool
- Large ergonomic seat with hygiene recess
- Height adjustable corrosion protected legs



Invacare Alize (with arms)

- Stable shower stool
- Large ergonomic seat surface with hygiene recess makes personal hygiene easier
- Fixed armrests
- Height-adjustable, corrosion-protected legs



Invacare Alize shower seat

- Ergonomic seat surface with hygiene recess
- Height-adjustable, corrosion-protected legs
- Height adjustable 500 - 560 mm



3 year warranty


Aquatec Galaxy

- Folding, height-adjustable shower chair, perfect for travelling and compact for storage
- Ergonomic seat surface, can be supplied with or without hygiene recess
- Soft backrest for additional support


Technical data

Order numbers	Invacare Alize stool	Invacare® Alize stool c/w arms	Invacare Alize chair	Aquatec Galaxy
Shower stool with hygiene recess	H2120/1	H2240/1		
Shower stool with armrests and hygiene recess		H2240/1		
Shower chair with backrest and armrests, with hygiene recess			H2480/1	1470805
Dimensions				
Seat width	410 mm (A)	410 mm (A)	410 mm (A)	400 mm (A)
Seat depth	400 mm (B)	400 mm (B)	400 mm (B)	345 mm (B)
Seat height	470 - 550 mm (C)	470 - 550 mm (C)	500 - 560 mm (C)	470 - 620 mm (C)
Total width	590 mm (D)	590 mm (D)	590 mm (D)	580 mm (D)
Total depth	480 mm (E)	480 mm (E)	480 mm (E)	460 mm (E)
Total height	550 mm (C)	550 mm (F)	940 mm (F)	840 - 990 mm (F)
Graduation, height adjustment				
Height to armrests				630 - 780 mm (G)
Width between the armrests		450 mm (G)	450 mm (G)	510 mm (H)
Width, hygiene recess	140 mm	140 mm	140 mm	130 mm
Depth, hygiene recess	170 mm	170 mm	170 mm	140 mm
More info				
Weight	4.5kg	5.7 kg	7.5kg	7.5 kg
Load capacity	120kg (19 stone)	120 kg (19 stone)	120kg (19 stone)	100 kg (16 stone)
Colour	White	White	White	White
Warranty	2 years	2 years	2 years	3 years


Invacare **Alize**

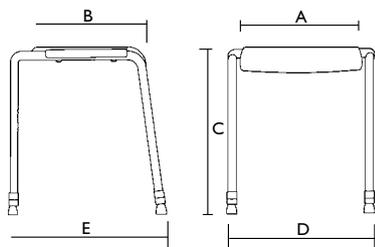

Invacare **Alize with arms**

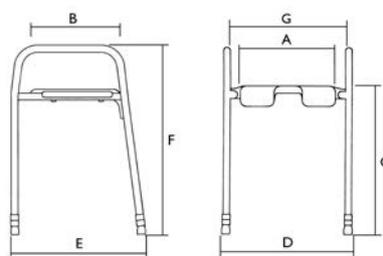

Invacare **Alize shower seat** / Aquatec **Galaxy**

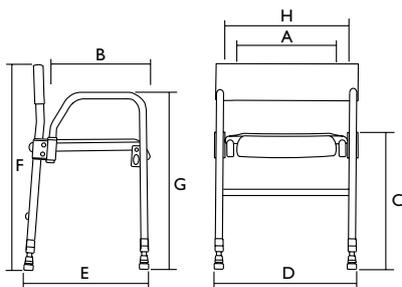



Aquatec **Pico**

- Modular shower chair system
- Four configurations
- Textured seat surface for safe sitting
- Large hygiene recess, either at the front or back, due to two mounting options for the backrest
- Height-adjustable legs
- Rubber feet provide a secure stand
- Simple and easy assembly, no tools required
- Load capacity of up to 160 kg (25 stone)

Options and features



Pico



Pico with armrests



Pico with backrest

3 year warranty



Technical data

Order numbers	Aquatec Pico
Shower stool	9.20.001
Shower stool with armrests	9.21.001
Shower chair with backrest	9.23.001
Shower chair complete	9.22.001
Backrest	16469
Armrests (pair)	16468

Dimensions	
Seat width	430 mm (A)
Seat depth	420 mm (B)
Seat height	425 - 575 mm (C)
Total width	575 mm (D)
Total depth	520 mm (E)
Total height	825 - 975 mm (F)
Height to armrests	605 - 755 mm (G)
Width between the armrests	430 mm (H)
Width, hygiene recess	180 mm
Depth, hygiene recess	145 mm
Graduation, height adjustment	25 mm

More info	
Weight, shower stool	3,2 kg
Weight, shower stool with armrests	4,1 kg
Weight, shower chair with backrest	4,0 kg
Weight, shower chair complete	4,8 kg
Load capacity	160 kg! (25 stone)
Colour	White
Warranty	3 years

Aquatec **Pico**

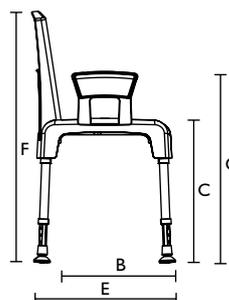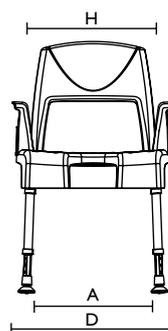



Aquatec **Pico 3 in 1**

- Stationary, height adjustable multifunctional shower chair
- Versatile: shower chair, commode or commode support frame
- Easy handling: toilet pan can be removed from all four sides
- Absolutely corrosion-free
- Armrests, backrest, toilet pan and toilet pan holder are included
- Seat height adjustable in 6 steps from 425 – 575 mm
- Very high ILoad capacity of 160 kg (25 stone)

Options and features



Pico 3 in 1 with one armrest



Pico 3 in 1 as commode support frame



Technical data

Order numbers	Aquatec Pico 3 in 1
Shower chair complete	1525887
Universal soft seat	1526032

Dimensions	
Seat width	430 mm (A)
Seat depth	420 mm (B)
Seat height	425 - 575 mm (C)
Total width	575 mm (D)
Total depth	520 mm (E)
Total height	825 - 975 mm (F)
Height to armrests	605 - 755 mm (G)
Width between the armrests	430 mm (H)
Graduation, height adjustment	25 mm
Height, ground to lower seat edge:	
with holder for toilet pan	340 - 490 mm
without holder for toilet pan	375 - 525 mm
Width, cut out in seat plate	215 mm
Depth, cut out in seat plate	250 mm
Graduation, height adjustment	25 mm

More info	
Weight, shower chair complete	4,8 kg
Load capacity	160 kg! (25 stone)
Colour	White
Warranty	3 years

Aquatec **Pico 3 in 1**

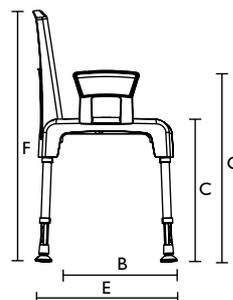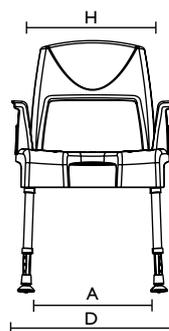

Pico 3 in 1 with universal soft seat

Invacare® | **Folding shower seats, bathboard, bathmat**



Invacare **R8802**

- Space-saving: mounted on the wall and can be folded when not in use
- Anatomically shaped seat surface for greater comfort when seated
- Simple and easy assembly (3 screws only)*



Invacare **R8804-60**

- Space-saving: mounted on the wall and can be folded when not in use
- Anatomically shaped seat surface for greater comfort when seated
- Armrests for secure, lateral hold
- Transfer is made easier by individual fold-away armrests

*Screws are not included.



3 year warranty

Aquatec® **Sansibar**

- Space-saving: mounted on the wall and can be folded when not in use
- Large hygiene recess makes personal hygiene easy
- Height-adjustable support legs for individual adjustment and additional stability
- Compensates for floor unevenness through finely adjustable feet

AQUATEC®



Yes, you can.®



Invacare **Marina**



Invacare **H190**

- Attractively designed bathtub board
- Non-slip seat surface
- Hand grip provides additional safety
- Practical shower head holder on both sides
- Individually adjustable using adjustment elements
- Also available as an XL version

- High-quality bathtub mat
- Held securely in the bathtub by numerous small suction cups on the underside of the mat
- Non-slip surface for a secure hold while bathing or showering
- Extra long: covers the entire base of the bathtub



Technical data

Order number	Invacare R8802	Aquatec® Sansibar	Invacare R8804-60	Invacare Marina / Marina XL	Invacare H190
Folding shower seat	R8802	9.10.101	R8804-60		
Bathtub board with hand grip				1525940	
Bathtub board XL				1471594	
Bathtub mat (6 pack)					H190



Dimensions					
Seat width	380 mm (A)	400 mm (A)	380 mm (A)	–	410 mm (A)
Seat depth	370 mm (B)	410 mm (B)	370 mm (B)	–	810 mm (B)
Seat height	–	460 - 530 ¹ mm (C)	440 - 550 mm (C)	–	10 mm (C)
Total width	385 mm (C)	400 mm (D)	545 mm (D)	330 mm (A)	410 mm
Total depth	470 mm (D)	460 mm (E)	470 mm (E)	690 / 740 mm (B)	810 mm
Total height	180 mm (E)	530 - 600 mm (F)	770 - 880 mm (F)	160 mm (C)	10 mm
Width, hand grip	–	–	–	180 mm (D)	–
Adjustment range, inner dimension	–	–	–	390 - 630 / 395 - 630 mm (E)	–
Graduation, height adjustment	–	25 mm	–	–	–
Stepless height adjustment	–	–	0 - 110 mm	–	–
Height to armrests	–	–	700 - 810 mm (G)	–	–
Height, seat to backrest	–	–	150 mm (H)	–	–
Height, backrest	–	–	180 mm (I)	–	–
Width between the armrests	–	–	460 mm (J)	–	–
Depth, folded	105 mm	110 mm	115 mm	–	–
Height, folded	520 mm	540 - 610 mm	1080 - 1190 mm	–	–

More info					
Weight	4.8 kg	2.6 kg	9 kg	1.6 kg	1.7 kg
Load capacity	135 kg (21 stone)	130 kg (19 stone)	135 kg (21 stone)	150 kg (24 stone)	130 kg (19 stone)
Colour	White	White	White	White	Lavender blue
Warranty	3 years	2 years	2 years	2 years	2 years


Invacare® **R8802**

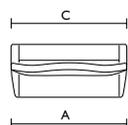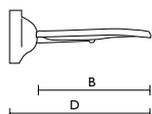


Aquatec® **Sansibar**

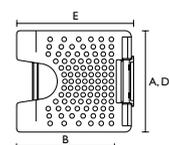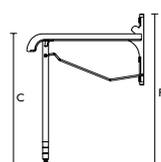

Invacare **R8804-60**

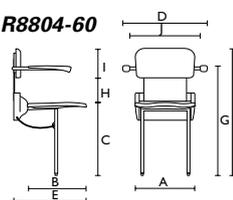

Invacare **Marina**

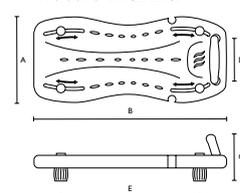

Invacare® | Toilet seat raisers



Invacare **Finesse**

- Toilet seat raiser made of soft, skin-friendly material
- Raises the toilet seat by 90 mm
- Fits on the toilet easily and safely, no assembly required
- Large hygiene recess at the front for easy personal hygiene



Invacare **Maris**

- Toilet seat raiser
- Very easy to attach to the toilet
- Smooth surface for easy cleaning
- Large hygiene recess (250 x 190 mm)



Aquatec® **90**

- Toilet seat raiser
- Raises the toilet seat by 100 mm
- Attached safely to the toilet with two side clamps
- Large turning knobs for quick attachment and detachment
- Particularly high load capacity: 225 kg (35 stone)



**High weight limit
225kg (35 stone)**

3 year warranty



AQUATEC®




Yes, you can.®


Technical data

Order number	Invacare H304	Invacare H301	Aquatec® 90
Toilet seat raiser	1470982		
Toilet seat raiser with lid			8.01.001
Toilet seat raiser without lid		1522262	8.00.001
Dimensions			
Seat width	385 mm (A)	365 mm (A)	365 mm (A)
Seat depth	430 mm (B)	385 mm (B)	395 mm (B)
Seat height	90 mm (C)	–	–
Seat height front	–	125 mm (C)	105 mm (C)
Seat height rear	–	135 mm (D)	115 mm (D)
Total width	400 mm	380 mm	–
Total width (with lid)	–	–	380 mm (E)
Total width (without lid)	–	–	375 mm
Total depth	435 mm (F)	400 mm (F)	410 mm (F)
Total height	135 mm (G)	185 mm (G)	–
Total height (with lid)	–	–	135 mm
Total height (without lid)	–	–	125 mm (G)
Width of opening	225 mm (H)	190 mm (H)	210 mm (H)
Depth of opening	275 mm (I)	250 mm (I)	270 mm (I)
More info			
Weight	2.2kg	0,9 kg	–
Weight (with lid)	–	–	1,6 kg
Weight (without lid)	–	–	1,1 kg
Load capacity	125 kg (19.5 stone)	120 kg (19 stone)	225 kg! (35.5 stone)
Colour	Grey	White	White
Warranty	2 years	2 years	3 years


Invacare **H304**

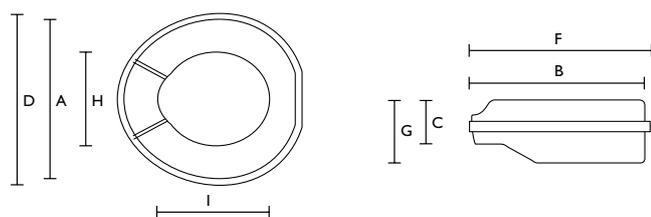

Invacare **H301**

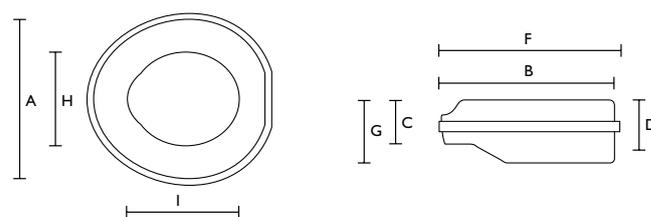

Aquatec® **90**

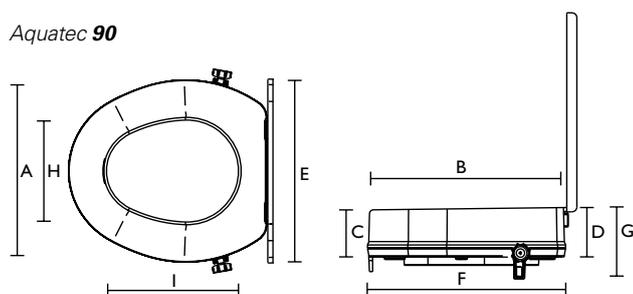



3 year warranty

Aquatec® 900

- Securely attached toilet seat raiser
- Individual adjustment to the user at any time, no tools required:
 - three different seat heights: 2, 4 and 6"
 - three different angles, lower at the front
- Fold-way, ergonomic armrests for additional support when sitting down and standing up
- Available with or without armrests.
- Seat can be folded up in all setting positions – ease of cleaning



Aquatec 90000

- Toilet seat raiser
- Available in three different seat heights: 2, 4, 6"
- Swing-away, ergonomic armrests for additional support when sitting down and standing up
- Innovative clip system: The lid and toilet seat can be removed from the base; no tools are required


Technical data

Order numbers	Aquatec® 900	Aquatec 90000
Toilet seat raiser with armrests	10128-10	8.31.101 (seat height 100 mm) 8.31.111 (seat height 60 mm) 8.31.121 (seat height 20 mm)
Toilet seat raiser without armrests	10129-10	
Accessories: Adapter for forward-inclined assembly		1470825


Dimensions		
Seat width	360 mm (A)	365 mm (A)
Seat depth	425 mm (B)	390 mm (B)
Seat height 1	60 mm (C)	20 / 60 / 100 mm (C)
Seat height 2	100 mm (D)	–
Seat height 3	150 mm (E)	–
Total width with armrests	590 mm (F)	540 mm (D)
Total width without armrests	430 mm (G)	380 mm (E)
Total depth	490 mm (H)	475 mm (F)
Total height	–	–
Total height with armrests	380 mm (I)	230 / 270 / 310 mm (G)
Total height without armrests	180 mm (J)	30 / 70 / 110 mm (H)
Height, seat to armrests	200 mm (K)	200 mm (I)
Width between the armrests	510 mm (L)	450 mm (J)
Width of the splash guard	205 mm (M)	–
Depth of the splash guard	260 mm (N)	–
Width of opening	225 mm (O)	205 mm (K)
Depth of opening	285 mm (P)	275 mm (L)
Length, fit 1	420 mm (Q)	–
Length, fit 2	440 mm (R)	–

More info		
Weight	–	–
Weight with armrests	5.0 kg	4.4 kg (SH:100 mm) 4 kg (SH:60 mm) 3.6 kg (SH:20 mm)
Weight without armrests	2.7 kg	–
Load capacity	120 kg (19 stone)	150 kg (23 stone)
Colour	White	White
Warranty	3 years	3 years

Aquatec 900

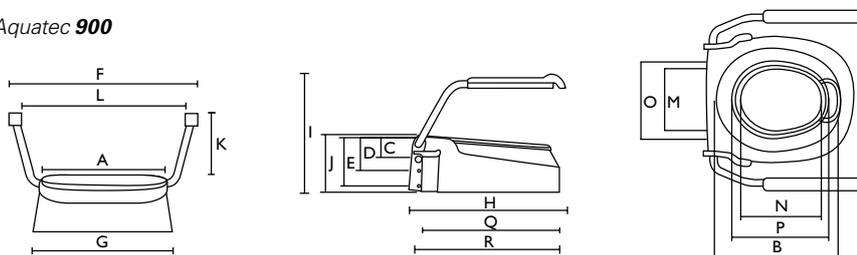

Aquatec 90000

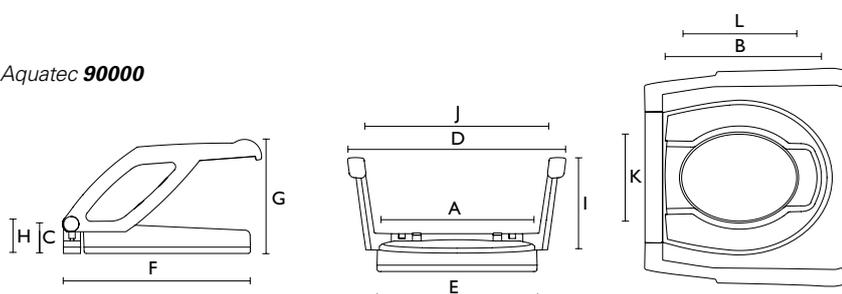

Invacare® / DOLOMITE | **Rollators**





Invacare® / DOLOMITE | **Rollators**

Walking Rollators

<i>Invacare P429/2 Delta</i>	40-41
<i>Invacare P452E/3 Banjo</i>	42-43
DOLOMITE JAZZ	44-45
DOLOMITE LEGACY	46-47
DOLOMITE SYMPHONY	48-49
DOLOMITE MAXI+	50-51
DOLOMITE ALPHA	52-53
DOLOMITE STEP UP	54-57

Invacare® | **Delta**



Lighter and smoother - wherever you want to go!

The new edition of the *Invacare Delta* has been improved to make it lighter and easier to use plus it has some additional features for added convenience.

This lightweight rollator now features a new style of wheel so that its easier to manoeuvre and runs smoothly even over rough surfaces. It has a handy basket, bag and tray; providing more space to carry items.



Yes, you can.®



Every day use

Whether you're out for the day or just popping to the shops, the *Invacare Delta* now comes with a shopping basket and tray providing additional space so you can carry so much more. The basket and bag are large enough to carry a variety of items plus the sturdy tray is practical for use indoors or out.

Sturdy, smooth support for longer walks

With easy to use push-to-lock brakes, the *Delta* is both comfortable to operate and use. It is equipped with new wheels that have been tried and tested to withstand the tough demands of everyday use. In fact our tests* have proven that they are durable and one of the best wheels on the market.

Features



A practical basket to use when out and about.



A sturdy plastic tray to help you around the house.



Tray helps you indoor.

* Tested 200000 cycles using 80% of the user weight.

Technical data



670 mm



600 mm



460 mm



810 - 940 mm

Delta



190 mm



6 kg



135 kg

Delta

Colour



Burgundy

Product Code

Delta 1515533



Invacare® | **Banjo**



A Modern Classic

.....
*Invacare **Banjo*** incorporates all the best features of a classic, quality four wheel walker with a fresh design and additional adjustment options. It has been built to perform well outdoors and fold easily for storage when not in use.





Simplicity and Comfort

The anatomical handle grips are comfortable and offer a secure grip for braking, whilst the height of the handle has an improved range from 710 mm to 980 mm. The safety catch has also been upgraded so that the **Banjo** folds easily through a smart and simple locking system.

Multiple Carrying Options

The *Invacare Banjo* has a hard wearing plastic detachable shelf, a crutch and cane holder, and also spacious basket all fitted as standard.

Features



Tray
Made of impact resistant plastic.



Practical basket
Max load: 9kg.



Crutch holder
Keep a cane or crutch on the rollator.

Technical data

Banjo	600 mm	690 mm	440 mm	710 - 980 mm	610 mm

Banjo	390 mm	190 mm	9 kg	130 kg	260 mm

Colour



Silver

Product Code

Banjo 1452442



DOLOMITE | JAZZ



An innovation in rollator design

DOLOMITE JAZZ combines modern design features with exceptional functionality. The clever cross folding design mechanism means it is quick and simple to fold away for storage or transportation. It offers various height and handle adjustments combined with a built-in brake system to provide stability, comfort and safety for a wide range of individuals.





Easy to operate, simple to store

Made of aluminium, the lightweight **DOLOMITE JAZZ** folds simply into a compact, upright position allowing for easy storage. The spacious basket is made of a durable and flexible material that folds with the rollator. It offers an easy to use brake system that is integrated into the frame ensuring a sleek and stylish design.

Adaptable to individual needs

The **DOLOMITE JAZZ** is available in two heights with easily adjustable handles that can be set in precisely the right position to suit individual needs. The seat has been designed for comfort and a back support comes as standard.

NEW! Slow Down Brake - A safety feature that adds an adjustable resistance to one or both wheels as required.

Features and Options



Comfortable backrest

A comfortable and safe support to lean against. Mounted in the handle bars.



Cane holder

Easy to fit, easy to use. No tools.
Part no 1512912



Basket Cover

Stylish and simple to fit, gives the basket weather protection and privacy. Comes in two colours - red or silver.
Part No 1524818 (Red)
Part No 1524819 (Silver).



Slow-down brake

A safety feature that prevents the walker rolling away. Perfect for gait training and one handed solutions.
Part no 13041-21

Technical data

JAZZ 600	620 mm	670 mm	440 mm	750 - 970 mm	600 mm
JAZZ 520	620 mm	670 mm	440 mm	670 - 810 mm	520 mm
JAZZ 600	390 mm	200 mm	7,8 kg	130 kg	260 mm
JAZZ 520	390 mm	200 mm	7,2 kg	130 kg	260 mm

Colour



Charcoal Grey

Product Code

JAZZ 600 **1507823**
JAZZ 520 **1507775**

Including tray, basket and back-strap.



DOLOMITE | LEGACY



The Original

The **DOLOMITE LEGACY** has everything you would look for in a quality four wheeled walker - reliability, stability and durability.

The **LEGACY** is a simple tried and trusted design that has been the inspiration for many similar walkers over the years.

A backstrap, basket and tray all come as standard on the **LEGACY**.





Open frame design

The spacious walking area within the rollator frame is the signature feature of **DOLOMITE** products. The flip-up seat creates space and makes it easier for the individual to access objects and working surfaces. The simplicity of the rear fork design reduces the possibility of obstruction to the individual's legs and feet. The comfortable, anti-slip seat is structured with a rounded rear edge and generous seating area.

Simple folding mechanism

An easily accessible handle makes the **DOLOMITE LEGACY** simple to fold for storage or transportation.

Features and Options



Backrest

Gives secure seating. Tray and basket helps transporting items.



Under Seat Bag

Hangs under the seat to create additional storage. Part No 12523



Slow Down Brakes

Gives an additional wheel resistance which slows down one or both wheels. Part No 13041-21

Technical data

Legacy 520	600 mm	660 mm	460 mm	770 - 850 mm	540 mm
Legacy 600	600 mm	710 mm	460 mm	840 - 930 mm	620 mm

Legacy 520	400 mm	200 mm	7,1 kg	125 kg
Legacy 600	400 mm	200 mm	7,4 kg	125 kg

Colour



Blue

Product Code

Legacy 450 1507482
Legacy 520 1507372
Legacy 600 1507356



DOLOMITE | *SYMPHONY*



Extra adaptable and lightweight

The **DOLOMITE SYMPHONY** is the lightweight option in our range of rollators. As well as the standard **DOLOMITE** features, the **SYMPHONY** has a lightweight aluminium frame, with a total weight of only 6.1kg - 6.8 kg. It also comes with an extensive range of adaptations for those more specific requirements.

Included as standard on the **SYMPHONY** - backstrap, basket and tray.



Light, stylish and durable

The aluminium frame provides a low product weight, whilst retaining stability and manoeuvrability. Its rustproof properties ensure that the walker retains its good looks come rain or shine.

Easier to use brakes

The long handles provide good leverage on the brakes, making them easier to operate for those with weak hand strength.



Features and Options



Oxygen cylinder holder
Part No. 12682-25



Slow Down Brakes
Gives an additional wheel resistance which slows down one or both wheels.
Part No 13041-21



Under Seat Bag
Hangs under the seat to create additional storage.
Part No 12523

Technical data

Symphony 450	600 mm	590 mm	460 mm	630 - 780 mm	470 mm
Symphony 520	600 mm	620 mm	460 mm	700 - 850 mm	540 mm
Symphony 600	600 mm	700 mm	460 mm	780 - 930 mm	620 mm

Symphony 450	370 mm	200 mm	6,1 kg	100 kg
Symphony 520	370 mm	200 mm	6,6 kg	100 kg
Symphony 600	370 mm	200 mm	6,8 kg	100 kg

Colour



Silver

Product Code

Symphony 450 1511173
Symphony 520 1507267
Symphony 600 1507170



DOLOMITE | *MAXI+*



Stronger and smarter

The **DOLOMITE MAXI+** is a heavy duty walker that has been developed to meet the needs of larger individuals. The generous, reinforced frame has a weight capacity of 200 kg (32st).

A backstrap, basket and tray all come as standard on the **MAXI+**.





Extra space

The **DOLOMITE MAXI+** features a wider seat than conventional walkers; the additional 100 mm of space between the handles provides more comfort.

Adaptable to individual needs

The **DOLOMITE MAXI+** is available in two frame heights with adjustments possible within each frame. It is recommended that individuals utilize the higher seat height when having difficulty moving from a sitting to a standing position. Further adaptation is possible with a selection of accessories from the range.

Features and Options



Comfortable backrest



Slow-down brake

A safety feature preventing the walker from rolling away. Part No 13041-21.



Cane holder

Makes the cane easily accessible without getting in the way. Part No 12606.

Technical data

Maxi 55	690 mm	650 mm	530 mm	720 - 890 mm	560 mm
Maxi 65	690 mm	740 mm	530 mm	810 - 980 mm	650 mm

Maxi 55	470 mm	200 mm	9,6 kg	200 kg	
Maxi 65	470 mm	200 mm	9,9 kg	200 kg	

Colour



Blue

Product Code

Maxi 55 **1511171**
Maxi 65 **1511172**



DOLOMITE | ALPHA



Customised Walking Solutions

The **DOLOMITE ALPHA** is a four wheeled walker with a comprehensive range of options that can provide additional support for differing conditions from arthritis, orthopaedic needs, post-operative care to gait training. In particular the handle and brake selection enables the frame to be adapted to individual requirements.



Infinite adjustable handle

The basic version of **DOLOMITE ALPHA** has the necessary height, width and depth adjustments needed to provide individuals requiring extra support with a good walking position. The rounded armrests offer secure and comfortable support for the forearms. (See picture over)

Advanced Customisation

The **DOLOMITE ALPHA ADVANCED** is unique in the opportunities it presents for multi-directional adjustments afforded by the ball-and-socket design. The angle of the handles, the operating direction of the brakes and height of the brake handle can all be set as required. This configurability can offer pressure management to individuals who may be suffering from high levels of discomfort. (See picture far left)



Features and Options



Platform support

Large rest area with platform support.



Adapt Handle with one-hand-brake



One-hand-brake

Operate both brakes with one hand.

Technical data



Basic	650 mm	700 mm	510 mm	950 - 1240 mm	620 mm
Advanced	650 mm	700 mm	510 mm	950 - 1240 mm	620 mm
Basic Platform	650 mm	700 mm	510 mm	1030 - 1240 mm	620 mm
Advanced Platform	650 mm	700 mm	510 mm	1030 - 1240 mm	620 mm



Basic	430 mm	200 mm	10,9 kg	125 kg
Advanced	430 mm	200 mm	11,4 kg	125 kg
Basic Platform	90 - 270 mm	200 mm	12,4 kg	125 kg
Advanced Platform	100 - 300 mm	200 mm	12,9 kg	125 kg



Colour



Blue

Product Code

ALPHA BASIC	12230-35-21
ALPHA ADVANCED	12220-35-21
ALPHA BASIC with platform	12230-02-35-21
ALPHA ADVANCED with platform	12220-02-35-21



DOLOMITE | *STEP UP*



The Optimum Rehabilitation Choice

The **DOLOMITE STEP UP** is a modern platform walker which maximises gait training techniques for rehabilitation purposes. It gives the patient confidence to start walking again following surgery or recovering from an injury. It is also suited to helping those who have long term problems with their walking. The **STEP UP** is available in three height adjustable versions - manual, hydraulic and electric. Assembly of accessories is tool free.





Ergonomic design

The single central column allows the carer easy access to help the patient and the simple adjustments mean the **STEP UP** can be set to suit specific needs. This means excellent manoeuvrability whilst continuously promoting an upright posture. The base is raised at the front to allow a longer step, while the lowered rear end helps with access under beds, etc. The unique fixed angle of the arm supports give a 25° angle between the shoulder point and wrist, providing the most natural position, maximising upper body strength and keeping elbows close. The platform is tilted 5° giving further support, and the 110° angled hand grip allows for a natural wrist position.

Fully Adjustable


The **STEP UP** is fully adjustable. And enables for the adjustment in the width of the base and arm supports, the height of the central column and the depth and angle of hand grips. The arm pads have been made using a soft material for pressure relief and are securely attached with microfilm. The arm pads are simple to remove for cleaning or replacing.

Transfer and access

The **STEP UP** is well thought out for transfer situations. The lowered rear end helps to reverse under the bed, whilst

the generous width allows wheelchairs to access the frame.

Features and Options

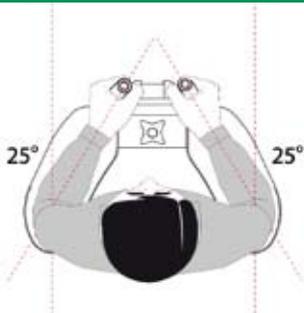

Release the central knob to adjust the width of arm supports.

The width of the legs on the **STEP UP** can be adjusted.

Accessories and options

Quick and easy assembly in track or prepared recess.



Practical Basket
Keeps belongings within easy reach.
Part No 1531390



Tray
To carry those small items.
Part No 1531391



Elbow Stop
Secures elbows from dropping off the side.
Part No 1531448



Swivel Castor with directional brakes
Helps focus on walking instead of direction.
Part No 1531408



Brake System
For those capable of operating a brake function.
Hand-brake kit
Part No 1537656
One-hand brake
Part No 1531450



Oxygen tank holder
Securely mount an oxygen tank on the walker.
Part No 1531406



Separate IV Holder
Separate IV-stand makes life complicated. Not anymore.
Part No 1531407



Stand assist incl. sling
For prescribed standing periods
Part No 1531449 (standing plate + knee support + sling)
Standing plate only
Part No 1535330



Ball grips
Suitable for those who have difficulty closing their hand.

Others

Small castor 100 mm diameter.

Models



Manual
Pull out the pin, move to required height, release the pin.



Hydraulic
Pull up the lever to adjust to required height with hydraulic assistance.



Electric
Using the controller, set the height with the electric motor.

Technical data



Step Up Manual	580/670/760 mm	700/790/880 mm	870 - 1320 mm	340 - 470 mm	620 - 750 mm
Step Up Hydraulic	580/670/760 mm	700/790/880 mm	870 - 1320 mm	340 - 470 mm	620 - 750 mm
Step Up Electric	580/670/760 mm	700/790/880 mm	870 - 1320 mm	340 - 470 mm	620 - 750 mm



Step Up Manual	810 mm	16 kg	150 kg	125 mm	215 mm
Step Up Hydraulic	810 mm	18 kg	150 kg	125 mm	215 mm
Step Up Electric	810 mm	23 kg	150 kg	125 mm	215 mm

Colour



RAL 1035

Article code

Step Up Manual 1521677
Step Up Hydraulic 1521678
Step Up Electric 1521679

Invacare® | **Walking Aids**




 Invacare® | **Walking Aids**

Walking frames	<i>Invacare Asteria</i>	60-61
	<i>Invacare Brio</i>	60-61
	<i>Invacare Aventia</i>	60-61
	<i>Invacare P435B Escort</i>	60-61
	<i>Invacare Foria</i>	60-61
Crutches	<i>Invacare Crutches</i>	62-63
Canes	<i>Invacare Canes</i>	64-65
Paediatric Walking Aids range	<i>Invacare Pluto</i>	66
	<i>Invacare Fun Kids Crutch</i>	67
	<i>Invacare Fixi and Fox Rollator</i>	67
	<i>Invacare Walk-On walking frame</i>	67

Invacare® | *Walking Frames*



Stable, Simple and Safe

The *Invacare* range of **walking frames** cover every option you could need from a walker. This includes fixed frames, folding frames and reciprocal frames. And different types of wheel options are also available.

They are designed to be lightweight, give maximum support, and to fold away easily for storage or transportation.





Invacare Escort

A medium size domestic aluminium frame which is height adjustable. Wheel kit also available.



0435CEB
Fixed Wheels



Invacare Brio

The **Brio** is a dual release aluminium walking frame which is height adjustable and folds easily for storage or transport.



Invacare Asteria

The **Asteria** is a stable, lightweight and user friendly folding walking frame. The central button makes it easy to open and close the frame.



Invacare Aventia

A folding reciprocal frame which moves forward in step with you! The **Aventia** is made of lightweight aluminium and is height adjustable.



Invacare Foria

A height adjustable folding aluminium frame which has handles at two different levels. This will help someone to rise from a seated position to a standing position.

Accessories for Brio and Asteria



6270E
Fixed wheels.



6265E
Swivel wheels with push down brake.



6267E
Swivel wheels.



6015
Pouch - **Brio** only.

Technical data

	Escort	Brio	Asteria	Aventia	Foria
Code	P435B	6291E-A 6291E-JR	1515569	P402	P421B
Max width (mm)	620	600	600	650	570
Max depth (mm)	400	430	450	440	500
Height (mm)	770 - 850	770 - 950 JR 640 - 820	770 - 950	780-940	790 - 920
Weight (kg)	2,5	2,5	2,8	2,5 - 3	2,5
Weight Capacity (kg)	130	136	135	100	130
Color	Grey	Grey	Grey	Grey	Grey



Invacare® | **Crutches**



From standard crutches to special needs

.....
Invacare **Crutches** comes in all shapes in sizes. So whether you are looking for a basic value crutch, an everyday disposable crutch, or something more colourful, we have the solution.

You will also find children's and bariatric crutches in our range.



Invacare Elbow Crutch



Double adjustable **elbow crutch** with ergonomic handle, anti-rattle screw, height adjustment guide. Available in three sizes : Adult, Tall and Junior.

Patient Height guide :
 Adult Crutch : 5'2" – 5'11"
 Tall Crutch : 5'11" – 6'8"
 Junior Crutch : 4'5" – 5'2"

Invacare Magic Twin



Single adjustable open cuff crutch with ergonomic handle and anti-rattle screw. **Magic Twin** crutches will self support by interlocking.

Available in
 Red
 Blue
 Black

Invacare Big 150



Heavy duty crutch with ergonomic handle and height adjustable open cuff. Fixed length which can be cut to suit patient's required measurements.

Invacare Fun Kids



Suited for children, colourfully designed closed cuff crutches which are double adjustable with ergonomic handle.

Invacare Safe In



Stylish double adjustable crutch with ergonomic handle and closed safety cuff. Features anti rattle screw and comes in two colours.

Available in
 Red
 Blue

Invacare Anatom Soft



Stylish double adjustable crutch with soft anatomically shaped handle for extra comfort and closed safety cuff. Features anti rattle screw and comes in two colours.

Available in
 Red
 Blue

Invacare Big Foot



Non sliding heavy duty rubber ferrule with large base and flexible shape to allow flat contact at different angles. Fits 18.5mm – 20.5mm Part No. 1540666

Invacare Crutch Belt

Attach the end of this stretchy belt to your crutches, and you are free to release your arms without putting down the crutches. Part No. 1540664



Technical data

	Elbow Crutch	Magic Twin	Big 150	Fun Kids	Safe In	Anatom Soft
Code	8153-A (Adult) 8153-T (Tall) 8153-J (Junior)	1540650 (Red) 1540651 (Blue) 1540652 (Black)	1529868	1540661	1540657 (Blue) 1540658 (Red)	1540659 (Blue) 1540660 (Red)
Height Floor to Handle	737-965mm (A) 813-1041mm (T) 584-813mm (J)	750-980mm	1080mm (Cut to length)	490-715mm	750-980mm	750-980mm
Height Handle to Cuff	267-343mm (A) 279-356mm (T) 267-343mm (A)	240mm	240- 265mm	175-225mm	240-265mm	240-265mm
Weight Capacity	137kg	130kg	150kg	100kg	120kg	120kg
Colours	Chrome	Red, Blue, Black	Black	Multi	Red, Blue	Red, Blue





Functionality, Style and Quality

Our range of canes and tripods are made to meet the demands and requirements of everyday life.

Designed and manufactured to the highest standards, we aim to provide the latest innovation as well as making sure you get the adjustability and style you have come to expect from Invacare.



Invacare **Prompto**



Adjustable aluminium cane with a fisher style handle. Available in left and right hand versions.

Invacare **Moviale Range**



Chestnut: **252** **253** **257** Beech: **254** **255** **270**

A classical range of wooden canes with different style handles to suit individual tastes. Lengths 900mm and 1000mm, can be cut to required size.

Invacare **Simplex**



A height adjustable aluminium cane with derby style handle.

Invacare **Foldeo**

Height adjustable aluminium cane which folds down to handy size and fits into its own pouch.



Invacare **Auxil**



Aluminium tripod available in two sizes (adult and junior). The **Auxil** is height adjustable.

O300 Cane holder



Create a place to hold a cane.

O400 Caneo



Hang your cane on the edge of a table, from a belt or pocket, etc.

O310 Dragonne



Cane wrist strap.

Technical data

	Prompto	Moviale	Moviale	Simplex	Foldeo	Auxil
Code	P264A (Right) P265A (Left)	P252 P253 P254	P255 P257 P270	1540625	P260	P328 (Adult) P328J (Junior)
Height Floor to Handle	790-1020mm	900mm	1000mm	730-955mm	830-930mm	700-920mm (A) 520-680mm (J)
Weight Capacity	100kg	100kg	100kg	100kg	100kg	100kg (P328) 100kg (P328J)



Invacare® | *Paediatric Walkers*



* Pictured with seat accessory

Invacare Pluto - Children's Reverse Walker

The **Pluto** reverse walker provides independence of movement for children. It's rear support design helps create an upright posture when standing and walking.

The **Pluto** promotes the therapeutic exercises needed to learn natural movements and increase confidence of balance. The automatic braking system in the back wheels can be activated to stop the walker from unintentionally rolling backwards.

It has steerable front wheels and comes in two frame sizes, both height adjustable and with further adjustments possible to the handgrips.

Accessories

1540676 Pluto Seat (Small)

1542470 Pluto Seat (Medium)

1540677 Front wheels with steerlock



Yes, you can.®



Invacare **Fun Kids**

Suited for children, colourfully designed open cuff crutches which are double adjustable with an ergonomic handle.



Invacare **Fixi and Fox**

The **Fixi and Fox** are aluminium rollators designed for children. They have 5 fixed height adjustments and fold down to a convenient size for storage or transportation. The dual action brakes are easy to operate and push to lock. Practical, fun and safe, the **Fixi and Fox** weigh only 5.2kg and 5.3kg respectively.

Accessories

1542469 Seat for **Fixi** or **Fox**.



Invacare **Walk On**

The **Walk On** is a children's aluminium walking frame with non-slipping handles for safer use. The **Walk On** has five fixed height adjustments and easily foldable.

It is available in two sizes, small and extra small, both come with and without wheels as required.

Technical data

	Pluto XS	Pluto S	Fun Kids	Walk On	Walk On Wheeled	Fixi	Fox
Code	1540674	1540675	1540661	1540668 (XS) 1540669 (S)	1540670 (XS) 1540671 (S)	1540672	1540673
Height (mm)	420-725	500-805	–	605 – 690 (XS) 690 – 775 (S)	605 – 690 (XS) 690 – 775 (S)	600-690	690 – 775
Width (mm)	660	680	–	540	540	575	575
Depth (mm)	750	800	–	440 (XS) 500 (S)	440 (XS) 500 (S)	625	670
Floor to handle (mm)	–	–	490-715	–	–	–	–
Handle to Cuff (mm)	–	–	175-225	–	–	–	–
Seat Height (mm)	–	–	–	–	–	425	500
Product Weight	6.2kg	6.6kg	–	2.7 kg (XS) 3.0 kg (S)	2.7 kg (XS) 3.0 kg (S)	5.2 kg	5.3 kg
Weight Capacity	60kg	80kg	100kg	150 kg	150 kg	100 kg	100 kg

Invacare® | **Hoists and Slings**



INVACARE
Yes, you can.®

Invacare® | **Hoists and Slings**

Lifter

<i>Invacare</i> Birdie® Family	70-73
<i>Invacare</i> Jasmine 	74-75
<i>Invacare</i> Roze 	76-77
<i>Invacare</i> Reliant 350 	78-79
<i>Invacare</i> Robin	80-83
<i>Invacare</i> EC-Track	84-87
<i>Invacare</i> Slings	88-95

Invacare® | **Birdie® and Birdie® Compact**



5 year warranty*



Transportability and patient comfort go hand-in-hand

With maximum space for the user, both the *Invacare Birdie* and the *Birdie Compact* offer a truly comfortable lift and transfer to or from beds, chairs or even the floor. The lifters are designed to ensure that folding and unfolding can be carried out easily and without the need for tools. In addition, the lifter can be dismantled into two parts if needed – also without tools.

INVACARE
Yes, you can.®



Total user comfort

The high boom maximizes space for the client even when the hoist is in the highest position. This gives an ergonomic lift and transfer as the client can be easily rotated in the hoist. This makes placing the client into a bed/wheelchair much easier for the carer. The optimization of space creates a reduction in risk of injury as the client's knees will not be close to the actuator of the hoist.

Ease of handling: simple to fold, unfold and dismantle without tools

Handling the lifter has been made as smooth and easy as possible as the hoist can easily be folded and unfolded without the need for tools. The **Birdie** takes up minimal space when stored and is easy to push and transport. If required, the mast and boom can be easily detached from the base (no tools required).

TÜV-certification for maximum security

For total security, the hoists come with a snap-lock system, thus minimizing the risk of accident. The lifters are CE-certified in line with the Medical Device Directive 93/42/ECC. To ensure

the hoists comply with EN ISO 10535 (the European standard for hoists) the hoist has been rigorously tested both by Invacare and the German TÜV institute, an accredited external body.

Features and Options



Long boom – high user comfort

The user can easily be rotated and positioned. No risk of injuries to the knees.



Easy to fold and dismantle

Ease of transport and storage. No tools required when folding / unfolding or when dismantling the hoist.



Snap lock

Optimal security to avoid risk of accidents.



Emergency stop and emergency lowering

Manual or electric emergency lowering (compact version features only electric lowering). Available with fixed or detachable battery. Off board battery charger is available.

Features and Options



Click on spreader bar
Easy to change the spreader bar.



Manual leg spread
Simply step on the pedal to adjust the leg spread.



Handle for leg spread
Handle can easily be attached to operate the manual leg spread.



Electrical leg spread
Available for *Invacare Birdie* (option).



2 point quick release spreader bars
Easy to attach and detach. Available widths: 350 mm, 450 mm, 550 mm.



Push bar and hand control
Superior ergonomics for carers.



Castors
100 mm or 75 mm castors available. (braking castors are shown)



Scale
If required, a scale is available approved. According to the medical standard EN 45501, certified and sealed in order to be used to medicate.



4 point quick release spreader bars
Easy to attach and detach. Available widths: 450 mm and 550 mm.

Technical data

					
Birdie	1150 mm	655 mm	42 kg	465 - 1685 mm	100 mm
Birdie Compact	1000 mm	530 mm	36 kg	515 - 1585 mm	100 mm

					
Birdie	1040 mm	170 kg	40 cycles	2,9 Ah	1400 mm
Birdie Compact	890 mm	150 kg	40 cycles	2,9 Ah	1070 mm

All measurements are taken with 75 mm castors and 2 point spreaderbar.
For 100 mm castors add 15 mm for height measures and 20 mm for width.

Frame colour



Light grey /
Dark grey

*Warranty:

Actuator: 5 years
Frame: 3 years
Battery: 1 year

Invacare® | **Jasmine**

NEW



5 year warranty*



Flexibility and comfort – an ideal solution for Acute and Residential care

The *Invacare Jasmine* Mobile Lifter is a flexible solution that offers a high level of comfort for all clients. A versatile lifter, the *Jasmine* provides safe patient handling up to 32 stone (200kg) weight capacity. With an expansive, wide-opening base and high lifting mechanism, the *Jasmine* offers easy manoeuvrability and positioning for a wide range of lifts.

With the latest control box technology, the *Jasmine* offers a 'Soft Start' and 'Soft Stop' function that provides a more comfortable and dignified transfer.



Yes, you can.®



Features and Options

Flexibility and maneuverability

The *Invacare Jasmine* is exceptionally mobile and easy to transport. The wide-opening base allows easy access to wheelchair or armchair, whilst the high lifting mechanism is essential for transfers to and from elevated surfaces. Flexible, and with optimized maneuverability, the *Jasmine* covers the needs of both the carer and client.

Safety first, with an intelligent control box

The *Jasmine* is fitted with the control box with diagnostic system: an intelligent service lifetime monitoring system that highlights the battery status and service condition of the actuator to ensure safe patient handling.



Spreader bar options

The *Jasmine* can be ordered with 2 or 4-point spreader bars. In widths: 450 mm and 550 mm.



Scale

A scale is available. Approved according to medical standard EN 45501, certified and sealed.



Linak® Jumbo™ Care

The intelligent monitoring system includes a service reminder, emergency operation, audio battery alarm, and battery indicator.

Technical data

Invacare <i>Jasmine</i> 1340 mm	690 mm	1050 mm	540 - 1920 mm	120 mm
Invacare <i>Jasmine</i> 1530 mm	200 kg	43 kg	40 cycles	

Product codes:

- CJ1541963 - Jasmine Lifter - 2 pnt 450mm spreader bar
- CJ1541964 - Jasmine Lifter - 2 pnt 550mm spreader bar
- CJ1541965 - Jasmine Lifter - 4 pnt 450mm spreader bar
- CJ1541966 - Jasmine Lifter - 2 pnt 550mm spreader bar

Frame colours



Moderate white

Blue dusk

***Warranty:**

- Actuator 5 Years
- Chassis 3 Years
- Battery 1 Year



NEW



5 year warranty*



Versatile and safe - an ideal solution for long-term care facilities

The *Invacare Roze* stand assist offers an alternative to a passive lift when a degree of weight bearing ability remains. Ideal for rehabilitation, degenerative conditions and clients whose condition varies from day to day, the *Roze* can provide dignity and a greater degree of independence. With the latest control box technology, the *Roze* offers a 'Soft Start' and 'Soft Stop' function that provides a more comfortable transfer. Combining a range of adjustable features with a lightweight aluminum construction, the *Roze* is equally at home in Residential and Acute Care environments.



Yes, you can.®



Features and Options

A versatile stand assist

The *Invacare Roze* is a versatile stand assist lifter that supports the transfer of individuals who can only bear a minimal amount of their own weight. Ideal for transfers from a seated to seated or seated to standing position, the small footprint of the *Roze* and wide leg opening offers a flexible and comfortable transfer.

Safety first, with an intelligent control box

The *Roze* is fitted with the control box with diagnostic system: an intelligent service lifetime monitoring system that highlights the battery status and service condition of the actuator to ensure safe patient handling.



Height-adjustable leg support

The padded leg support can be adjusted for different user needs.



Removable foot plate

The footplate is horizontal and can be removed for standing/walking a patient during rehabilitation.



Linak® Jumbo™ Care

This intelligent monitoring system includes a service reminder, emergency operation, audio battery alarm, and battery indicator.

Technical data

Product code: 1156019



Invacare Roze 940 - 1610 mm



1050 mm



660 mm



1110 mm



1380 mm



Invacare Roze 120 mm



100-200 lifts per full charge



53 kg



Max. 200 kg

Frame colours

***Warranty:**



Moderate white



Blue dusk

Actuator 5 Years
Chassis 3 Years
Battery 1 Year



Invacare® | **Reliant 350 Stand Assist**



Fresh new look and design improvements for ease of use



For safe and comfortable assisted transfers

The *Invacare Reliant 350 Stand Assist* is a compact and maneuverable mobile lifter that provides safe and comfortable assisted transfers for those with limited mobility or rehabilitation needs. The *Invacare Reliant 350* can be used to transfer a client that cannot transfer unaided, but has some level of weight bearing ability, from one seated position to another. The robust design with a maximum weight limit of 158 kg (25 Stone) allows for additional support whilst transferring, preventing unnecessary strain on both patient and carer.



Yes, you can.®




Adaptable and Easy to Use

Reliant 350 features a multi-adjustable lift boom and a selection of sling mounting points to provide an adjustable and ergonomic lift angle. Large rear and small front castors help maneuver the **Reliant 350** into small or confined spaces and the base has a lever adjustment mechanism for added stability. In addition, it has a generous non slip footplate providing a stable surface for the feet, as well as a height-adjustable padded knee guard.

Safe and Reliable


The **Reliant 350** provides a range of safety features including an emergency stop button, emergency mechanical descent system, as well as an anti-entrapment device which stops the boom from lowering if it meets resistance. The battery recharge indicator is visual and audible. For transport and storage, the **Reliant 350** can be easily dismantled.

Features and Options

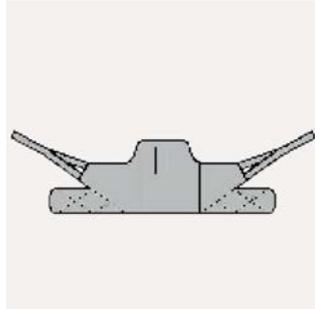

Stand Assist Sling

Designed to be used with a **Stand Assist** lifter, this sling features the latest innovative fabric for improved comfort and handling. Available in sizes; S, M, L & XL.

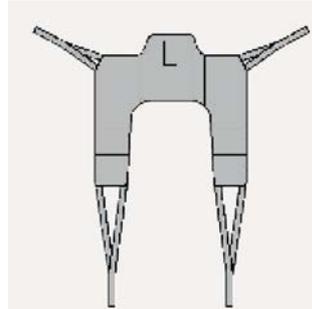

Transfer Stand Assist Sling


Designed for **Stand Assist** lifters, for those who require additional leg support, this sling provides a supportive and comfortable transfer. Available in sizes; S, M, L & XL.



Independent Charger


An off-board charger and spare battery allow for 24 hour use.

Technical data

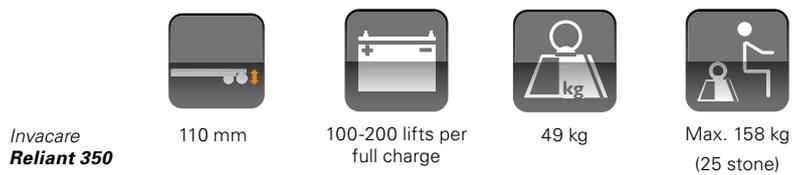

Frame colour

Product Code: RPS350-1E

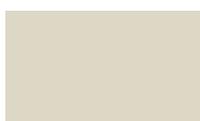

Beige



Invacare® | **Robin**®



3 year warranty



Moving forward with dignity and care

The *Invacare Robin* ceiling hoist provides a unique and innovative method of transfer. The two strap design has no obstructive spreader bar; offering a safe and dignified transfer with optimal client and carer contact.





Comfort

The two strap design gives unsurpassed stability during the transfer. The smooth soft start function makes positioning more gentle, and the two speed motor reduces the amount of time required for the transfer. The strap only configuration eliminates the need for a spreader bar of any size, adapting automatically to the weight and shape of the client.

Transfer with dignity


Personal contact is important. The absence of a spreader bar provides the client with a sense of freedom reducing the fear of being close to obtrusive parts of the lifting equipment. This enhances the transfer experience by introducing a unique dignified method of moving people.

High safety and easy to work with

The *Invacare Robin* has been conceived with safety in mind. In the event of an emergency situation, lowering, raising and stop functions are standard features built into the motor unit. The lifting mechanism offers an overload security

system, self adjusting strap inlet and centrifugal braking. Our patented strap system enables the carer to attach the sling to the hoist with the minimum of effort and the added security of a simple locking mechanism.

Features and Options

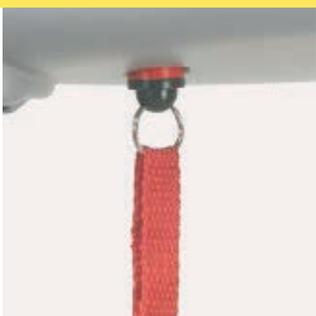

Emergency stop and lowering



Strap inlet self adjusts automatically to size of client



Hook self locks when weighted



Cable disengages to avoid damage if unit is pulled along via hand control

Features and Options



Robin
The unique strap system without the use of a spreader bar.



Robin Mover
With an integrated traverse motor.



Blue light illuminates when motor is in use, and a signal light flashes when battery is low



Charge motor through hand control when not in use



External emergency battery



Control- buttons on motor unit
For emergency use.



A scale which is easy to mount and use



Flexiscope
If a specific situation requires a spreader bar.



Stretcher, that offers unique stability without a leveller



① **Hand control**
Identified by symbols.
② **Robin Mover hand control**
With symbols for up/down and direction.

Technical data

					
Robin	470 mm	250 mm	190 mm	2,5 meter	30 mm/s (85 kg), 40 mm/s (0 kg)
Robin Mover	470 mm	250 mm	310 mm	2,5 meter	30 mm/s (85 kg), 40 mm/s (0 kg)

		
Robin	7 kg	200 kg (32 stone)
Robin Mover	9,5 kg	200 kg (32 stone)

Colour



white/grey

Warranty: 3 years

except on battery and parts of normal wear and tear (e.g. hoist tape)

Lifts per charge; Aprox 40. (user weight 85 kg, 500 mm up/down)

Voltage between battery/hoist; 24V DC

Electrical connection of charger: 230V AC 50 Hz

Invacare® | **EC-Track™**




Discrete, flexible and safe

.....
Invacare EC Track system completes the versatility of the *Invacare Robin*. To meet the varied transfer needs, we offer a wide range of tracking options for different kind of environments and types of installation. The **EC Track** has a discrete modern design with a low height and long span, ensuring that the complexity of the various transfer needs are met. **EC Track** is constructed from light aluminium material and is easy to install.




Yes, you can.®

Features and Options

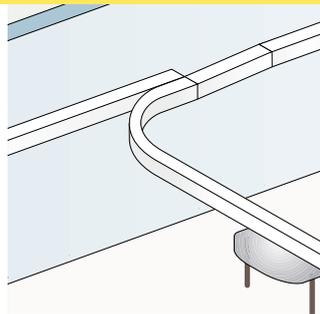


Track switch

Makes it possible to change direction between two tracks.

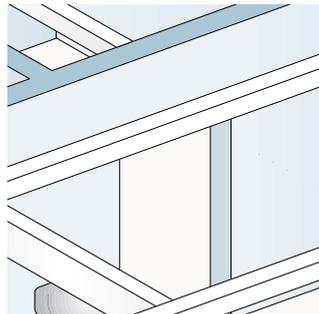


Room to room transfer

With two *Invacare Robin* motors and a room-to-room trolley, it is possible to transfer the user from room to room.

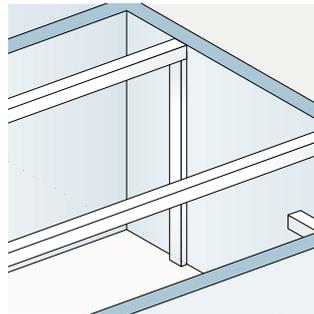



Wall support

The tracks can be installed in places where the wall does not have sufficient carrying capacity.

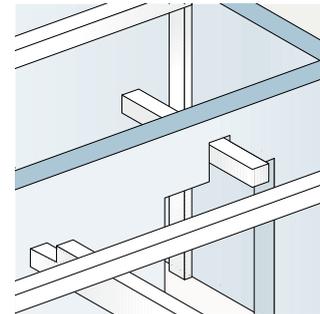


Transit coupling

To transfer the user from one room to another, using only one *Robin* motor.

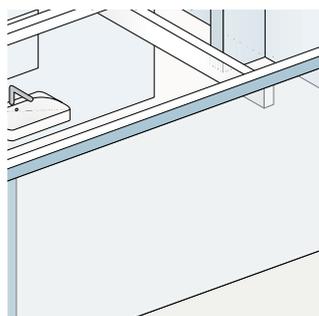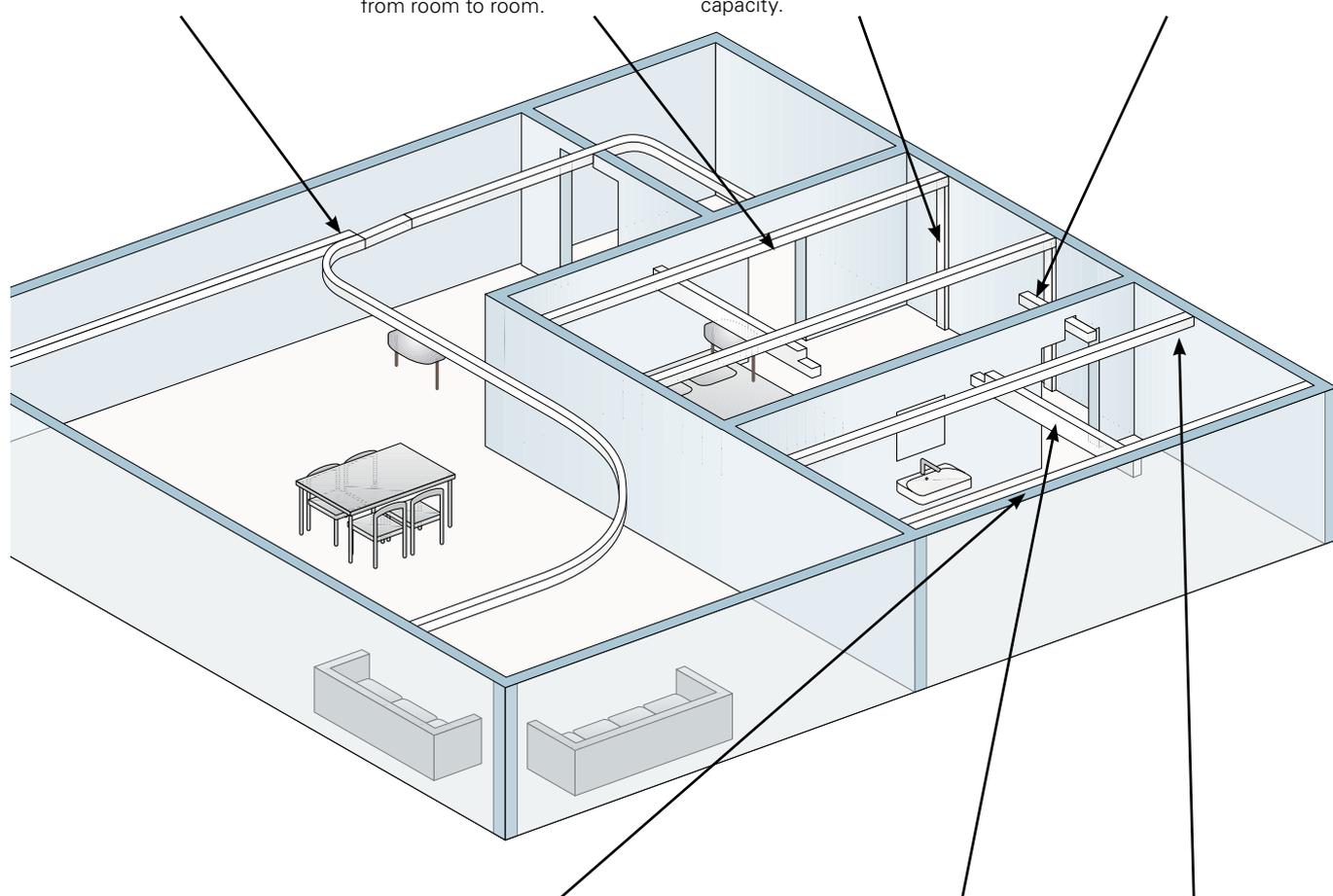


Wall bracket

To install the track along the wall.

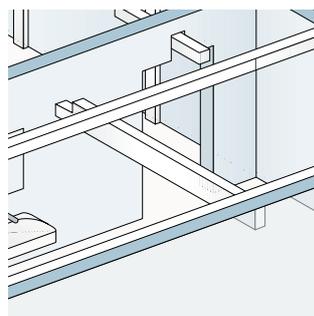

Traverse systems

Permits an unlimited number of lifting places.

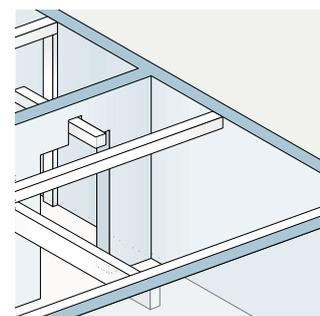

Wall bracket

To install the track in a free span. The tracks are mounted at the ends.

Features and Options



Invacare EC Track
Is available in 3 different profiles.



Curves
Are available in 30, 45, 60 and 90 degrees as standard. Only S-Profiles.



The **Invacare Robin** trolley



Invacare Robin
Connects to the track, without the use of tools.



Standard bracket
Available for all track options.



The Quick bracket
Single point installation (unavailable in the UK).



Wall bracket
An alternative to ceiling mounting.



Wall support
To support the wall-mounted track



Traverse trolley



Track switch
When a change in direction is required during movement of the patient.



Transit coupling
To transfer from one room to another, using only one **Robin** motor.



EC - Gantry

Tracks on „legs“- no need for installing tracks in wall or ceiling. EC - Gantry is very easy to install and dismount and once mounted it is very easy to raise and lower the gantry without using tools.

The EC - Gantry is available in two versions:

1. With castors - height adjustable from 1950 - 2300 mm
2. Without castors - height adjustable from 2160 - 2510 mm



EC -Gantry Traverse

Room covering tracks on „legs“- no need for installing tracks in wall or ceiling.

EC - Gantry Traverse is very easy to install and dismount and once mounted it is very easy to raise and lower the gantry without using tools.

EC- Gantry is height adjustable from 2180 to 2530 mm.

Technical data




Invacare **EC-Track** 200 kg

Frame colour



white

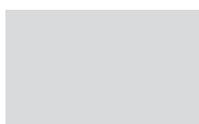

aluminium

Invacare® | **Slings**




Optimal positioning and ease of application

Invacare's new range of slings provides a safe and comfortable transfer solution for many types of physical disability. This selection of high quality slings has been developed with the collaboration of Occupational and Physiotherapists, to ensure optimal fit, support for the client and ease of application for the carer.

Most slings are available in five sizes with the option of fabric types. A range of adjustments enable the slings to meet the specific needs of the individual.



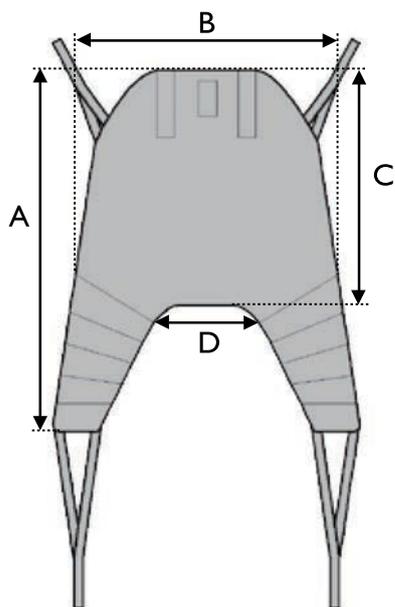


Invacare Universal High Sling

Supports the whole body and the head giving the client a comfortable and slightly reclined position.

The **Universal High Sling** is typically used for those who have limited control of the head, trunk or hips – or for transfers to and from a lying position. The **Universal High Sling** can also be used for seated to seated transfers. This sling can also be used without the head support, by simply removing the plastic rigidisers.

Size	A	B	C	D	Product number / Material		
					Polyester fabric	Net fabric	Spacer fabric
XS	1040	740	610	340	1483776	1485950	1485955
S	1080	840	730	380	1483777	1485951	1485956
M	1420	1000	830	400	1483778	1485952	1485957
L	1550	1100	920	510	1483779	1485953	1485958
XL	1680	1150	920	550	1485949	1485954	1485959

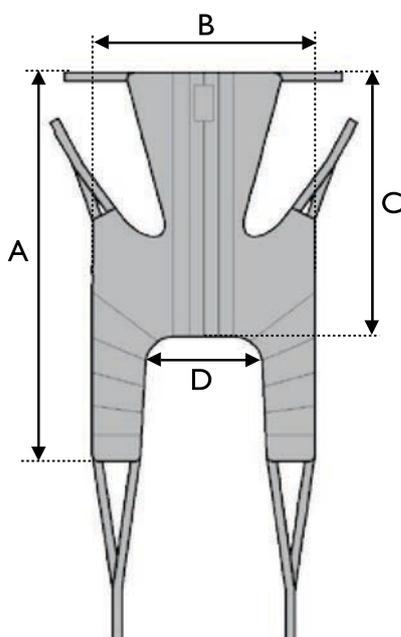


Invacare Universal High Plus Sling

Supports the whole body including the head giving the client a comfortable and slightly reclined position. The **Universal High Plus sling** also allows clients keep their arms outside the sling. It is typically used for those who have limited control of the head, trunk or hips – or for transfers to and from a lying position. The **Universal High sling** can also be used for seated to seated transfers. This model has a stiffened fabric in the back section, to ensure optimal weight distribution and comfort for the client.




Size	A	B	C	D	Product number / Material
					Polyester fabric
S	1380	840	880	340	1485960
M	1510	950	930	450	1485961
L	1580	1100	970	560	1485962
XL	1590	1190	970	630	1485963

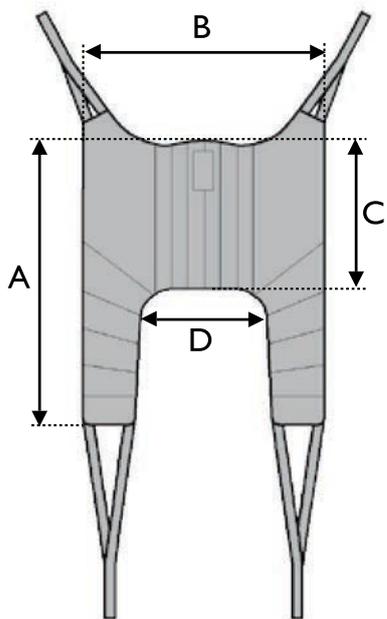

Invacare **Universal Standard Sling**

Supports the whole body; including the shoulders. This sling should only be used for those who have good head control. The **Universal Standard sling** is often preferred as an all-round sling for the elderly, as it supports the entire body and gives the client a feeling of comfort as well as an unobstructed view. This sling can be used for seated to seated transfers and is easy to apply. This model also has a stiffened fabric in the back section, to ensure optimal weight distribution and comfort for the client.


Size	A	B	C	D	Product number / Material		
					Polyester fabric	Net fabric	Spacer fabric
S	1000	920	510	410	1485882	1478411	1485363
M	1090	1050	560	480	1485883	1478412	1485364
L	1180	1090	580	560	1485884	1478413	1485365
XL	1240	1250	585	710	1485885	1478414	1485366

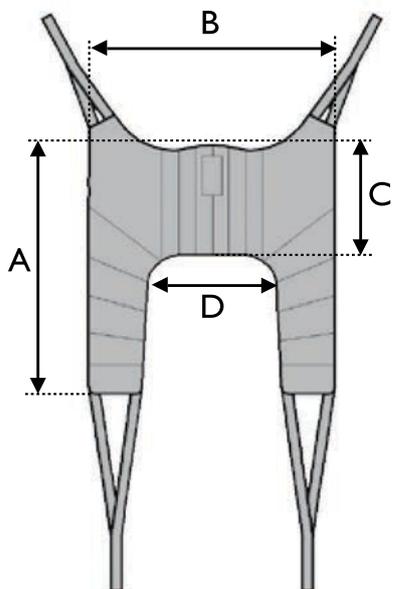

Invacare **Universal Low Sling**

Supports the hips and lower back. The **Universal Low sling** can be used for those who have some trunk control and full head control. It is used for seated to seated transfers and is easy to apply. This sling also has stiffened back section to ensure optimal weight distribution and comfort for the client.


Size	A	B	C	D	Product number / Material
					Polyester fabric
S	850	800	350	410	1491893
M	920	940	375	480	1491894
L	950	1075	405	575	1451067

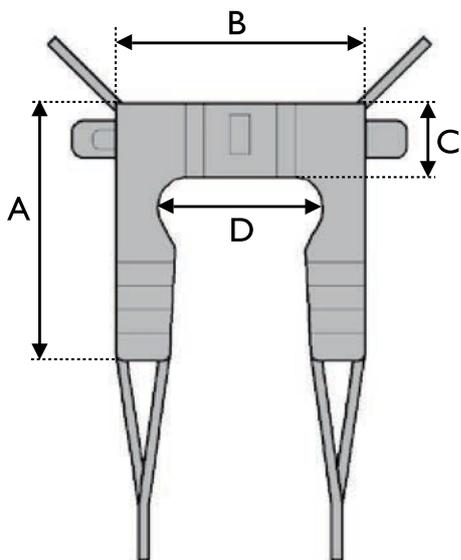

Invacare Dress Toileting Sling

Within the Invacare sling range there are two toilet transfer slings available; with or without head support. Both models come with a belt for increased support around the body and for safe positioning. The **Dress Toileting Sling** is suitable for those who have limited mobility. The model with head support is designed for clients who also have reduced head control. The **Dress Toileting Sling** is suitable for those who want to be able to dress and undress during the toilet transfer. It can be used for seated to seated transfers and is very easy to apply.

Size	A	B	C	D	Product number / Material	
					Toileting polyester fabric	Toileting High polyester fabric
S	870	710	200	430	1478440	1478444
M	950	800	230	490	1478441	1478445
L	1020	900	230	600	1478442	1478446
XL	1100	1080	270	760	1478443	1484780

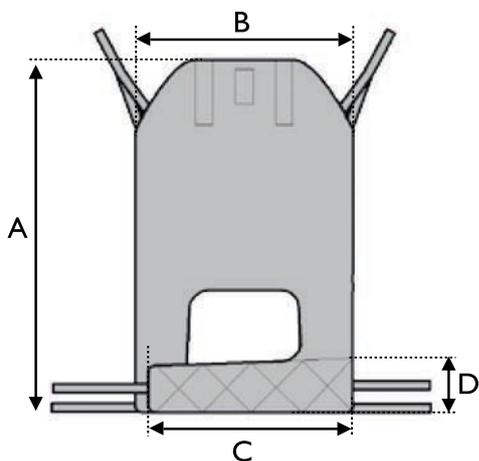


Invacare Easy-fit Sling

Supports the whole body and the head, and is therefore suited for clients with the most challenging physical disabilities*. The specially designed separate leg supports are easy to apply, and the sling affords a good hip position, which allows the client to sit relatively upright. The **Easy-fit sling** is also well-suited for amputees. It is easy to apply in both sitting and lying positions.


Size	A	B	C	D	Product number / Material	
					Polyester fabric	Spacer fabric
XS	1100	710	595	170	1451101	1451106
S	1130	760	640	180	1451102	1451107
M	1200	800	675	190	1451103	1451108
L	1270	860	720	210	1451104	1451109
XL	1390	940	780	220	1451105	1451110

* Following risk assessment

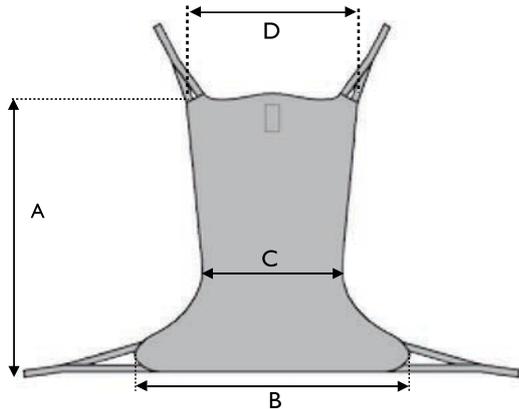


Invacare Comfort Sling, Comfort Sling with Toilet Opening and Comfort Standard Sling

This sling supports the whole body and the head. The **Invacare Comfort Sling** is ideally suited to those who have chronic pain. The sling is also useful for above knee amputees, or for those who due to postural position or pain prefer to remain in a seated position. The **Comfort Sling** is applied while the client is lying down, and is used to transfer the client into a seated position. The **Comfort Sling with Toilet Opening** is applied in the same way and is used for toileting.

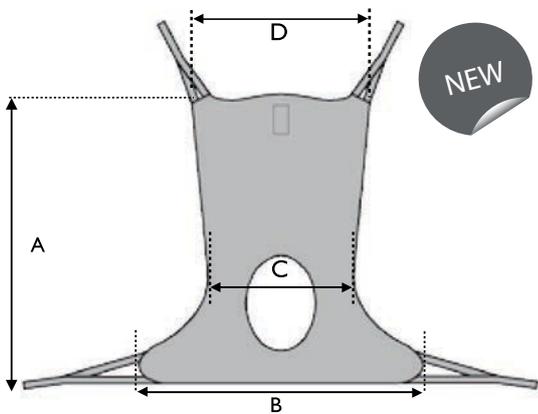

Invacare **Comfort Sling** with and without Toileting Opening


Size	A	B	C	D	Product number / Material			
					Comfort		Comfort Toileting	
					Net fabric	Spacer fabric	Polyester fabric	Net fabric
XS	1040	970	510	960	1485967	1485948	1451091	1451095
S	1150	1050	550	1020	1485968	1485554	1451092	1451097
M	1395	1170	610	1160	1485969	1485777	1451093	1451098
L	1445	1290	675	1280	1485970	1485778	1451094	1451099
XL	1460	1330	745	1470	1485970	1451090	1451095	1451100

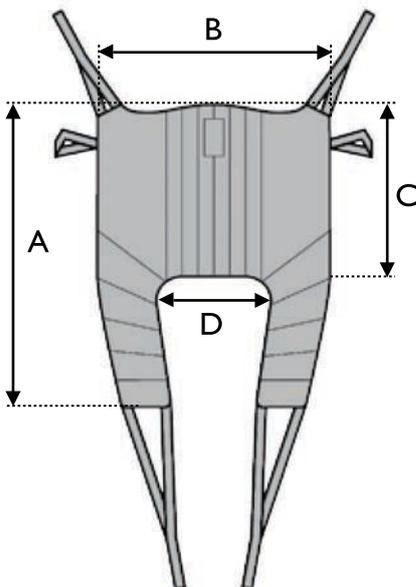

Invacare Amputee Sling

This model is specially designed for the safe transfer of clients with single or double amputations. The sling affords the client comfortable support in a slightly reclined position. It is important that the client has good head control.

Size	A	B	C	D	Product number / Material
					Polyester fabric
S	1100	930	490	320	1485367
M	1160	970	520	360	1485368
L	1360	1230	570	430	1485369
XL	1400	1350	570	480	1485370

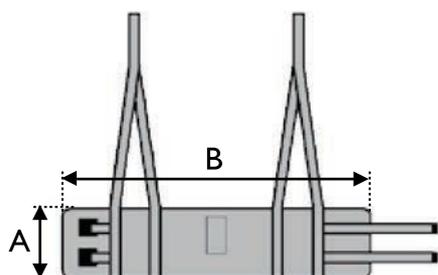


Invacare Standing Transfer Vest

The **Standing Transfer Vest** is applied around the trunk. The client's feet must be on the floor and the client must also have some standing ability. The **Standing Transfer Vest** is used for those with reduced balance and/or walking ability. The vest is especially useful for toileting. It also reduces the risk of injury to the carer.

A groin support is available as an accessory to the **Standing Transfer Vest**. This helps keep the vest in place.

Invacare Standing Transfer Vest

Size	A	B	Product number / Material	
			Polyester fabric	Groin Band
XS	150	750	1485404	–
S	150	850	1485405	1492662
M	180	950	1485406	–
L	180	1050	1484683	1492663
XL	180	1200	1484684	–

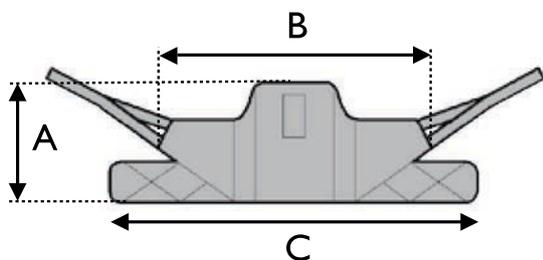


Invacare Stand Assist Sling

Designed for use with the **Reliant 350 lifter**, but compatible with other types of **Stand Assist lifter***, the **Invacare Stand Assist Sling** offers a comfortable and safe transfer for clients who have some weight bearing ability, but cannot stand independently. In addition, the client should be coherent, have good body and head control and the ability to grip the lift boom during the transfer. The sling is well padded and comes with non-slip material on the back to ensure a comfortable transfer position.

Size	A	B	C	Product number / Material
				Polyester fabric
S	260	840	950	1484685
M	330	895	1170	1484686
L	360	950	1270	1484687
XL	360	1020	1810	1484688

May be suitable for use with **Stand Assist Lifters** that are compatible with fabric loop slings. Following risk assessment.

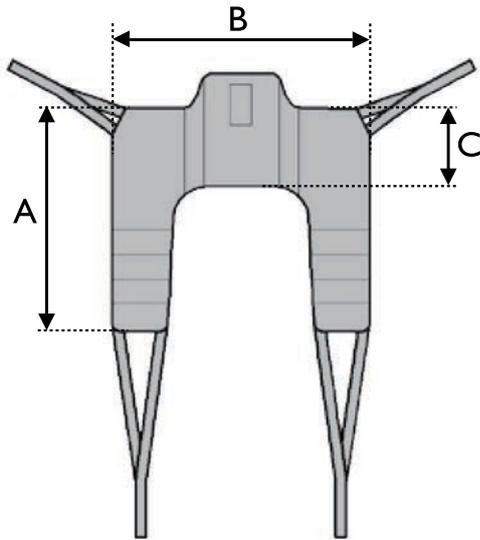

Invacare Transfer Stand Assist Sling

The *Invacare Transfer Stand Assist Sling* can be used with the **Reliant 350** lifter and enables a simple seated transfer. This sling is used for those with good body and head control who are able to grip the lift boom.

Size	A	B	C	Product number / Material
				Polyester fabric
S	860	840	240	1484689
M	1020	895	260	1484690
L	1180	950	330	1484691
XL	1230	1020	350	1485117

Sling Accessories



Fleece padding for extra comfort



Leg insert, Foam. For **Universal Standard, Universal High** and **High Plus** in Net Material



Extension loops - in order to give more positioning possibilities



Head support stiffeners - for Slings with head support. **Easy Fit, Universal High**

Sling Types	Fabric Matrix			Size Matrix				
	Solid	Net	Spacer	XS	S	M	L	XL
Universal Low	●				●	●	●	
Universal Standard	●	●	●		●	●	●	●
Universal High	●	●	●	●	●	●	●	●
Universal High Plus	●	● <small>NEW</small>	● <small>NEW</small>		●	●	●	● ¹
Dress Toileting	●				●	●	●	●
Dress Toileting High	●				●	●	●	●
Easy-fit	●		●	●	●	●	●	●
Amputee	●				●	●	●	●
Comfort		●	●	●	●	●	●	●
Comfort Toilet	●	●		●	●	●	●	●
Stand Assist	●				●	●	●	●
Transfer Stand Assist	●				●	●	●	●

^{1/} Only available in solid material

Fabrics



Solid Polyester

A standard material that is easy to apply – slides smoothly and dries quickly after washing.



Net*

Used for bathing or if the user will stay seated on the sling for a period of time.




Spacer*

A slightly elastic fabric that adapts to the user's body shape. Well suited for those with skin integrity issues.

* It is important to carry out a risk assessment with reference to skin integrity for those who be seated on a sling for long periods.

Size guide

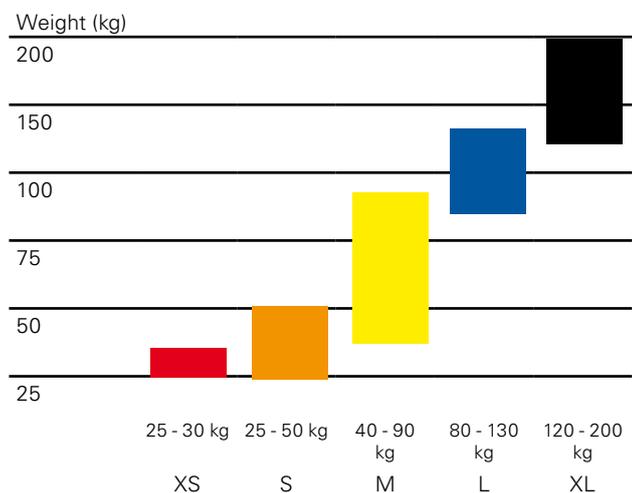

1 year warranty



200 kg**
(32 stone)



95°



All slings should be tested and adjusted to the individual user to ensure optimal comfort and safety during transfers.

** Except for stand assist, Transfer stand assist & Standing transfer vest in combination with groin band: the safe working load is 159 kg (25.5 stone).



Invacare® | **Beds**



Invacare® | **Beds**

Beds	Invacare Etude® Basic	98-101
	Invacare Etude® Medley Family	102-105
	Invacare Etude® Classic	106-109
	Invacare SB 755	110-115
	Invacare Octave®	116-117
Beds Accessories	Invacare Beds Accessories	118-127

Invacare® | *Etude® Basic*



Warranty
5 year actuators
3 year bed frame

The Tried and Trusted Home Care Bed

The *Invacare Etude Basic* is a high quality, height adjustable, electric profiling bed that is designed to meet the needs of the modern Community environment. Complete with a four section profiling mattress platform and combined with a dual height range of 330 - 730 or 400 - 800 mm, the *Etude Basic* provides increased comfort for the client and an enhanced working environment for the carer.



Yes, you can.®



Dignity & Independence

The *Invacare Etude's* extensive height range makes getting in and out of bed easier. The bed is supplied complete with cotsides as standard and can be fully adjusted to suit patient needs using the ergonomic hand control, which features light touch push buttons. A range of accessories are available such as support handles to provide additional assistance when moving in and around the bed. In addition an optional 'Soft Control' can be used to synchronise both the head section and knee break.

Easy Handling & Storage

Portability and storage are key benefits of the *Etude*. Its simplistic design allows easy assembly and dismantling without tools and the process can be safely and comfortably managed by one person (each section weighs less than 25 kg). The bed ends can be easily attached to the transport unit making the *Etude* lightweight and manoeuvrable for handling and storage.

Approved Strength & Stability

Probably one of the most robust Community beds on the market, the *Etude* satisfies all the requirements laid down by the European Standard for Nursing Beds EN1970 and carries the official blue stamp from TÜV (German Safety Standards Institute) for its strength, safety and stability. External

approval is complemented by Invacare's stringent demand for consistent high quality throughout. Internal testing, quality control and a continuous improvement environment ensures excellence in the *Etude's* functionality, quality and safety. Invacare is ISO 9001 and ISO 13485 accredited.

Features and Options



Invacare Scala Basic 2
Foldable steel side rail. Compliant with the new standard for beds IEC 60601-2-52.



Invacare Scala Basic Plus 2
Foldable steel side rail with plastic inserts at the mounting brackets. Compliant with the new standard for beds IEC 60601-2-52.



Invacare Scala Medium 2
Foldable steel side rail for higher mattresses. Compliant with the new standard for beds IEC 60601-2-52.



Invacare Scala Decubi 2
Foldable steel side rail for extra high mattresses. Compliant with the new standard for beds IEC 60601-2-52.

Features and Options



Transport unit

For easy storage and enhanced portability. The **Etude Basic** can be delivered fitted to a transporter kit. (Available on request)



Mattress Extension Support

Allows the mattress platform to be extended by 100 mm or 200 mm to suit many size requirements.



Double Height Positions on Bed Ends

Designed to suit the carers environment and equipment needs, the built-in, double-brackets offer the choice of height ranges of 330 - 730 mm or 400 - 800 mm.




Safe Support in Bed

Lifting pole with safe, adjustable grip handle. Can be attached to either left or right side.



Support Handles

Extra help for the patient getting in and out of bed. Available in 5 different sizes.

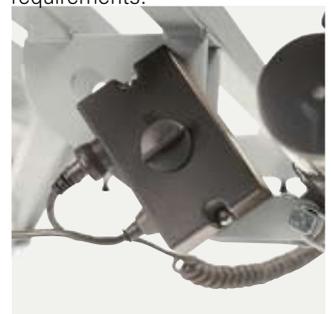

ACP-Box

By using the ACP-box, all handset functions can be blocked. This ensures that clients cannot inadvertently operate the bed whilst the carer is not present.



Well-Fitted Cable Network

Cables are designed to be hung on the bed frame - as a result, there are no loose cables or opportunities for cables to get tangled or trapped under the bed. N.B. Fitted with UK plug as std.



Curled cables for easy handling

The use of curled cables helps to avoid leads becoming trapped within sections whilst profiling.



Easy Maintenance Motors

Motors can be dis-engaged without the need for tools; cabling between the control unit and the motor is easy to change.

Please see *Invacare Bed Accessories* brochure for more accessories.

Technical data

Etude Basic	Width outside: 920 mm Width inside: 900 mm	Length outside: 2220 mm Length inside: 2080 mm Matress: 2000 mm	Height adjustable: 330-730 or 400-800 mm	Matress support dimensions: 4-sectioned: 850-240-350-560 mm	0 - 70° Backrest angle	0 - 18° Thigh angle

Etude Basic	0 - 25° Legrest lift	11° Tilt option*	4-sectioned: 78 kg	SWL 180 kg / 28.5 stone Max. user weight 145 kg /23 stone

* Please notice that the tilt function is only practised for users under medical supervision.

NB/ The bed must not be used by patients under 12 years of age, or by patients with body size equivalent to an average 12 year old or smaller.

For mattress compatibility, please refer to the user manual which can be found at www.invacare.co.uk

Frame colour and wood work



Light Grey

Beech

Products, Options and Features can be subject to change.

Invacare® | *Etude® Medley Range*



3 year warranty



Cost effective patient care

The **Etude Medley** from Invacare offers the perfect balance of budget, functionality and quality requirements. The four-sectioned profiling mattress support allows for a range of height combinations, providing an ideal position for the carer and comfort for the user.

The latest addition to the **Etude Medley** range is the **Etude Medley Low** which lowers to a 220 mm height to maximise safety for vulnerable clients who are prone to unexpected bed leaving. A specially designed bed end fascia then enables both the **Etude Medley** and **Etude Medley Low** to be transformed into the **Invacare Etude Medley Select**. This stylish bed is available in a beech finish that gives the type of homely feel that will integrate into most community and Residential Care environments. All beds are easy to transport and assemble.



Yes, you can.®

Improved functionality and efficiency

The *Invacare Etude Medley* has been specifically designed for quick transportation, to meet the needs of the transient working environments of Community Care provision. The **Medley** can be supplied on a transport kit so that the bed can easily be moved and installed. In addition, any assembling can be done without the use of tools or technicians.

Height is no longer an issue

The flexible design of the **Etude Medley** allows the bed to be mounted at a low height or a higher position, creating a range of working heights for the carer and transfer heights for the patient. This is facilitated by the unique feature of the bed end, which enables the mattress support to be mounted at both upper and lower positions.



Options and features



Double height positions on bed ends

Adjustable height range of 330 - 730 mm or 400 - 800 mm are available for the **Etude Medley**.



Four section profiling mattress support

Etude Medley comes with a full electric mattress support platform to suit individual patient requirements.



Full length side rails

The **Etude Medley** is available with wooden or steel full length side rails. Three quarter length metal side rails are also available in a range of heights.

Technical data

Etude Medley Width outside: 1020 mm Width inside: 900 mm	Length outside: 2220 mm Length inside: 2080 mm (+ 200 mm extensions)	Height adjustable: 330 - 730 or 400 - 800 mm	Mattress support dimensions 4 sectioned: 750 - 350 - 340 - 560 mm	0 - 70° Backrest angle	Thigh angle: 0 - 29° Legrest lift: 0 - 20°

Etude Medley 11° Tilt option*	69 kg	SWL 180 kg (28.5 stone)

* Please notice that the tilt function is only practised for users under medical supervision.

NB/ The bed must not be used by patients under 12 years of age, or by patients with body size equivalent to an average 12 year old or smaller.

Frame colour and wood work



Light Grey



Beech

For mattress compatibility, please refer to the user manual which can be found at www.invacare.co.uk



Options and features

Low height to maximise safety

The height of the **Etude Medley Low** bed offers increased security for clients enabling them to sleep with real peace of mind. The **Etude Medley Low** sits to the floor, (220 mm in height) which combined with a floor mattress, means that if a fall does take place any harm is reduced. Carers also have the option to raise the bed to a more suitable height, when necessary.

The **Etude Medley Low** makes it easier for clients to get in and out of bed, providing more independence.

The **Medley Low** is dual height adjustable, offering a 220 - 620 and 290 - 690 mm height option.



Easy to transport

The **Etude Medley** and **Etude Medley Low** are easy to assemble and dismantle and can be transported using transporter kits.



Floor mattress helps to prevent injuries

In the event of a client falling from the bed, a foldable mattress can be placed on the floor and attached to the **Medley Low** for added safety.

The part number for the foldable mattress is 1494608.

Technical data

Etude Medley Low




Width outside: 1020 mm
Width inside: 900 mm



Length outside: 2220 mm
Length inside: 2080 mm
(+ 200 mm extensions)



Height adjustable: 220 - 620 or 290 - 690 mm

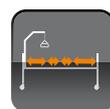

Mattress support dimensions 4 sectioned: 750 - 350 - 340 - 560 mm



0 - 70° Backrest angle



Thigh angle: 0 - 29°
Legrest lift: 0 - 20°

Etude Medley Low



11° Tilt option*



69 kg



SWL 180 kg (28.5 stone)

* Please notice that the tilt function is only practised for users under medical supervision.

NB/ The bed must not be used by patients under 12 years of age, or by patients with body size equivalent to an average 12 year old or smaller.

Frame colour and wood work



Light Grey



Beech

For mattress compatibility, please refer to the user manual which can be found at www.invacare.co.uk

Homely design

The **Etude Medley Select** offers a more attractive, homely design for **Etude Medley** and **Etude Medley Low**. Also available as a retro-fit, this aesthetic bed end cover fits with precision onto the bed actuators and can be removed at any time.

The leg supports are equipped with a Rastofix to enable individual positioning for the legs. The **Etude Medley Select** can be dismantled and stored easily using a transport adapter for moving and storage.




Options and features



Select design for Medley Low
The homely design of **Medley Select** is also available with the **Medley Low** model.



Lifting pole
Can be adjusted in both height and depth.

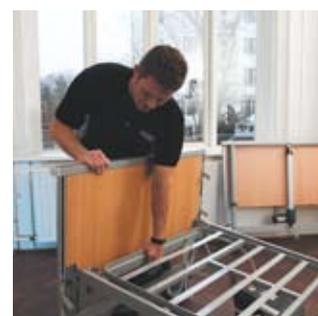

Medley can easily be dismantled by one person without using tools.

Technical data

Etude Medley Select	Width outside: 1050 mm Width inside: 900 mm	Length outside: 2250 mm Length inside: 2080 mm (+ 200 mm extensions)	Height adjustable: 330 - 730 or 400 - 800 mm Medley Low: 220 - 620 or 290 - 690 mm	Mattress support dimensions 4 sectioned: 750 - 350 - 340 - 560 mm	0 - 70° Backrest angle	Thigh angle: 0 - 29° Legrest lift: 0 - 20°

Etude Medley Select	11° Tilt option*	76 kg	SWL 180 kg (28.5 stone)

* Please notice that the tilt function is only practised for users under medical supervision.

NB/ The bed must not be used by patients under 12 years of age, or by patients with body size equivalent to an average 12 year old or smaller.

Frame colour and wood work



Light Grey

Beech

For mattress compatibility, please refer to the user manual which can be found at www.invacare.co.uk



Invacare® | *Etude® Classic - Extra Low*



Warranty
5 year actuators
3 year bed frame

Exceptionally Low - for safe nursing care

Based on the principles of comfort and safety, the *Invacare Etude Classic - Extra Low* boasts a solid, highly stable design. This versatile Community and Residential Care bed excels with an impressive carrying capacity of 220 kg (34.5 stone) and a height adjustment range of 200 - 600 mm to offer increased client safety when transferring whilst providing an enhanced environment for manual handling.



Yes, you can.®



Dignity & Independence

The *Invacare Etude Classic - Extra Low's* extensive lifting range makes getting in and out of bed easier. Not only does it empower the client in and around the bed, it also ensures optimal conditions for the provision of care. With professional assessment, vulnerable clients can be nursed at a low height, without the need for side rails*, which may make this bed suitable for those who leave their beds at night. The extensive lifting range and other flexible accessories makes day to day care significantly easier for everyone.

Easy Handling & Storage

Portability and compact storage constitute key features of the *Etude Classic - Extra Low's* design. In one simple lifting action the mattress support comes away from the bed ends and can be separated manually by releasing a simple mechanism. Having attached the bed ends to the transport unit, the *Etude* becomes lightweight (section weights less than 25 kg) and portable, letting the castors take the strain.

Quality, Strength and Stability

The *Etude Classic - Extra Low* satisfies all the requirements of the European Standard for Nursing Beds EN1970 and carries the official blue stamp from the TÜV (German Safety Standards Institute) for its strength, safety and stability. External approval is complemented by Invacare 's stringent demand for consistent high quality throughout the

manufacturing process. An internal testing, quality control and a continuous improvement environment ensure excellence in the *Etude Classic - Extra Low's* functionality, quality and safety. Invacare is ISO9001 and ISO13485 accredited.

* Following detailed clinical assessment of individual client and carer risk

Features and Options



Disassembled by one person without tools



Curled cables for easy handling

The use of curly leads has been strictly limited to functions requiring stretch capability.



Easy Maintenance Motors

Motors can be disengaged without the need for tools, and cabling between the control unit and the motor is easy to change.

Features and Options



Well-Fitted Cable Network

Cables are designed to be hung the bed frame - as a result, there are no loose cables.



Hand control

The colours makes it easy to benefit from all the positioning options.



Scandinavian Style

Britt III side rail offers a modern style to compliment the beech finish bed ends and helps complete the natural look.



Safe support in bed

Lifting pole with safe adjustable grip handle. Can be attached to either right or left side of bed.



Support rails

Supporting rails provide added safety and support when getting out of bed. No additional components or fittings are required.



Transport unit

For easy storage and enhanced portability.




Individuality & Customisation

Mattress support extension facilitates optimal length adjustment. Choose between 100 mm or 200 mm extension.



ACP-Box

By using the ACP-box, all handset functions can be blocked.

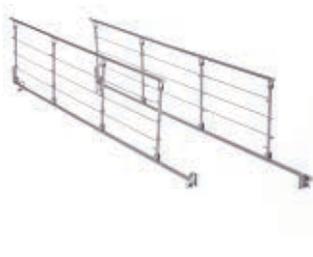

Invacare Scala Basic 2

Foldable steel side rail. Compliant with the new standard for beds IEC 60601-2-52.



Invacare Scala Basic Plus 2

Foldable steel side rail with plastic inserts at the mounting brackets. Compliant with the new standard for beds IEC 60601-2-52.



Invacare Scala Medium 2

Foldable steel side rail for higher mattresses. Compliant with the new standard for beds IEC 60601-2-52.



Invacare Scala Decubi 2

Foldable steel side rail for extra high mattresses. Compliant with the new standard for beds IEC 60601-2-52.

Please see *Invacare Bed Accessories* brochure for more accessories.

Technical data

Etude Classic	Width outside: 1030 mm	Length outside: 2400 mm	Height adjustable: 200 - 600 mm*	Mattress support dimensions: 4-sectioned: 850 - 240 - 350 - 560 mm	0 - 70° Backrest angle	0 - 20° Thigh angle
Extra Low	Width inside: 900 mm	Length inside: 2000 mm (+200 mm extension)				

Etude Classic	0 - 25° Legrest lift	11° Tilt option**	92 kg	SWL 220 kg (34.5 stone)	Bariatric use
Extra Low					

- * Use of mobile hoists must always be carried out when bed is in high position
- ** Please note that the tilt function is only practised for clients whilst under medical supervision.

NB/ The bed must not be used by patients under 12 years of age, or by patients with body size equivalent to an average 12 year old or smaller.

For mattress compatibility, please refer to the user manual which can be found at www.invacare.co.uk

Frame colour and wood work



Products, Options and Features can be subject to change.

Invacare® | **SB 755**



3 year warranty

When comfort comes first

Combining comfort, ease of handling and safety, the *Invacare SB 755* is ideal for all homecare and institutional environments.

The *SB 755* offers optimum comfort with an ideal seating positioning to minimise shear and friction. Carers can adjust the bed to meet individual user requirements.

Thanks to a unique technical feature the *SB 755* can be adapted to taller users as both the head and foot of the bed can be extended.

All these features provide real benefits to users and care-givers in terms of comfort and general well-being.



Yes, you can.®



Your comfort is our focus

The **SB 755** comes with several features to ensure an excellent seating position and to reduce shear and friction during movement. These include an improved auto-regression system, optimised length of mattress platform sections and a leg section that can be lowered below horizontal.

Invacare strives to cater for all patients and so offers a length extension of 100 mm at the head and foot of this model to meet the needs of taller users.

Ease of handling

No part of the **SB 755** weighs more than 25 kg and dismantlable scissors are available as an option.

Mounting and dismantling the bed has been made particularly easy to allow for greater ease of transport and storage. The bed can be stored on its base and rolled directly from van to bedroom. Assembly and transportation can be carried out by one person quickly and with minimum effort.

Safety is key


As always at Invacare, safety is our key concern. Side rails in both the standard bed length and the extendable version meet the increased strength demands of

the new IEC 60601-2-52_FDIS standard for beds, without compromising the ergonomics or working conditions of the carers.

Features and Options



Auto regression – improved functionality. Increased comfort when activating backrest, minimising shear/friction forces.

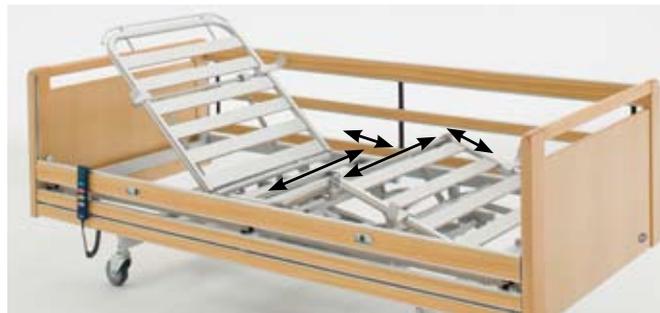

Optimised dimensions of mattress platform according to anthropometrical measurement, for increased comfort.



Extendable mattress platform. Additional 100 mm at head and foot. Bed easily extends with 50, 100, 150 or 200 mm to suit user needs.

Features and Options



Line Side rail

Full-length side rail meets the tough demands of strength and safety in the expected forthcoming standard for beds (IEC 60601-2-52_FDIS). The Line side rail is lightweight and easy to use.



Line Extendable side rail

Full-length side rail, which can be adjusted incrementally from 2000-2200 mm.



Removable panel

On the Victoria bed end, a removable panel allows easy access from the head and foot of the bed even with the full-length side rails mounted.



SB 755 Detachable



Bed ends can be dismantled without tools.

**Need to move the bed from one place to another?
No lift?
Narrow staircase?
...No problem!**

The **SB 755 Detachable** is designed to overcome all these challenges. This bed is so easy to mount and dismantle that it can be carried out by one person; our unique dismantlable system is designed for minimum effort - no need to worry about losing screws or parts! Designed with parts that weigh no more than 25kg each, the **SB 755 Detachable** is a novel design that is unique in the market. Such easy handling ensures that the bed is easy to carry – to whatever floor.



The seat section opens and the two mattress support parts are divided. The mattress platform is detached from the base with a quick release mechanism.



Mattress platform is effectively released and unlocked in the seat section.



The scissor is released from the base with the quick release handle.



Activate the quick release and the scissor is detached. The handle makes it very easy to dismantle.



Bed parts, mattress and accessories are stored directly on the base with special hooks; storage and transport are facilitated and the bed takes up less space.

Features and Options



Victoria bed end

Detachable panel - Available in beech, cherry and white beech. Compatible with the Line side rail range.



Sophie bed end

Available in beech, cherry and white beech. Compatible with the Line side rail range.



Emma bed end

Detachable - Available in beech, cherry and white beech. Compatible with the steel side rail range.



Susanne bed end

Available in beech, cherry and white beech. Compatible with the Line side rail range.



Anita bed end

Available in beech, cherry and white beech. Compatible with the Line side rail range.



Several handcontrols:

Invacare **Soft Control** or standard with "green" button function.



Lifting pole

The handle can be adjusted both in height and depth.



Swivel support handle

The swivel support handle provides good support for the client getting to a standing position.



Invacare **Scala Basic 2**

Foldable steel side rail.*



Invacare **Scala Basic Plus 2**

Foldable steel side rail with plastic inserts at the mounting brackets.*



Invacare **Scala Medium 2**

Foldable steel side rail for higher mattresses.*



Invacare **Scala Decubi 2**

Foldable steel side rail for extra high mattresses.*

*Compliant with the new standard for beds IEC 60601-2-52.



Cover for wooden side rail

Padded cover for wooden side rails. Cover remains in place when side rail is lowered.



Mattress support

Mattress support with aluminium slats.



Central brake

Available with a wheel diameter of 100, 125 or 150 mm.



Lowering kit

Using castors of 75 mm creates a lower height of 250 mm.

Please see **Invacare Bed Accessories** brochure for more accessories.

Technical data

SB 755	Width inside: 850 / 900 mm Width outside: 950 / 1010 mm	Length outside: 2130 mm Length inside: 2080 mm Mattress: 2000 mm	Height adjustable: 350-820 mm Total height: 1600-2070 mm	Mattress support dimensions: 800 - 260 - 260 - 680 mm	Under bed clearance: 170 mm	0° - 73° Backrest angle

SB 755	0° - 28° Thigh angle	0° - 16° Minimum angle below horizontal	0° - 18° Full legrest lift	25° (Angle below frame) 14° (Angle above frame)	SWL 235 kg (37 stone) Max. user weight 200 kg (31.5 stone)	Heaviest part of the product: 25 kg	Bariatric use

NB/ The bed must not be used by patients under 12 years of age, or by patients with body size equivalent to an average 12 year old or smaller.

For mattress compatibility, please refer to the user manual which can be found at www.invacare.co.uk

Frame colour and wood work



Aluminium grey

Beech

Invacare® | *Octave*®



3 year warranty

Setting New Standards in Bariatric Care

With a safe working load up to 60 Stone (385 kgs), the *Invacare Octave* is a robust yet stylish profiling bed solution that assists in the care and comfort of bariatric clients. With a fully electrical, 4-sectioned profiling mattress support as standard, the *Octave* provides optimal comfort and support to the client, as well as an ideal working environment for the carer. The frameless construction of the mattress platform allows improved side access to the client, and the bed ends can be detached independently of the side rails for quick and easy access to both head and feet.



Yes, you can.®



Features and Options

Great comfort and stylish design

The *Invacare Octave* is designed specifically with the bariatric client in mind. The bed frame extends to the edges of the bed base for enhanced support and ease of transfer. An easy to use hand control allows the client to operate the bed's height and profiling function. The *Octave's* heavy duty motors enable a smooth and almost silent operation. The *Octave* is a substantial bed, however, its design includes all the well-known qualities of Invacare's Scandinavian design.

Delivery and storage

To make delivery as easy as possible, the *Octave* is delivered on a transportation rack. Dismounting and assembly of the *Octave* couldn't be simpler, despite the size of the bed. Furthermore, the transportation rack makes it possible to economise on storage space and makes the *Octave* very easy to move between locations.



Support handle
For increased support whilst bed leaving.



Lifting pole
The handle can be adjusted in both height and depth, to offer support whilst re-positioning.



Extendable bedend
The bedend is extendable to offer mattress platform extension of 100 or 200 mm.

Technical data

<i>Invacare Octave</i>					
	Width outside: 1190 mm	Length outside: 2210 mm Length inside: 2070 mm Mattress: 2000 mm	Height adjustable: 400 - 800 mm Total height: 1720 - 2120 mm	Mattress support dimension: 810 - 340 - 330 - 520 mm (4-sectioned)	0° to 57° Backrest angle
	0° to 30° Thigh angle	0 to 30° Legrest lift	Steel 180 kg Wood: 199 kg	SWL 385 kg / 60 stone Max user weight: 350 kg / 56 stone	

Frame colour



Silver

Beech

NB:

The bed must not be used by patients under 12 years of age, or by patients with body size equivalent to an average 12 year old or smaller.

Invacare® | *Bed Accessories*



Supporting your needs

The *Invacare* bed range features a wide range of accessories which gives you the possibility to customise your bed both in terms of design and taste. This allows perfect fittings to personal environment with supporting accessories, tailored to individual client needs.

In the catalogue you will also find compatibility references for beds, side rails and bed ends.



Bed ends



Victoria bed end

Detachable bed end for **SB 755** in beech, white beech and cherry. Compatible with Line full length side rails.



Susanne bed end

Available in beech, white beech and cherry. Standard/low version for **SB 755**, compatible with Line full length side rails.



Anita bed end

Available in beech and cherry. Standard/low version for **SB 755**, compatible with Line full length side rails.



Sanne bed end

Available in beech. High/low version for **SB 750** and compatible with full length side rails.



Piggy bed end

Detachable bed end in beech. Available for **SB 750** and compatible with steel side rails.



Oda bed end

Detachable bed end in beech. Available for **SB 750** and compatible with steel side rails.



Flex bed end

Available in beech. High/low version for **SB 750** and compatible with full length side rails.

Side rails



Line

Full length aluminium side rail that meets the forthcoming bed safety standard (IEC 60601-2-52_FDIS). Available in beech, white beech and cherry.




Line Extendable

Full length aluminium side rail that meets the forthcoming bed safety standard (IEC 60601-2-52_FDIS). Available in beech, white beech and cherry. Adjustable length 2000-2200 mm.



Britt IV

Full length wooden side rail. Available in beech, white beech and cherry.

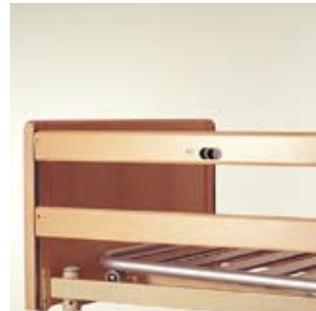

Britt III

Full length wooden side rail. Available in beech, white beech and cherry. Available in short or extended version.

Side rails



Lilly III
Full length wooden side rail with three bars. Available in beech and cherry, also in an extended version.



Eclipse
Metal side rail for **Etude**.



Invacare Scala Basic 2
Foldable steel side rail. Compliant with the forthcoming standard for beds IEC 60601-2-52.



Invacare Scala Basic Plus 2
Foldable steel side rail with plastic inserts at the mounting brackets. Compliant with the forthcoming standard for beds IEC 60601-2-52.



Invacare Scala Medium 2
Foldable steel side rail for higher mattresses. Compliant with the forthcoming standard for beds IEC 60601-2-52.



Invacare Scala Decubi 2
Foldable steel side rail for extra high mattresses. Compliant with the forthcoming standard for beds IEC 60601-2-52.



Wood side rail for Medley
Full length wooden side rail for **Invacare Etude Medley** with two bars. Available in beech only.



Steel side rail for Medley
Full length steel side rail for **Invacare Etude Medley**.



Side rail height extender (+150 mm)
for Line, Britt III, Britt IV and **Invacare Etude Medley** full length wooden side rails.

Side rail covers



Cover for Line side rail

Padded covers for Line side rail. Cover remains in place when side rail is lowered. Available in padded, extra padded and extended + 100 mm and + 200 mm.



Net cover for Line side rail

Net covers for Line side rail. Cover remains in place when side rail is lowered. Available as standard and extended + 100 mm and + 200 mm.



Cover for wooden side rail

Padded cover for wooden side rails. Cover remains in place when side rail is lowered. Available in standard, extended, short and net.



Cover for metal side rail

Padded cover for metal side rails. Cover remains in place when side rail is lowered. Available in standard, extended, and extra padded.



Safety cover*

For wooden side rail. Recommended for restless users.

Support handles



Support handle

A very stable support handle for getting in and out of bed:

- H250 x W300 mm
- H250 x W800 mm
- H400 x W300 mm
- H400 x W500 mm
- H400 x W950 mm

* When using side rails and the according safety cover, it is essential to ensure correct fitting otherwise there is a danger of entrapment or suffocation between the mattress support, side rail and bed end.

Support handles



Swivel support handle

The swivel support handle makes it possible to swing the legs out over the handle before the grip is folded out. Can be used alone or in pairs on the **SB 755** and **SB750**.



Support handle for *Octave*®

Makes getting in and out of bed much easier for users. H400 x W300 mm.



Sit-up support

Help for getting in and out of bed. For **SB 755** and **SB 750**.

Lifting poles



Lifting pole

Catching handle can be adjusted in both height and depth.



Swivel lifting pole

The lifting pole swings away to the side of the bed. Catching handle can be adjusted in both height and depth.



Lifting pole

Catching handle can be adjusted in both height and depth. For use with **Octave**® beds.



Free-standing lifting pole

Height adjustable, catching handle can also be adjusted.

Other accessories



Distance bumper

Prevents walls from being damaged when is the bed is moved. Can be placed on the base frame at the end of the bed or at the side of the bed.



Bumper wheels

The bumper wheel can be mounted horizontally or vertically. Available in three different sizes.



Castors


Different sizes of castors. Comes in 75, 80, 100, 125, and 150 mm.



Lowering kit

By using 75 mm wheels **SB 755** and **SB 750** can be lowered to 250 mm.

Other accessories

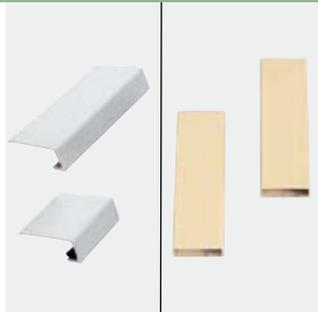

Cover tubes

To be mounted on the bed frame to cover the tubes when the bed is extended. For **SB 755** and **SB 750** in different sizes and colors.



Limitation crossbar

Prevents the leg support from dipping below horizontal position. Can be placed on the base frame. Available in ivory and aluminium. For **SB 755** and **SB 750**.



Transportation fitting for Etude®

Fittings for easy transport of the **Etude** beds.



Extension parts for Etude®

Available in 100 mm or 200 mm extension.



Extension for Octave®

Extension for **Octave** for greater flexibility. 100 mm to 200 mm extension.



Wooden adapters – Octave bed

Wooden adapters for mounting **Octave** side rail, comes in beech and cherry.



Side rail lock

Strap for locking the side rail in its lowest position. Can be used with any side rail.



Wooden front panel

To be mounted on the side of the bed. In combination with Flex, Sanne, Victoria, Sophie, Susanne, Anita bed ends.




Wooden front panel

To be mounted on the side of the bed. In combination with Piggy and Oda bed ends.




Transportation rack

Can be used **SB 755** and **SB 750**. For transporting or storing an assembled bed.

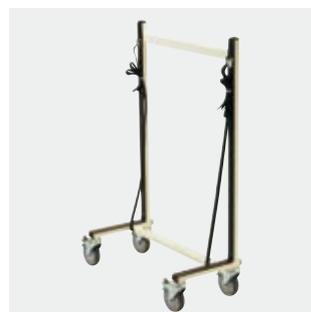

Transportation rack

To be used for transporting a bed.

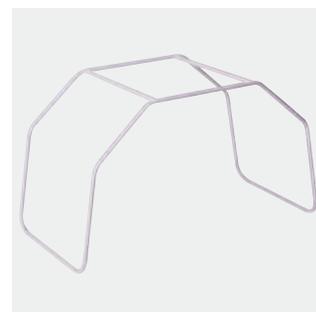

Bed Cradle

Bed cradle can be placed on the bed end to alleviate pressure from the duvet on the legs/feet.

Other accessories



Rastofix – handle
For easy adjustment of the leg section.



IV Rack
① Free-standing IV rack. Height-adjustable, castors for easy movement.
② Can be placed in the lifting pole tube or on the bed frame.



IV Rack mount fittings
Fittings for mounting the IV rack either in the lifting pole tube or on the bed frame.

Electrical accessories



Hand control
Hand control HB80 and HL80 with green button function. Each function of the HL80 can be individually locked. Available for **SB 755**.



Hand control
Hand control HB70. With green button function for user comfort. Available for **SB 755** and **SB 750**.



Soft control
Function-based hand control with easy activated buttons. Available for beds with 4 sections.



Hand control for Octave®
Hand control HL85. Each function on the HL85 can be individually locked.



Hand control for Alegio® standard
Dewert hand control for **Alegio** with green button function.



Hand control for Alegio with locking function
Dewert hand control for **Alegio**, can be locked individually.



Holder for hand control
Flexible holder for hand control, easier access for users. Available for soft control, HB70, HB80, HL80 and HL85 and Dewert hand control. Comes in aluminium and ivory.



T-cable
Can be used when two hand controls are needed.

Electrical accessories



Battery back-up

To be used in case of power failure. The bed can be delivered prepared for battery back up. The battery back up will be installed afterwards. Available for **SB 755** and **SB 750**.



Extension cable

2.5 metre extension cable for Linak hand control. Can be used for **SB 755**, **SB 750**, **Etude®** and **Octave®**.



ACP box

Gives healthcare personnel the ability to lock the function of one or more motors. Can be used for **SB 755**, **SB 750**, **Etude** and **Octave**.



ACP box

Gives the healthcare personnel the ability to lock the function of one or more motors. Can be used with **Alegio**.



Field power supply

Can be used with **SB 755** to adjust the height or any other motor when no power is available. For service personel only, not to be used with a user in the bed.

Furniture



Bed tray

Stable table/tray, fits directly on all beds with full length wooden side rails.



Pausa

Height-adjustable tilt table with plate that can be adjusted from 720 - 1160 mm in height. Max load: 10 kg.



Kauma

Height-adjustable table with plate that can be adjusted from 710 - 1135 mm in height. Max load: 10 kg.



Sonnet Combi

Bedside table in beech. All horizontal surfaces are in plastic laminate. Practical top shelf, one drawer and a cupboard. The flap can be adjusted from 310 - 930 mm in height.

Furniture



Back Support

Manual backrest support, with possible head support. Can be adjusted from a 45 - 60 degree angle.

Pressure Reducing Accessories



Softform® Flexipad

Highly versatile pads filled with silicone gel, used to position and support limbs, joints, the torso or the head.



Softform Heelpad

Enables vulnerable heels to be immersed into independent gel sacs, offering significant pressure reduction and comfort to clients.



Softform Odstock Wedge

Assists in the reduction of interface pressures on the sacrum and heels and ideal for improving posture for those with back discomfort.



Propad® Leg Trough

Elevates the heels and allows them to be suspended over the edge of the trough, providing pressure relief for the vulnerable heel area.



Flo-tech® Back Support

Provides a high degree of stability and improved posture, resulting in enhanced protection from pressure ulcer development



Flo-tech Lumbar Pad

Designed for use on existing chairs or wheelchairs and can improve posture, comfort and support.

Invacare® | **Mattresses**



INVACARE
Yes, you can.®



 Invacare® | **Mattresses**
Mattresses

<i>Invacare Softform® Premier</i>	130-131
<i>Invacare Softform® Excel</i>	132-133
<i>Invacare Softform® Premier MaxiGlide</i>	134-137
<i>Invacare Softform® Premier Active</i>	138-141
<i>Invacare Softform® Premier Visco</i>	142-143
<i>Invacare Softform® Odstock Wedge</i>	144-145
<i>Invacare Softform® Heelpad</i>	146-147
<i>Invacare Softform® Flexipad</i>	148-149
<i>Invacare Propad® Overlay</i>	150-151
<i>Invacare Propad® Premier Overlay</i>	152-153
<i>Invacare Propad® Visco Overlay</i>	154-155
<i>Invacare Propad® Leg Trough</i>	156-157

Invacare® | **Softform® Premier**



Warranty
Foam 10 Years
Cover 5 Years

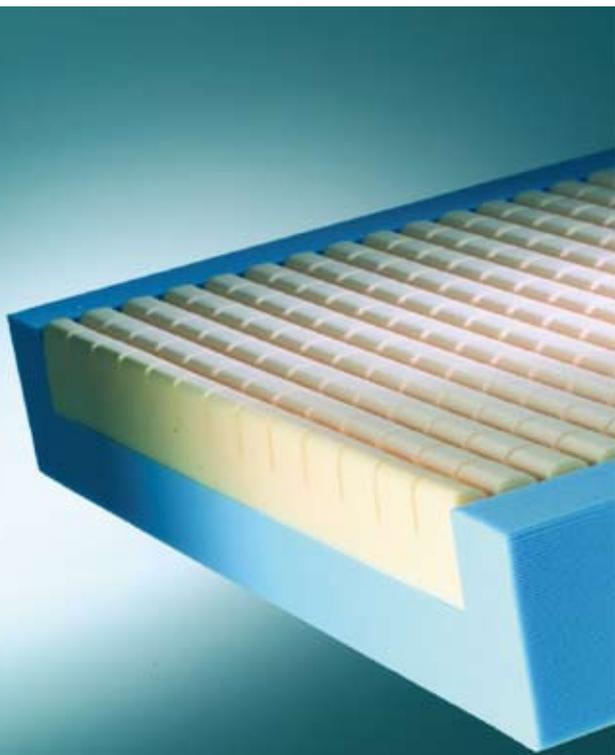

Clinically Proven - Pressure reducing mattress

The *Invacare Softform Premier* is a high specification replacement mattress designed to meet the demands of the Community and Acute healthcare environment. Practical and durable, this mattress delivers exceptional levels of comfort and pressure reduction. The **Softform Premier** is suitable for those considered at 'Very High Risk'* of developing pressure ulcers.



Yes, you can.®



Effective pressure ulcer prevention

The *Invacare Softform Premier* features a unique contoured insert pad, constructed from 'keyhole cut' castellated foam cells, which move independently to provide excellent patient comfort, reduce shear and friction and optimise pressure reduction. The mattress insert has a non glued construction, which allows excellent conformity when profiled on electric bed frames.

Cover is resistant to fluid ingress


The outer cover is a multi stretch, waterproof, polyurethane fabric with high frequency welded seams; and a full length flap concealing the zip. These features combine to effectively reduce the potential of fluid entering into the mattress. The base of the mattress features a toughened PU coating to prolong the longevity of the mattress.

Features



Does not require turning or rotating

This mattress does not require turning or rotating, significantly minimising manual handling issues.

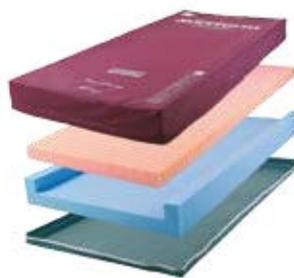

Single piece foam "U" core

Facilitates patient transfer by providing effective side-wall support. Robust foam base prevents the mattress from bottoming out.



Full length flap conceals zip

The full length flap protects the zip. This reduces the possible risk of fluid ingress, preventing contamination of the foam core.

Technical data



Softform Premier Single Size (660 - 900mm)
Double Size (910 - 1370 mm)
King Size (1380 - 1525 mm)



Length: 1410 - 2100 mm



Height: 152 mm



Weight: 13 kg



Max. User Weight: 247.6 kg (39 Stone)



Softform Premier Cleaning: 80°

Manufactured to European Standards: BS EN 597-1 Cig, BS EN 597-2 Match.
Foam Interior: BS 5852: 1990 Crib Source 5; Cover BS7175: 1989 Crib Source 5.

Colour



Fuchsia (cover)

* Essential nursing care is pivotal in pressure ulcer prevention. This mattress will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action based planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. A range of assessment scales can be used as a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.



Invacare® | **Softform® Excel**



Warranty
Foam 2 Years
Cover 2 Years



Pressure reducing mattress

The *Invacare Softform Excel* provides a high level of comfort. The pressure reducing mattress core is constructed from high specification foam which provides even weight distribution for greater comfort and pressure re-distribution. With a robust foam base it prevents the mattress from bottoming out and is suitable for those considered 'At Risk'* of developing pressure ulcers.





Features and Options

Effective pressure ulcer prevention

The core of the *Invacare Softform Excel* is constructed from a series of 'keyhole cut' cells, which move independently to evenly distribute weight, promoting greater comfort and facilitating pressure ulcer prevention. The robust base features a firm layer of high density foam to provide substantial support.

Multi stretch PU waterproof cover

The cover fabric provides excellent recovery and durability, and helps to reduce shear and friction forces. The base of the mattress is a tough coated polyurethane (PU) which protects the mattress when used with electric profiling bed frames.




Zip around three sides

Allows easy inspection of the foam core.



Rotate-only mattress

Rotate-only turning cycle will increase product longevity. Non turn reduces the amount of manual handling required.

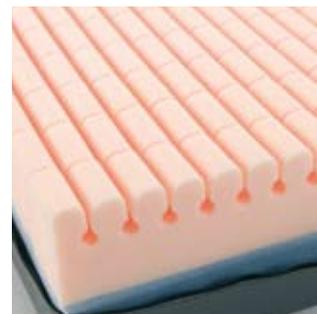

Strong, durable construction

The *Invacare Softform Excel* is an established choice of healthcare professionals.

Technical data

Softform Excel	Standard width: 800, 855, 880, 900 mm	Standard length: 1950 mm	Height: 152 mm	Weight: 10 kg	Max. Client Weight: 140 kg (22 stone)	Cleaning: 80°

Manufactured to European Standards:
BS7177: Crib 5

Colour



Blue (PU)

* Essential nursing care is pivotal in pressure ulcer prevention. This mattress will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action based planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. A range of assessment scales can be used as a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.



Invacare® | *Softform® Premier MaxiGlide*



Warranty
Foam 8 Years
Cover 4 Years

Designed to reduce shear and friction

The Invacare **Softform Premier MaxiGlide** mattress offers a new and innovative glide mechanism that significantly reduces shear and friction forces when used with an electric profiling bed. The patented design of the mattress allows the top surface to provide excellent pressure reducing properties along with a high level of comfort and support, whilst the base of the mattress uniquely separates to conform to the bedframe. The **Softform Premier MaxiGlide** is suitable for use with those who are considered to be at 'Very High Risk'* of developing pressure ulcers.



Yes, you can.®



Effective pressure ulcer prevention

The **Softform Premier MaxiGlide** features a deep-contoured insert pad constructed from 'keyhole' cut castellated foam cells, which move independently to maximise the body contact area and optimise pressure reduction.

Reduces the need for manual handling

The profiling mattress maximises the benefits of an electric profiling bed, delivering improved positioning whilst reducing the requirement for repositioning of the client.


Multi stretch cover

The mattress cover is multi stretch and conforms to the deep contoured castellations on the upper surface of the mattress, minimising rucking and further reducing the potential for tissue damage.

Features



Smooth sliding surface

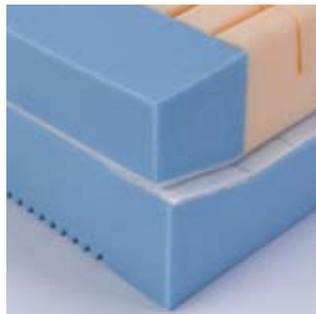

Top surface remains with patient




Bottom surface profiles with the bedframe

Features



Patented design

The "Glide" mattress system is only available from Invacare and approved distributors.

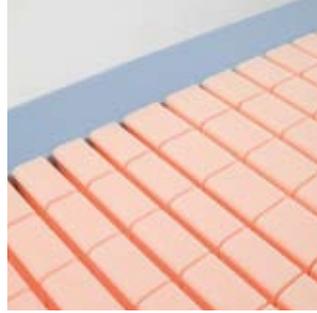

Modular Construction

Enables easy replacement of all parts, increasing product longevity.



Pressure Ulcer Prevention

The mattress provides an excellent pressure reducing surface, designed to maximise the body contact area and prevent the build up of pressure ulcers*.



Cover is resistant to fluid ingress

Features include two-way stretch, vapour-permeable and waterproof polyurethane cover with high frequency welded seams and a full length flap concealing the zip.



Fully detachable top cover and toughened polyurethane base

The toughened PU coated base protects the mattress against metal bed frames. The top cover is easily removed for laundering and foam inspection.



Does not require turning or rotating

Mattress does not require turning or rotating, significantly reducing manual handling risk.

Technical data

**Soffform
Premier
MaxiGlide**



Std Width:
800, 855, 880,
900 mm



Std Length:
1950-2100 mm



Height:
152 mm



Weight:
13.5 kg



Max. Client Weight:
247.6kg / 39 stone



Cleaning: 80° C
Refer to
instructions
for use and relevant
infection control
guidelines

Manufactured to European Standards:
Interior Foam: BS 5852 Crib 5
Cover: BS 7175 Crib 5

Colour



Baltic Blue

* Essential nursing care is pivotal in pressure ulcer prevention. This mattress will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action based planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. A range of assessment scales can be used as a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.

Invacare® | *Softform® Premier Active*




Warranty
Foam 10 Years
Cover 5 Years
Pump 2 years
Insert 2 years

Static Dynamic pressure reducing mattress

The unique patented design of the *Invacare Softform Premier Active* features an alternating air insert beneath the castellated foam of the tried and trusted *Softform Premier* mattress. The mattress retains its properties as a 'High/Very High Risk'* static surface, but should a patient require stepping up to a dynamic surface, a discreet pump can be fitted to the air insert, transforming the mattress into an alternating surface delivering additional levels of pressure relief. There is no requirement to transfer the client, thereby minimising any discomfort.



Yes, you can.®

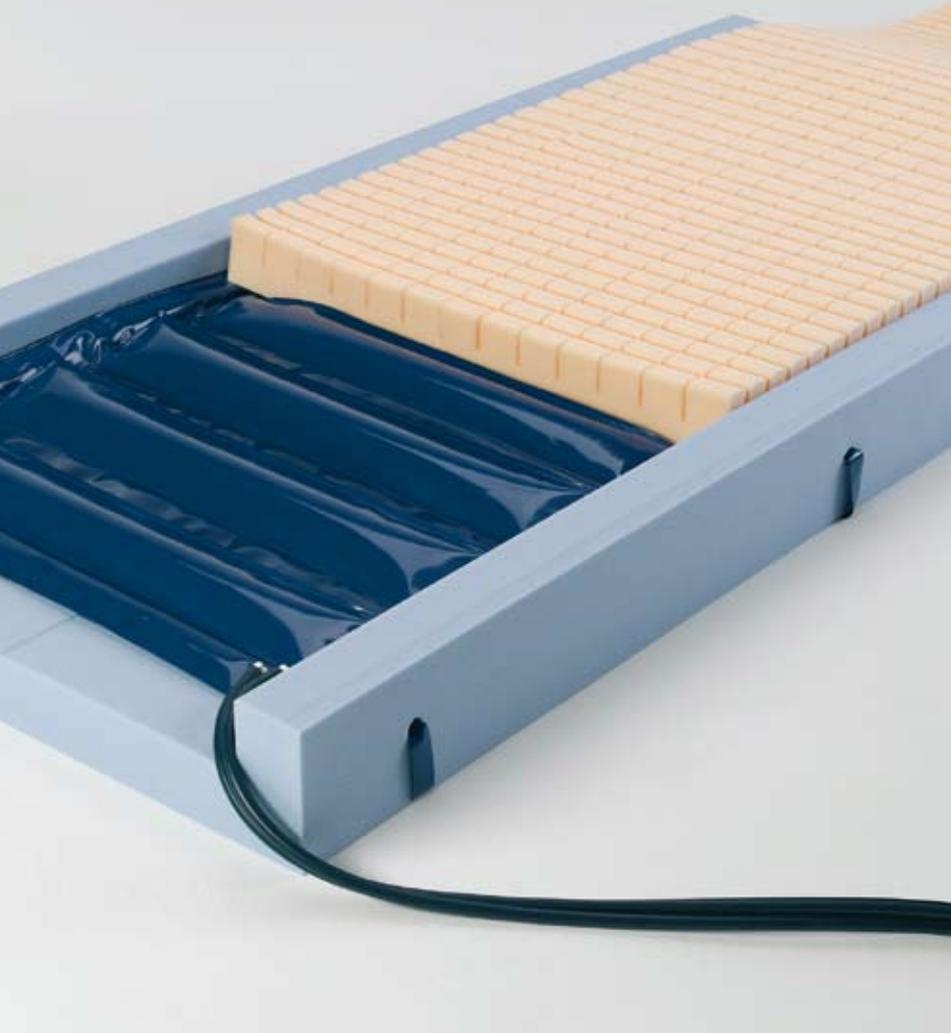

Static Dynamic mattress

The air insert contained within the mattress foam construction allows stepping up to and down between static and dynamic alternating systems with minimum patient disruption. The air cells alternate on a ten minute cycle beneath the foam insert, delivering a very low impact at the point of change, thereby delivering high levels of comfort in addition to improved pressure reduction.

Pump automatically assesses patient weight

Sensory software within the pump can gauge the patients weight, enabling pressure to be maintained accordingly at an appropriate level. The pump is also fully alarmed in case of power failure or fault.

Multi stretch cover

The mattress cover is multi stretch and conforms to the deep contoured castellations on the upper surface of the mattress, minimising "rucking" and further reducing the potential for tissue damage.

Features



Air insert within Softform Premier mattress foam construction

Allows stepping up and down between static and dynamic alternating systems with minimum patient disruption.



Inflation of air cells alternates on a 10 minute cycle

Delivers an extremely comfortable mattress surface, while reducing interface pressures on any one area of the body.



High risk static mattress

In the event of a power failure, the patient is still resting on a 'High Risk' pressure reducing surface, minimising the potential for pressure damage occurring.

Features



The *Invacare Softform Premier* is a leading pressure-reducing static mattress.



To step up a patient, retract the pipes from the mattress.



Fit pipes to the Active pump.



Simply switch on, no additional adjustments required.



The client is now on an alternating surface.

Technical data



Ranging from
800 - 900 mm



Ranging from
1810 x 2100 mm



152 mm



13 kg



247.6 kg
(39 stone)



Cleaning: 80°

**Softform
Premier Active**

Manufactured to European Standards: BS EN 597-1 Cig, BS EN 597-2 Match, BS7177:1996 Crib Source 5. Tested to pass ISO 1419:1995.

Colour




Fuchsia (cover)

* Essential nursing care is pivotal in pressure ulcer prevention. This mattress will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action based planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. A range of assessment scales can be used as a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.

Invacare® | **Softform® Premier Visco**



Warranty
Foam 4 Years
Cover 4 Years

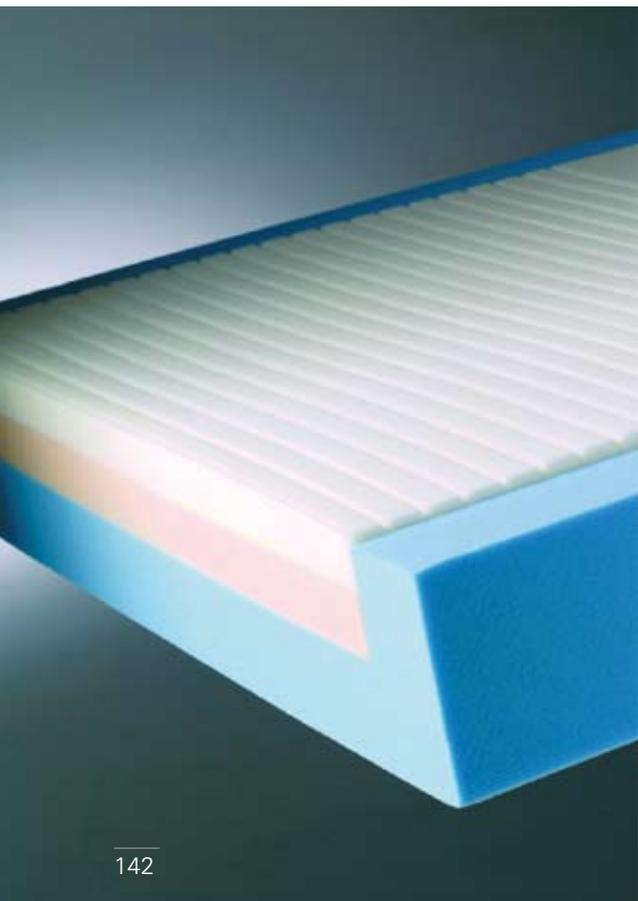

Pressure reducing visco elastic foam mattress

The *Invacare Softform Premier Visco* is a pressure reducing mattress suitable for those considered at 'High Risk'* of developing pressure ulcers. The mattress features a high density visco elastic (memory) foam insert, which is temperature-sensitive to the body, allowing it to conform to body shape. This provides exceptional comfort and facilitates pressure reduction by minimising peak pressures.





Features and Options

Comfort and pressure reduction combined

The *Invacare Softform Premier Visco* is constructed from a range of superior quality materials. The lower 'U' shaped layer is made up of robust high density foam, to prevent the mattress bottoming out, and to improve longevity. The central core includes a high density, visco elastic (memory) foam top section which warms to conform to body shape, maximising the body contact area to give good pressure reduction.

Exceptional comfort

By conforming to the individual body shape of the client, the high quality visco elastic foam featured in this mattress offers exceptional comfort and support.



Cover is resistant to fluid ingress

The multi stretch, waterproof polyurethane cover has high frequency welded seams and has a full length flap concealing the zip. The cover also has a toughened PU coated base to prolong the longevity of the mattress.



Modular Construction

Enables easy replacement of parts, increasing product longevity. Insert pad can be interchanged with regular castellated foam, if required.



Does not require turning

Mattress does not require turning, significantly reducing amount of manual handling issues.

Technical data



Softform Premier Visco

Single Size (660 - 900 mm)
Double Size (910 - 1370 mm)
King Size (1380 - 1525 mm)



Length: 1410 - 2080 mm



Height: 152 mm



Weight: 13 kg



Max. User Weight: 247.6 kg (39 stone)



Softform Premier Visco

Cleaning: 80°

Manufactured to European Standards: BS EN 597-1 Cig, BS EN 597-2 Match, BS7177: 1996 Crib Source 5

Colour



Fuchsia

* Essential nursing care is pivotal in pressure ulcer prevention. This mattress will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action based planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. A range of assessment scales can be used as a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.



Invacare® | *Softform® Odstock Wedge*



Pressure reducing wedge

The *Invacare Softform Odstock Wedge* was developed in partnership with clinicians to assist in the reduction of interface pressures on the sacrum and heels without significantly increasing pressure on the thighs or calves. It is designed to improve and maintain patient posture in the recumbent position, and to reduce the tendency for the client to slip down the bed, thereby reducing shearing forces on the sacrum. The **Odstock Wedge** is ideal for improving posture for those with back discomfort.



Yes, you can.®



Features and Options

Reduces pressure on the heels

Clinical research has effectively demonstrated that the degree of knee flexion necessary to achieve substantial levels of pressure relief is relatively small. The *Invacare Softform Odstock Wedge* was designed by clinicians to reduce pressure on the heels through re-positioning the hip and elevating the knee.

Reduces lumbar lordosis

By tilting the pelvis, the **Odstock Wedge** also reduces lumbar lordosis – desirable following back surgery and ideal for improving the sleeping posture of patients with back problems.



High frequency welded seams

Ensures the core is protected from fluid ingress, increasing product integrity.



High quality construction

The multi stretch cover, is waterproof and vapour-permeable. The foam core is constructed from, high density foam.



Designed to be used in conjunction with the *Invacare Heelpad*

Used in conjunction, these two products combine to significantly reduce pressure on the vulnerable heel area.

Technical data



Width: 570 mm



Length: 865 mm



Height: 115 mm



Weight: 1.5 kg



Cleaning: 80°

Softform Odstock Wedge

Refer to instructions for use and relevant infection control guidelines.

Manufactured to European Standards:
Interior Foam: BS 5852: Crib 5,
Cover: BS7175: Crib 5

Colours



Fuchsia (Standard)



Blue (Option)

Warranty: 1 year



Invacare® | **Softform® Heelpad**



Pressure Reducing Heelpad

The *Invacare Softform Heelpad* is designed to reduce pressure on the vulnerable heel area. The controlled volume of gel in the sacs allows pressure displacement and even weight distribution. Pressure on a patient's heels are kept to a minimum when the heels are in contact with the gel sacs. The *Invacare Softform Heelpad* was designed to be used in conjunction with the *Invacare Odstock Wedge*, but can also be used independently as an effective pressure reducing device.





Features and Options

Designed to reduce pressure on the heels

The *Invacare Softform Heelpad* enables vulnerable heels to be immersed into independent gel sacs, offering significant pressure reduction and comfort to clients.

Enhanced combined performance

The *Softform Heelpad* offers enhanced performance, when used in conjunction with the *Softform Odstock Wedge*, by providing improved support for the lower limbs, and pressure reduction for the heels.



Single and Double Heelpads available

Single and Double formats available to meet specific patient requirements.



Heelpads supplied with three platilon covers

Platilon protective covers are disposable, and help to avoid contamination of the heelpad.



Designed to be used in conjunction with the Invacare Odstock Wedge

Used in conjunction, these two products combine to significantly reduce pressure on the vulnerable heel area.

Technical data

Softform Heelpad					
	Width: 620 mm	Length: 480 mm	Height: 50 mm, 32 mm	Weight: 3.58 kg	Cleaning: 80°

Manufactured to European Standards:
 Interior Foam: BS 5852: 1982 Crib 5, BSEN 597-1,2
 Cover: BS 7175: 1989 Crib 5

Colour



Fuchsia

Warranty 1 Year



Invacare® | **Softform® Flexipad**



A range of versatile support pads

Invacare Flexipads are highly versatile pads filled with silicone gel, and can be used to position and support limbs, joints, the torso or the head.



Yes, you can.®

Improves positioning

Invacare **Flexipads** can be used to position limbs for treatment, whilst recovering from an injury, or to assist in maintaining a more comfortable position to alleviate pain.

Improves support

Easy to use and adjust to position, **Softform Flexipads**, can provide supplementary short-term support for patients experiencing difficulty in maintaining posture.



Features and Options



Single Flexipad

Fitted with a strap for fastening to chairs, the **Invacare Flexipad** increases the surface area under elbows, arm or back.



Double Flexipad

Two gel sacs with a foam core – highly versatile in a number of applications, including supporting large bony prominences.



Flexipad Ultra

A lighter, more flexible format, available with or without straps.

Technical data



Width: 220 mm



Length: 330 mm



Height: 50 mm, 32 mm



Weight: 0.8 kg



Cleaning: 80°

Softform Flexipad

Manufactured to European Standards:
Cover: BS 7175: 1989 (Crib 5)

Colours



Blue



Black Smooth

Warranty 1 Year



Invacare® | *Propad®* Mattress Overlay



3 Year warranty



The Original Mattress Overlay

Widely used by healthcare professionals in Hospitals, Nursing Homes and in the Community, the *Invacare Propad Mattress Overlay* has proved itself to be an essential element of a successful pressure care programme. Designed to sit on top of a conventional mattress, the *Propad Mattress Overlay* provides comfort and pressure area care to those considered at 'High Risk'* of developing pressure ulcers.



Yes, you can.®

Effective Pressure Ulcer Prevention


The foam core has independent, geometrically-cut surface cells which significantly reduce shear and friction, whilst increasing the overall support surface area, thereby reducing peak pressure.

Even Weight Distribution

The **Propad Overlay** has larger, more supportive cells at the centre and smaller cells at each end. This enables the mattress overlay to achieve even weight distribution by conforming to the shape and form of the body and cradle bony prominences.



Features and Options

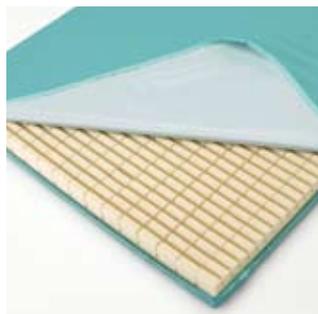

Cover helps to prevent a build up of shear and friction forces

Cover is multi stretch, water-resistant and vapour-permeable.



Available in a range of sizes

Single - 1875 x 876 x 85 mm,
Double - 1875 x 1372 x 85 mm,
King Size - 1980 x 1524 x 85 mm.



Covers have an anti-slip base

Ensures safe use of the **Propad Overlay** and helps to prevent slippage.

Technical data



Std sizes:
1875 x 876 x 85 mm ¹
1875 x 1372 x 85 mm ²
1980 x 1524 x 85 mm ³

Additional sizes available on request



Height:
85 mm



Weight: 5.5 kg



Max User Weight:
108 kg
(17 stone)



Cleaning: 80°

Propad Mattress Overlay

- 1/ Single
- 2/ Double
- 3/ King Size

Manufactured to European Standards: BS 7177 Crib 5

Warranty 3 Years

Colour



Green (PU)

* Essential nursing care is pivotal in pressure ulcer prevention. This mattress will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action based planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. A range of assessment scales can be used as a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.



Invacare® | **Propad® Premier
Mattress Overlay**



3 Year warranty



The clinically effective and easy to carry overlay

Created specifically for community application, the *Invacare Propad Premier Mattress Overlay* has been designed for use in conjunction with conventional mattresses to deliver exceptional comfort and pressure reduction. The *Propad Premier Mattress Overlay* is suitable for those considered at 'High Risk'* of developing pressure ulcers, by providing an excellent pressure reducing surface.



Yes, you can.®



Features and Options

Effective Pressure Ulcer Prevention

The *Invacare Propad Premier Mattress Overlay* has a foam core which has independent, geometrically-cut surface cells. This significantly reduces shear and friction whilst increasing the support surface area, thereby reducing peak pressure.

Even Weight Distribution

The foam core has larger, more supportive cells at the centre and smaller cells at each end. This enables the mattress overlay to achieve even weight distribution by conforming to the shape and form of the body and cradles bony prominences.



Folds down to a quarter of its length

Facilitates easy storage and delivery. Carrying and retaining strap ensures ease of handling and transportation.



Cover helps to prevent a build up of shear and friction forces against the client

Cover is multi stretch, water-resistant and vapour-permeable. Covers feature an anti-slip base.



Size: 73.8 x 34.5 x 3.3/8" (1875 x 876 x 85 mm). Suitable for a standard single bed. Additional sizes available on request.

Technical data

Propad Premier Mattress Overlay					
	Std sizes: 73.8 x 34.5 x 3.3/8" (1875 x 876 x 85 mm) Additional sizes available on request	Height: 85 mm	Weight: 5.5 kg	Max Client Weight: 108 kg (17 stone)	Cleaning: 80°

Manufactured to European Standards: BS 7177 Crib 5

Warranty 3 Years

Colour



Green (PU)
(Standard)

* Essential nursing care is pivotal in pressure ulcer prevention. This mattress will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action based planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. A range of assessment scales can be used as a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.



Invacare® | **Propad® Visco Overlay**



3 Year warranty



Dual layered overlay designed for comfort and support

The *Invacare Propad Visco Overlay* features a high density visco elastic (memory) foam core, which is temperature-sensitive, allowing it to conform to the shape of the body. This provides a support surface with exceptional levels of comfort, and pressure reduction. The *Propad Visco Overlay* is suitable for those considered at 'High Risk'* of developing pressure ulcers.



Yes, you can.®



Features and Options

Comfort and Conformity

The interior of the overlay has a top layer of high density visco elastic foam, which is temperature-sensitive. The foam warms and conforms to body shape, maximising the body contact area. This provides exceptional comfort and facilitates pressure reduction by minimising peak pressures.

High density supportive foam base

The dual layer construction, combines a high level of comfort with improved longevity, due to the supportive, robust foam base.



Cover helps to prevent a build up of shear and friction forces against the client

Cover is multi stretch, water-resistant and vapour-permeable.



Available in a range of sizes

The *Invacare Propad Visco Overlay* is available in 2 standard sizes:

Single: 1875 x 876 x 85 mm,
Double: 1875 x 1372 x 85 mm.
Additional sizes available on request.



Covers have an anti-slip base

Ensures safe use of the *Propad Visco Overlay* and helps prevent slippage.

Technical data



Size:
1875 x 876 x 85 mm ¹
1875 x 1372 x 85 mm ²

Additional sizes available on request



Height:
85 mm



Weight: 7.1 kg



Max Client Weight:
108 kg
(17 stone)



Cleaning: 80°

**Propad
Visco
Overlay**

1/ Single
2/ Double

Manufactured to European Standards: BS 7177 Crib 5

Warranty 3 Years

Colour



Blue (PU) (Standard)

* Essential nursing care is pivotal in pressure ulcer prevention. This mattress will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action based planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. A range of assessment scales can be used as a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.



Invacare® | **Propad® Leg Trough**



Shaped leg support - provides comfort and promotes recovery

The *Invacare Propad Leg Trough* has been designed in consultation with clinicians to provide maximum support for the whole leg. The **Leg Trough** elevates the heels and allows them to be suspended over the edge of the trough. This provides maximum pressure relief for the vulnerable heel area.





Features and Options

Provides comfort and support to promote recovery

The leg trough offers optimum support and protection for the whole leg, reducing pressure and promoting recovery.

Relieves pressure on vulnerable heels

Elevates the heels and allows them to be suspended over the edge of the trough, providing pressure relief for the vulnerable heel area.



Heavy duty water-resistant cover

A two-way stretch, vapour permeable and anti-static polyurethane (PU) water-resistant cover.



Cover includes a zip and is removable

The cover can be unzipped and removed from the foam core to allow laundering of the cover.



Adjustable side-squeeze clip buckles

The adjustable buckles enable appropriate and comfortable positioning of the leg.

Technical data



Propad Leg Trough

Length: 510, 610, 660 mm



Width: 110 mm



Height: 115 mm



Weight: 0.8 kg



Cleaning: 80°

Manufactured to European Standards:
Interior Foam: BS 5852: 1982 Crib 5, BSEN597-1,2
Cover: BS7175: 1989 Crib 5

Colour



Blue

Warranty 1 Year



Invacare | **Cushions**




 Invacare® | **Cushions**
Cushions

<i>Invacare Flo-tech® Lite</i>	160-161
<i>Invacare Flo-tech® Contour</i>	162-163
<i>Invacare Flo-tech® Lite Visco</i>	164-165
<i>Invacare Flo-tech® Contour Visco</i>	166-167
<i>Invacare Flo-tech® Plus</i>	168-169
<i>Invacare Flo-tech® Image</i>	170-171
<i>Invacare Flo-tech® Solution</i>	172-173
<i>Invacare Flo-tech® Back Support</i>	174-175
<i>Invacare Flo-tech® Lumbar Pad</i>	176-177
<i>Invacare Propad® Original</i>	178-179
<i>Invacare Propad® Low Profile</i>	180-181
<i>Invacare Propad® Premier</i>	182-183
<i>Invacare Propad® Visco</i>	184-185

Invacare® | **Flo-tech® Lite**



3 Year warranty



Slimline contoured cushion offering comfort and improved posture

The *Invacare Flo-tech Lite* is a slimline contoured cushion that offers improved stability and postural support. This high quality, moulded foam cushion is designed to significantly improve client comfort whilst sitting.



Yes, you can.®



Slimline foam base made from high quality, resilient moulded foam

The high quality, resilient foam is manufactured using a unique formula, which offers high levels of pressure reduction and improves product longevity. Moulding allows the **Flo-tech Lite** cushion to provide a 3D form, which cannot be emulated in cut foam. The pommel at the front of the cushion separates the thighs, encouraging good posture and allowing air to circulate.

Contoured seating surface

Offers improved stability and postural support. In addition delivers a greater degree of weight distribution compared to that of a flat surface. The castellations of the cushion allow slight movement which in turn reduces shear and friction forces.

Features and Options



Water resistant cover

Cover is multi stretch, water-resistant and vapour-permeable. Cover features an anti-slip base and carry handle

Available in a wide range of sizes

To view the cushion size range, please visit www.invacare.co.uk.

Technical data



A range of 61 sizes



Height: 78 mm / 3"



Weight: 1.5 kg



Max. User Weight: 127 kg (20 stone)



Cleaning: 80°

Flo-tech Lite

Manufactured to European Standards: BS EN 597-1 Cig, BS EN597-2 Match, and BS 7176 (Crib 5)

Colours



Light Green



Black Smooth



Black Wipedown



Black Towelling



Black Spacer Fabric

Warranty 3 Years



Invacare® | **Flo-tech® Contour**



3 Year warranty



Designed for comfort and improved posture


The surface of the *Invacare Flo-tech Contour* cushion is contoured into sections. This unique design improves posture and increases the support surface area, effectively dispersing weight away from specific pressure points. The combination of improved posture and increased support surface area will contribute significantly to pressure ulcer prevention. This deep cushion is suitable for heavier clients (up to 30 stone), and those considered 'At Risk'* of developing pressure ulcers.



Yes, you can.®

Moulded foam cushion

Moulding allows **Flo-tech** cushions to provide a 3D form, which cannot be emulated in cut foam. Moulded foam is highly durable, increasing the product's longevity. The pommel at the front of the cushion separates the thighs, encouraging good posture and allows air to circulate.

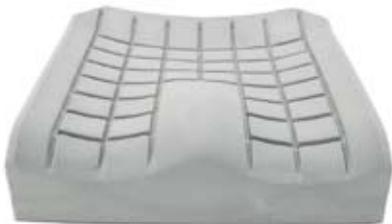

Contoured seating surface

Offers improved stability and postural support. The castellations of the cushion allow slight movement which in turn reduces pressure. The contoured surface helps to prevent the client sliding forward, therefore reducing the potential for shear and friction forces.

Features and Options



Water Resistant Cover

Cover is multi stretch, water-resistant and vapour-permeable. Cover features an anti-slip base.

Available in a wide range of sizes

To view the cushion size range, please visit www.invacare.co.uk.

Technical data



A range of 83 sizes



Height: 103 mm / 4"



Weight: 1.38 kg



Max. User Weight: 191kg (30 stone)



Cleaning: 80°

Flo-tech Contour

Manufactured to European Standards: BS EN 597-1 Cig, BS EN 597-2 Match, and BS 7176 (Crib 5)

Warranty 3 Years


Colours



Black Smooth



Black Wipedown

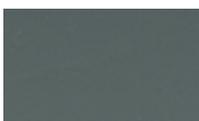

Black Towelling



Black Spacer Fabric



Kingfisher Green

* Essential nursing care is pivotal in pressure ulcer prevention. This mattress will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action based planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. A range of assessment scales can be used as a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.

Invacare® | **Flo-tech® Lite Visco**



3 Year warranty



Low profile visco elastic (memory) foam cushion offering exceptional comfort

The *Invacare Flo-tech Lite Visco* cushion is constructed from high density, visco elastic (memory) foam which is temperature sensitive, allowing it to conform to body shape. This gives exceptional levels of comfort whilst reducing peak pressures.



Yes, you can.®

Visco Elastic (memory) foam designed to maximise body contact area

The temperature-sensitive foam conforms to the body, maximising the contact area, delivering exceptional levels of comfort and optimising pressure reduction.

Contoured seating surface withommel and side support

Offers increased lateral support and pelvic stability to enhance the client's posture.




Features and Options



Waterproof resistant cover
Cover is multi stretch, water-resistant and vapour-permeable. Cover features an anti-slip base and carry handle.



Available in a wide range of sizes
To view the cushion size range, please visit www.invacare.co.uk.

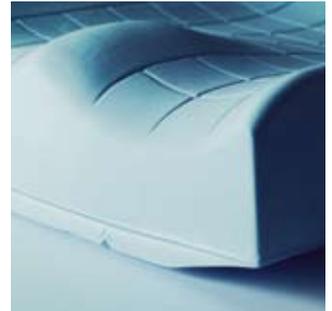

Three foam densities available
Available in three densities to accommodate a range of client weights:
Soft 5 stone (up to 32 kg),
Medium 8 stone (up to 51 kg)
and Firm 16 stone (up to 100kg).

Technical data


Flo-tech Lite Visco A range of 61 sizes	Height: 78 mm / 3"	Weight: 1.26 kg	Max. User Weight: 100 kg (16 stone)	Cleaning: 80°

Water Resistant
Vapour Permeable

Manufactured to European Standards:
Interior Foam: BS EN 597-1, BS EN 597-2, BS 5852:
1982 Crib 5
Cover: BS 7177: 1989 Crib 5

Colours

Warranty 3 Years

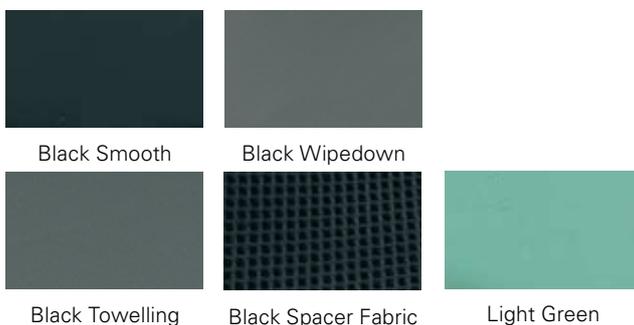

Invacare® | **Flo-tech® Contour Visco**



3 Year warranty



Conforms to the shape of the body

The *Invacare Flo-tech Contour Visco* cushion is moulded from visco-elastic (memory) foam, which is temperature-sensitive to the body. As the foam warms, the contoured surface area conforms to body shape, improving posture, increasing comfort and maximising the surface area in contact with the client, thereby optimising pressure reduction. The *Flo-tech Contour Visco* is suitable for those considered at 'High Risk'* of developing pressure ulcers.



Yes, you can.®

Deep moulded visco elastic (memory) foam cushion

The temperature-sensitive foam conforms to the client's body, maximising the contact area, delivering exceptional levels of comfort and optimising pressure reduction. The deep foam makes this cushion particularly suitable for heavier and taller clients.



Contoured seating surface withommel and side support

Offers increased lateral support and pelvic stability to enhance posture.

Features and Options



Water Resistant Cover

Cover is multi stretch, water-resistant and vapour-permeable. Cover features an anti-slip base and carry handle.

Available in a wide range of sizes

To view the cushion size range, please visit www.invacare.co.uk.

Four foam densities available

Available in four densities to accommodate a range of client weights -
 Soft: 8 stone (up to 51 kg),
 Medium: 12 stone (up to 76 kg),
 Firm: 16 stone (up to 102 kg) and
 Very Firm: 20 stone (up to 125 kg).

Technical data



Flo-tech Contour Visco
A range of 51 sizes



Height: 103 mm / 4"



Weight: 1.60 kg



Max. User Weight: 125 kg (20 stone)



Cleaning: 80°

Water Resistant
Vapour Permeable

Manufactured to European Standards:
 Interior Foam: BS EN 597-1, BS EN 597-2, BS 5852: 1982 Crib 5
 Cover: BS 7175: 1989 Crib 5

Warranty 3 Years

Colours



Black Smooth



Black Wipedown



Black Towelling



Black Spacer Fabric



Kingfisher Green

* Essential nursing care is pivotal in pressure ulcer prevention. This mattress will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action based planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. A range of assessment scales can be used as a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.



Invacare® | **Flo-tech® Plus**



3 Year warranty



Deep contoured cushion with gel sac

The *Invacare Flo-tech Plus* cushion incorporates a gel sac which is located towards the rear of the cushion. The positioning of the gel sac accommodates the Ischial Tuberosities, helping to significantly reduce pressure in this vulnerable area. This delivers superior levels of protection for those considered to be at 'High Risk'* of developing pressure ulcers.



Yes, you can.®



Features and Options

Deep contoured foam base with gel sac at rear of the cushion

The positioning of the gel sac provides increased immersion and protection to the Ischial Tuberosities which may provide protection for clients who are considered to be at 'High Risk' of developing pressure ulcers.

Contoured seating surface

Offers improved stability and postural support. Also delivers a greater degree of weight distribution compared to that of a flat surface. The castellations of the cushion allow slight movement which in turn reduces shear and friction forces. The deeper moulded foam offers greater protection to those who are considered clinically obese, or clients who would benefit from a higher sitting position.



Water Resistant Cover

Cover is multi stretch, water-resistant and vapour-permeable. Cover features an anti-slip base.

Available in a wide range of sizes

To view the cushion size range, please visit www.invacare.co.uk.

Technical data



A range of 80 sizes



Height: 103 mm / 4"



Weight: 2.5 kg



Max. User Weight: 191kg (30 stone)



Cleaning: 80°

Flo-tech Plus

Manufactured to European Standards: BS EN 597-1 Cig, BS EN 597-2 Match, and BS 7176 (Crib 5)

Warranty 3 Years

Colours



Fuchsia



Black Smooth



Black Wipedown



Black Towelling



Black Spacer Fabric

* Essential nursing care is pivotal in pressure ulcer prevention. This mattress will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action based planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. A range of assessment scales can be used as a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.



Invacare® | **Flo-tech® Image**



3 Year warranty

Slimline contoured cushion with gel sac

The *Invacare Flo-tech Image* cushion uses a combination of a slim, contoured foam with a large gel sac covering the sacral area. This increases even weight distribution and immersion of the Ischial Tuberosities. The design delivers a high degree of protection to clients who are considered to be at 'High Risk'* of developing pressure ulcers.



Yes, you can.®

Slimline contoured foam and gel combination cushion

Some clients prefer a slimline cushion, finding it more compatible with their daily living activities. The *Invacare Flo-tech Image* combines the benefit of a contoured foam upper surface with a large gel sac covering the ischial and sacral area, offering increased support and reducing the occurrence of peak pressures.

Contoured seating surface with gel sac

Offers improved stability and postural support. The moulded foam castellations and gel combine effectively to reduce shear and friction forces on tissue, and help to prevent the client sliding forward in the seat.



Features and Options



Water Resistant Cover

Cover is multi stretch, water-resistant and vapour-permeable. Cover features an anti-slip base and carry handle.

Available in a wide range of sizes

To view the cushion size range, please visit www.invacare.co.uk.

Technical data



A range of 50 sizes



Height: 78 mm / 3"



Weight: 1.82 kg



Max. User Weight: 140 kg (22 stone)



Cleaning: 80°

Flo-tech Image

Water Resistant
Vapour Permeable

Manufactured to European Standards:
BS EN 597-1 Cig, BS EN 597-2 Match, and BS 7176 (Crib 5)

Warranty 3 Years


Colours



Black Smooth



Black Wipedown

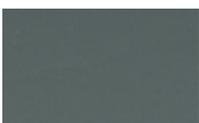

Black Towelling



Black Spacer Fabric



Mulberry

* Essential nursing care is pivotal in pressure ulcer prevention. This mattress will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action based planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. A range of assessment scales can be used as a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.



Invacare® | **Flo-tech® Solution**



3 Year warranty



Superior protection for 'Very High Risk' clients

The *Invacare Flo-tech Solution* cushion uses a combination of a slim, contoured foam base with a dual layer of gel sacs over the entire cushion surface. The dual gel sacs and the unique moulded foam shape of the *Flo-tech Solution* significantly increases the support surface area. This delivers a superior level of protection for those considered to be at 'Very High Risk'* of developing pressure ulcers.




Delivers superior protection from pressure ulcers

The *Invacare Flo-tech Solution's* design increases the contact area between the client and the foam cushion. The gel re-distributes pressure away from bony prominences, offering increased protection. The additional layer of gel is incorporated under the ischial and sacral area to further enhance the pressure reducing qualities of this cushion, making this cushion suitable for 'Very High Risk' clients.

Slimline contoured cushion

Some clients prefer a slimline cushion, finding it more compatible with their daily living activities. The low profile is particularly beneficial for the facilitation of side transfers. This may make this cushion more suitable for use with existing equipment.

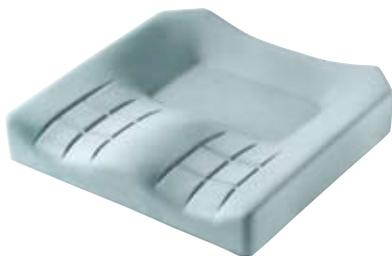

Features and Options



Water Resistant Cover

Cover is multi stretch, water-resistant and vapour-permeable. Cover features an anti-slip base and carry handle.

Available in a wide range of sizes

To view the cushion size range, please visit www.invacare.co.uk.

Invacare Flo-tech Deep Solution (Option)

Option of a 4" (10 cm) "Deep" foam base is also available in a range of sizes.

Technical data



A range of 44 sizes



Height: 78 mm and 102 mm (3" Standard, 4" Option)



Weight: 3.08 kg



Max. User Weight: 140 kg (22 stone)



Cleaning: 80°

Warranty 3 Years

Flo-tech Solution

Manufactured to European Standards: BS EN 597-1 Cig, BS EN 597-2 Match, and BS 7176 (Crib 5)

Colours



Black Smooth



Black Wipedown



Black Towelling



Black Spacer Fabric



Terracotta Swatch (Only available for welded wipe down covers)



Teal blue PU

* Essential nursing care is pivotal in pressure ulcer prevention. This mattress will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action based planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. A range of assessment scales can be used as a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.



Invacare® | *Flo-tech® Back Support*



3 Year warranty



The Flo-tech Back Support Range

The *Invacare Flo-tech Back Support* range has been designed to provide added stability to a wheelchair or standard armchair and it can be used in conjunction with any *Flo-tech* cushion.

The *Flo-tech Back Support* range provides excellent comfort and improved positioning, which in turn helps to improve or maintain independence whilst in the seated position. The unique supportive design offers additional support and may provide protection from pressure ulcer formation and spinal deformities.



Yes, you can.®



Features and Options

A highly protective back support

The **Flo-tech Back Support** provides a high degree of stability and provides improved posture resulting in enhanced protection from pressure ulcer development.

Added support of the lumbar region

The **Flo-tech Back Support** with Lumbar can be used for those who wish to have increased lateral and pelvic stability. The Visco Elastic (memory foam) Lumbar Support contours to the individual and can be easily adjusted and secured with velcro pads to accommodate a range of heights.



Standard cover is moisture resistant

The black towelling cover is two-way stretch, moisture-resistant and vapour-permeable, which helps to prevent a build up of shear and friction forces.



Extendable straps

Offers the client a securely fixed **Back Support** on a standard armchair or wheelchair. The straps are easy to adjust and can be fixed to a position that provides the client with the maximum comfort and support.



Complete Care - The Flo-tech Lumbar Pad

Designed for lumbar support, this pad provides increased lateral and pelvic stability, as well as helping to restore the natural lumbar curve.

Technical data



Flo-tech Back Support

Max User Weight: 165 kg




Cleaning: 80°

Back Support is manufactured to European Standards: BS 5852: 1990 Crib 5
 Back Support Lumbar is manufactured to European Standards: BS 5852: Crib 5

Colours



Black Smooth

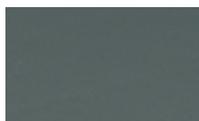

Black Towelling



Black Spacer Fabric

Warranty 3 Years



Invacare® | **Flo-tech® Lumbar Pad**



Improving posture, comfort and support

The *Invacare Flo-tech Lumbar Pad* is designed for use on existing chairs or wheelchairs and can improve posture, comfort and support. The *Lumbar Pad* is simple and easy to use. It is fitted with adjustable straps and comes complete with a black fabric cover.



Yes, you can.®

Visco Elastic (Memory) Foam

The **Lumbar Pad** is made up of two foams; Visco Elastic (Memory) foam on one side and resilient, high density foam on the other. This combination of foams provides the client with added support to the lower back. The resilient high density foam increases the support of the product whilst the Visco Elastic foam absorbs pressure from the client.



Features and Options



Adjustable Straps

Provides the client with improved comfort, adjustability and positioning.

Technical data



Flo-tech Lumbar Pad Size: 420 x 200 x 50 mm

Cleaning: 80°

The Lumbar Pad is manufactured to European Standards:
BS5852: 1990 Crib 5

Colour



Black Fabric

Warranty 1 Year

Invacare® | **Propad® Original Cushion**



Designed for comfort, support and pressure reduction

The *Invacare Propad Original* cushion has been designed to give comfort and pressure reduction. The *Propad Original* cushion is suitable for those considered 'At Risk'* of developing pressure ulcers.



Yes, you can.®



Features and Options

Deep profiled foam cushion, designed to deliver pressure care

The *Invacare Propad Original* Cushion is made from a layer of castellated foam. Each castellation moves independently to deliver enhanced support and pressure area care to the client, as well as good air circulation, reducing the potential for sweating and associated skin maceration.

Contoured seating surface

Offers improved stability and postural support; castellations reduce shear and friction forces on tissue.



Cover helps prevent a build up of shear and friction forces against the client

Cover is multi stretch, water-resistant and vapour-permeable. Cover features an anti-slip base.



Available in a wide range of sizes

To view the cushion size range, please visit www.invacare.co.uk.

Technical data



Wide range of sizes



Height: 102 mm / 4"



Weight: 1 kg



Max User Weight: 110kg / 17.3 stone



Cleaning: 80°

Propad Original Cushion


Manufactured to European Standards:
 Interior Foam: BS EN 597-1 Cig, BS EN597-2 Match, BS 5852 (Crib 5)
 Cover: BS 7175 (Crib 5)

Warranty 1 Year

Colours



Blue

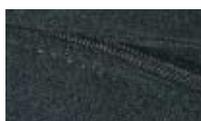

Black Towelling

* Essential nursing care is pivotal in pressure ulcer prevention. This product will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action based planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. A range of assessment scales can be used as a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.



Invacare® | ***Propad® Low Profile Cushion***



Comfort and support in a low profile

The *Invacare Propad Low Profile* cushion provides clients who require a lower seating position the option of having a lower height cushion. This cushion is suitable for those who find a low profile cushion more compatible with their daily living activities.



Yes, you can.®



Features and Options

Castellated Foam Cushion

The *Invacare Propad Low Profile cushion* is made from a layer of castellated foam, each castellation moves independently to deliver support and pressure area care to the client, as well as good air circulation, reducing the potential for sweating and associated skin maceration. This cushion is suitable for 'Low Risk' clients who require cushioning for comfort.

Reduces shear and friction forces

The castellations in the *Invacare Propad Low Profile cushion* reduce shear and friction forces on tissue.



Cover helps to prevent a build up of shear and friction forces

Cover is multi stretch, water-resistant and vapour permeable. Cover features an anti-slip base.



Available in a wide range of sizes

To view the cushion size range, please visit www.invacare.co.uk.

Technical data



A range of 102 sizes



Height: 65 mm / 2"



Weight: 0.8kg



Max User Weight: 76 kg (12 stone)



Cleaning: 80°


Propad Low Profile Cushion

Manufactured to European Standards: BS5852: Crib 5
BS7175: Crib 5

Colours



Blue PU

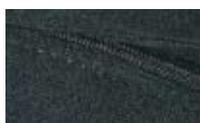

Black Towelling

Warranty 1 Year

Invacare® | **Propad® Premier Cushion**



3 Year warranty



Designed for comfort, support and pressure reduction

The *Invacare Propad Premier Cushion* has been designed to give comfort and pressure reduction. The cushion has a dual layer construction to ensure improved support as well as prolonged longevity. The **Propad Premier Cushion** is suitable for those considered 'At Risk'* of developing pressure ulcers.



Yes, you can.®



Features and Options

Dual-layer profiled foam cushion, designed to deliver enhanced pressure care

The lower layer of the *Invacare Propad Premier Cushion* is constructed from a resilient high density foam and is shaped to deliver increased stability and support for the client. The upper section of the cushion features a layer of castellated foam. Each castellation moves independently to deliver enhanced support and pressure area care to the client, as well as ensuring good air circulation, reducing the potential for sweating and associated skin maceration.



Cover helps to prevent a build up of shear and friction forces

Cover is multi stretch, water-resistant and vapour-permeable. Cover features an anti-slip base.

Available in a wide range of sizes

To view the cushion size range, please visit www.invacare.co.uk.

Technical data



Wide range of sizes



Height: 102 mm / 4"



Weight: 1 kg



Max User Weight: 110kg / 17.3 stone



Cleaning: 80°

Propad Premier Cushion

Manufactured to European Standards: BS EN 597-1 Cig, BS EN 597-2 Match, and BS 7175 (Crib 5)

Warranty 3 Years

Colours



Blue PU



Black Towelling



Black Smooth

* Essential nursing care is pivotal in pressure ulcer prevention. This product will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action based planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. A range of assessment scales can be used as a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.



Invacare® | **Propad® Visco Cushion**



Dual layered cushion - conforms to the shape of the client

The *Invacare Propad Visco cushion* is a dual layered cushion, using a combination of high density, temperature-sensitive, visco-elastic (memory) foam in conjunction with a highly resilient foam base, to provide added support and longevity. The **Propad Visco cushion** is suitable for those considered 'At Risk'* of developing pressure ulcers.



Yes, you can.®

Visco Elastic (Memory) Foam Top

The top foam layer of the cushion is made from high density visco elastic (memory) foam, which conforms to the shape of the body, combining exceptional comfort with a fully supported feeling.

Robust Foam Base

The high density foam base, provides added support, and prevents the cushion from bottoming out, ensuring longevity of the product.



Features and Options



Cover helps to prevent a build up of shear and friction forces against the client

Cover is multi stretch, water-resistant and vapour permeable. The cover also features an anti-slip base.

Available in a range of sizes

To view the cushion size range, please visit www.invacare.co.uk.

Technical data



A range of 3 sizes



Height: 102 mm / 4"



Weight: 1.4kg



Max User Weight: 110kg / 17.3 stone



Cleaning: 80°

Propad Visco Cushion

Manufactured to European Standards:
BS 5852: Crib 5
BS 7175: Crib 5

Warranty 1 Year

Colours




Blue PU



Black Smooth



Black Towelling

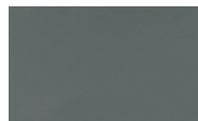

Black Wipedown

* Essential nursing care is pivotal in pressure ulcer prevention. This product will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action based planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. A range of assessment scales can be used as a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.

Technical Information

- Pressure Area Care Products

Product	Warranty	Fire Testing	Grade Ref & Colour
<i>Invacare Softform®</i>			
Premier Active	5 Years cover 10 Years foam 2 Years insert 2 Years pump	BS 5852: Crib 5	RX39/120 Pink RX39/200 Blue
Premier MaxiGlide	4 years cover 8 years foam	BS 5852: Crib 5 BS 7175: Crib 5	GR. RX 39/120 Pink GR. RX 39-200 Blue
Premier	5 Years cover 10 Years foam	BS 7177: Crib 5	Gr. 39/120 Pink Gr. RX 39/200 Blue
Premier Visco	4 Years *	BS 5852: Crib 5	Gr. RX 39/120 Pink Gr. RX 39/200 Blue Gr. VC5580 White
Excel	2 Years	BS 7177: Crib 5	Gr. RX 39/120 Pink Gr. RX 39/200 Blue
Heel Pad	1 Year	BS 5852: Crib 5	Gr. RX 33/125 White
Odstock Wedge	1 Year	BS 5852: Crib 5	Gr. RX 33/125 White
<i>Invacare® Propad®</i>			
Premier Overlay	3 Years**	BS 7177: Crib 5	Gr. RX36/125 Beige
Original Overlay	3 Years**	BS 7177: Crib 5	Gr. RX36/125 Beige
Visco Overlay	1 Year	BS 5852: Crib 5	Gr. VC5580 White Gr. RX36/125 Beige
Premier Cushion	1 Year	BS 5852: Crib 5 BS 7175: Crib 5	Gr. RX36/125 Beige Gr. RX39/200 Blue
Original Cushion	1 Year	BS 5852: Crib 5 BS 7175: Crib 5	Gr. RX36/125 Beige
Low Profile Cushion	1 Year	BS 5852: Crib 5 BS 7175: Crib 5	Gr. RX36/125 Beige
Visco Cushion	1 Year	BS 5852: Crib 5 BS 7175: Crib 5	VC5580 RX39/200 Blue
Leg Trough	1 Year	BS 5852: Crib 5	RX39/200 Blue
<i>Invacare Flo-tech®</i>			
Flo-tech® Cushions	3 Years**	BS 7176: 1995 Crib 5	Gr. FG 70/180 Grey
Flo-tech® Visco Cushion: Soft	3 Years**	BS 5852: 1990 Crib 5	Gr. LM 85/115 Blue
Flo-tech® Visco Cushion: Medium	3 Years**	BS 5852: 1990 Crib 5	Gr. LM 85/150 Blue
Flo-tech® Visco Cushion: Firm	3 Years**	BS 5852: 1990 Crib 5	Gr. LM 85/190 Blue
Flo-tech® Visco Cushion: Very Firm	3 Years**	BS 5852: 1990 Crib 5	Gr. LM 85/250 Blue
Flo-tech® Back Support	3 Years**	BS 5852: 1990 Crib 5	Gr. FG 70/180 Grey
Flo-tech® Back Support Lumbar	3 Years**	BS 5852: 1990 Crib 5	Gr. FG 70/180 Grey VC55085 White
Lumbar Pad	1 Year	BS 5852: 1990 Crib 5	Gr. RX 39/200 Blue VC55085 White



Yes, you can.®

Product	Nominal Density Range (kg/m)	Nominal Hardness Range (N)	Fatigue Class	Weight of Product (kg)
<i>Invacare Softform®</i>				
Premier Active	38 - 40	110 - 130	Very Severe	15
	38 - 40	175 - 215	Very Severe	
Premier MaxiGlide	38 - 40	105 - 135	Very Severe	13.5
	38 - 40	180 - 200	Very Severe	
Premier	38 - 40	105 - 135	Very Severe	13
	38 - 40	180 - 200	Very Severe	
Premier Visco	38 - 40	105 - 135	Very Severe	13
	38 - 40	180 - 200	Very Severe	
	55	80	N/A	
Excel	35 - 37	105 - 135	Very Severe	13
	38 - 40	180 - 200	Very Severe	
Heel Pad	32 - 34	110 - 140	Severe	3
Odstock Wedge	32 - 34	110 - 140	Severe	2,5
<i>Invacare® Propad®</i>				
Premier Overlay	35 - 37	110 - 140	Very Severe	5,5
Original Overlay	35 - 37	110 - 140	Very Severe	5,5
Visco Overlay	55	80	N/A	7,1
	35 - 37	110 - 140	Very Severe	
Premier Cushion	35 - 37	110 - 140	Very Severe	1
	38 - 40	180 - 200	Very Severe	
Original Cushion	35 - 37	110 - 140	Very Severe	1
Low Profile Cushion	35 - 37	110 - 140	Very Severe	1
Visco Cushion	55	80	N/A	1,4
	38 - 40	180 - 200	Very Severe	
Leg Trough	38 - 40	180 - 200	Very Severe	0,75
<i>Invacare Flo-tech®</i>				
Flo-tech® Cushions	65 - 75	170 - 190	Very Severe	
Flo-tech® Visco Cushion: Soft	68 - 78	90 - 120	N/A	
Flo-tech® Visco Cushion: Medium	68 - 78	120 - 150	N/A	
Flo-tech® Visco Cushion: Firm	68 - 78	155 - 195	N/A	
Flo-tech® Visco Cushion: Very Firm	68 - 78	195 - 235	N/A	
Flo-tech® Back Support	65 - 75	170 - 190	Very Severe	
Flo-tech® Back Support Lumbar	65 - 75	170 - 190	Very Severe	
	50 - 60	55 - 85	N/A	
Lumbar Pad	38 - 40	175 - 215	Very Severe	
	50 - 60	55 - 85	N/A	

Invacare® | Pressure Ulcer Risk Assessment

When assessing an individual for a pressure reducing cushion or mattress it is always advisable to carryout some form of risk assessment. This risk assessment can include a formal risk assessment scale such as the Norton scale (1962), the Braden scale (1985) and the Waterlow scale (1985). The Waterlow scale for example provides a guide category of potential risk for that individual, based on the final score.

Category	Waterlow Score
At Risk	10+
High Risk	15+
Very High Risk	20+

However, it is recommended that the results of a formal risk assessment are always considered alongside more informal clinical judgement before the total care package is decided, of which pressure management equipment is an important part. This ensures that factors that may not be included in the formal risk assessment are given sufficient consideration.

In addition, Invacare recommend that guidance on pressure ulcer risk assessment is sought from relevant Agencies and Panels as detail below.

References

Bergstrom N., Braden, B., Laguzza A., Holman, V. The Braden Scale for predicting pressure sore risk: reliability studies. Nurs Res 1987; 34(6);205-210.

European Pressure Ulcer Advisory Panel and National Pressure Ulcer Advisory Panel. Prevention and treatment of pressure ulcers: quick reference guide. Washington DC: National Pressure Ulcer Advisory Panel; 2009.

National Institute for Clinical Excellence. Inherited Clinical Guideline B. Pressure Ulcer Risk Assessment and Prevention. National Institute for Clinical Excellence, April 2001.

Norton, D., Exton Smith A.N., McLaren R. An investigation of geriatric nursing problems in hospital. London: National Corporation for the Care of Old People. 1962.

Waterlow, J. A risk assessment card. Nurs Times 1985; 81 (49);5115

Waterlow, J. The Waterlow pressure ulcer prevention manual. Taunton: Newtons. 2005.



	Product	Waterlow Score
Mattress	Softform Premier	20+
	Softform Premier MaxiGlide	20+
	Softform Premier Active	20+
	Softform Premier Visco	15+
	Softform Excel	10+
Overlay	Propad Overlay	15+
	Propad Premier Overlay	15+
	Propad Visco Overlay	15+
Cushion	Flo-tech Solution	20+
	Flo-tech Contour Visco	15+
	Flo-tech Plus	15+
	Flo-tech Image	15+
	Flo-tech Contour	10+
	Propad Original	10+
	Propad Premier	10+
	Propad Visco	10+
	Flo-tech Lite Visco	Up to 15
	Propad Low Profile	Up to 10
Flo-tech Lite	Up to 10	

Cleaning Instructions

General Cleaning – to prevent delamination, use a solution of neutral detergent and hot water, or alternatively a 0.1% hyperchloride solution. Apply with single use wipe. Following application, the cover should be rinsed with clean water (damp cloth) and thoroughly dried.

Removable contaminated covers should be laundered at 80°C with a dilute detergent solution. Please ensure covers are thoroughly rinsed and dried before placing back on the Invacare product.

(Invacare recommends that care staff should always actively liaise with a person qualified in infection Control regarding decontamination protocols.)

NB: Always ensure you read the user instructions prior to use, and inspect the interior and exterior of cushions, mattresses, overlays and accessories on a regular basis.

* Standard Configuration

** Reducing Warranty, please contact Invacare for full details.

Invacare® | Pictograms with Description

Width between castors	Height between floor and seat	Total outside width	Height from floor to headrest	Length from front castor to footrest	Maximum user weight	Tilt angle	Back rest recline angle
Weight of shower chair	Minimum & maximum seat height	Total chair height	Seat width without side flaps	Seat width with side flaps	Seat depth	Length of base	Width of base / footprint
Total length (backrest reclined)	Width of widest point on backrest	Backrest height	Tilt angle	Maximum user weight	Total weight of bathlift without battery	Weight of battery	Number of cycles
Inside width when hoist is closed	Inside width when hoist is open	Length of boom	Height of legs	Total product weight	Maximum user weight	Floor to boom height at highest position	Length of motor
Width of motor	Height of motor	Strap length	Strap length increments	Width outside: width inside	Width outside	Length outside: length inside	Height adjustable



Mattress support dimensions



Backrest angle



Overall width: lying surface:



Overall length: lying surface



Thigh angle



Legrest lift height



Tilt option



Thigh angle legrest lift



Legrest angle



High-Low: total height:



Under bed clearance



(Angle below frame)
(Angle above frame)



Minimum angle below horizontal



Full legrest lift



Total product weight



Anti-static wheels with independent braking



Various sizes



Height



Width



Depth



Height



Depth



Width



Various sizes



Width



Height



Depth



Length



Width



Height



Length



Height



Width



Width outside wheels



Width between wheels



Width between handles



Product weight



User weight limit



Depth



Overall height



Seat height



Width folded



Wheel diameter



Depth folded