

Invacare® | **Action<sup>1</sup><sub>NG</sub>**



### **Simple, robust and efficient**

The Invacare **Action<sup>1</sup><sub>NG</sub>** steel wheelchair is the latest addition to the **Action<sub>NG</sub>** Family. Although made of steel, this remains a lightweight wheelchair of exceptional quality. The rigidizer blocks, commonly used with the **Action<sub>NG</sub>** Family range, make the product extremely stable and, consequently, very easy to propel. The **Action<sup>1</sup><sub>NG</sub>** is available in fixed configurations for immediate delivery, but Invacare also offers a wide range of accessories to meet the differing clinical needs of users and carers for optimal, everyday, use.

So if you're looking for a high quality wheelchair that's practical, economical and convenient, then look no further than the **Action<sup>1</sup><sub>NG</sub>**.



**Yes, you can.®**



### Cross compatible to suit every need

The *Invacare Action<sup>1</sup><sub>NG</sub>* shares multiple components with other wheelchairs from the **Action<sub>NG</sub>** Family, including brakes, wheels, legrest assortment, cane holder and drip stand. That means you can adapt your **Action<sup>1</sup><sub>NG</sub>** to fully meet the user needs. Enabling the sharing of common components reduces the complexity of handling, prescribing and adjusting the product.

### Manoeuvrable, and easy to transport and store

The **Action<sup>1</sup><sub>NG</sub>** has been design for optimal manoeuvrability, whether being pushed or self-propelled, indoor or out. It even remains mobile in tight spaces or crowded areas because of its modest dimensions and compact frame. Plus, thanks to it's single crossbar system, compactness and quick release rear wheels axles, the **Action<sup>1</sup><sub>NG</sub>** is easy to handle, transport and store. Simply take off the rear wheels and fold it away.

### Simple to customize

The **Action<sup>1</sup><sub>NG</sub>** is available in both transit and self propulsion versions to accommodate individual needs. The wheelchair can easily be adjusted to

the user morphology. A wide range of accessories is available to enhance patient's comfort and customize the chair to their specific requirements.

## Features and Options



#### Anti-tippers

Provide extra stability on slopes or uneven ground.



#### Attendant hub brakes

Attendant operated hub brakes are available and include push to lock brakes.



#### Anti-barrier

Allow users to pass through narrow doors by removing the rear wheels.

**Features and Options**



**Transit version**

A transit version is available with 12" solid rear wheels for dependent users.



**Legrest**

Different types of legrest are available : amputee support ①, elevating legrest 0° to 80° ② and fixed 80° legrest ③, to meet every user needs.



**Seat belt**

A buckle belt is available to increase safety while seated.

**Dual Handrim version**

An hemiplegic version is available with Dual HR system. The dual handrim can be mounted on either the right or left side.



**Headrest**

Can be adjusted in both height and depth on the backbrace to bring more comfort to the user.



**Front castor**

2 types of front castor are available : 8" x1.3/8" solid ① and 8" x 2" solid ②.



**Armrest assortment**

2 types of swing away armrest are available for easy transfer : ① fixed height desk (240 mm) and ② fixed height full length (340 mm).



**Cane holder**


Hold the cane securely to the chair.

## Features and Options



### Brake lever

The folding lever of push to lock brakes makes transfer easier.

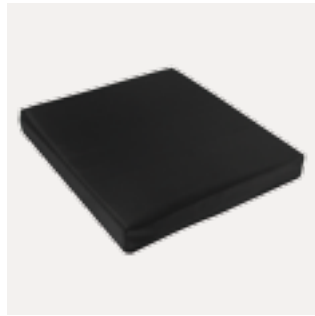

### Cushion

A 50 mm standard foam cushion with nylon and anti-slip cover is available.



### Drip stand

## Technical data

<b>Self propulsion</b>	405 / 455 / 505 mm	425 mm	3°	450 mm	7°	410 mm
<b>Transit</b>	405 / 455 / 505 mm	425 mm	3°	450 mm	7°	410 mm
<b>Self propulsion</b>	330 - 455 mm	240 mm	Seat width + 215 / 205 mm**	315 / 305 mm**	1045 / 1085 mm**	815 / 855 mm**
<b>Transit</b>	330 - 455 mm	240 mm	Seat width + 175 mm	265 mm	925 mm	695 mm
<b>Self propulsion</b>	22°	880 mm	127 kg	16,5 / 18,6* / 18,8 kg**	10,5 / 11,5* / 12,9 kg**	
<b>Transit</b>	22°	880 mm	127 kg	16,3 kg	14,6 kg	

\* With Dual Handrim  
\*\* With drum brakes

All product dimensions are in millimetres to apply EN 12183 standards. For more pre-sales information about this wheelchair, please consult the [www.invacare.eu.com](http://www.invacare.eu.com) website. It can be found under the section Low Active Wheelchairs / Action1NG / Sales & Images in the product catalogue. For information about typical user description, intended use, adjustments and transportation, please refer to the user guide.

Invacare Limited  
Pencoed Technology Park - Pencoed -  
Bridgend CF35 5AQ - United Kingdom  
Tel: +44 1 656 776 222  
Fax: +44 1656 776 220  
E-mail: [uk@invacare.com](mailto:uk@invacare.com)  
Sales Order E-mail:  
[ordersuk@invacare.com](mailto:ordersuk@invacare.com)  
[www.invacare.co.uk](http://www.invacare.co.uk)

Invacare Ireland Ltd.  
Unit 5, Seatown Business Campus  
Seatown Road, Swords,  
Co. Dublin - Ireland  
Tel: +353 (0)1 810 7084  
Fax: +353 (0)1 810 7085  
[ordersireland@invacare.com](mailto:ordersireland@invacare.com)  
[www.invacare.ie](http://www.invacare.ie)

© 2013 Invacare International Sàrl  
All rights reserved.

All information quoted is believed to be correct at time of print. Invacare reserves the right to alter product specifications without prior consultation.

A1NG - UK - 01/2013

## Frame colours



Night Blue



Black Nylon

Please notice colors may vary slightly from those displayed above.



**Yes, you can.®**