

DOLOMITE | *STEP UP*



The Optimum Rehabilitation Choice

The **DOLOMITE STEP UP** is a modern platform walker which maximises gait training techniques for rehabilitation purposes. It gives the patient confidence to start walking again following surgery or recovering from an injury. It is also suited to helping those who have long term problems with their walking. The **STEP UP** is available in three height adjustable versions - manual, hydraulic and electric. Assembly of accessories is tool free.





Ergonomic design

The single central column allows the carer easy access to help the patient and the simple adjustments mean the **STEP UP** can be set to suit specific needs. This means excellent manoeuvrability whilst continuously promoting an upright posture. The base is raised at the front to allow a longer step, while the lowered rear end helps with access under beds, etc. The unique fixed angle of the arm supports give a 25° angle between the shoulder point and wrist, providing the most natural position, maximising upper body strength and keeping elbows close. The platform is tilted 5° giving further support, and the 110° angled hand grip allows for a natural wrist position.

Fully Adjustable


The **STEP UP** is fully adjustable. And enables for the adjustment in the width of the base and arm supports, the height of the central column and the depth and angle of hand grips. The arm pads have been made using a soft material for pressure relief and are securely attached with microfilm. The arm pads are simple to remove for cleaning or replacing.

Transfer and access

The **STEP UP** is well thought out for transfer situations. The lowered rear end helps to reverse under the bed, whilst

the generous width allows wheelchairs to access the frame.

Features and Options

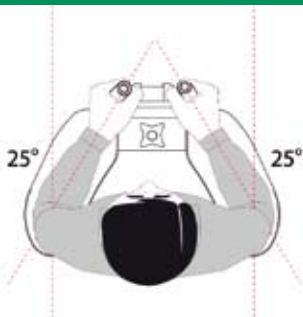

Release the central knob to adjust the width of arm supports.



The width of the legs on the **STEP UP** can be adjusted.

Accessories and options

Quick and easy assembly in track or prepared recess.

**Practical Basket**

Keeps belongings within easy reach.
Part No 1531390

**Tray**

To carry those small items.
Part No 1531391

**Elbow Stop**

Secures elbows from dropping off the side.
Part No 1531448

**Swivel Castor with directional brakes**

Helps focus on walking instead of direction.
Part No 1531408

**Brake System**

For those capable of operating a brake function.

Hand-brake kit

Part No 1537656

One-hand brake

Part No 1531450

**Oxygen tank holder**

Securely mount an oxygen tank on the walker.
Part No 1531406

**Seperate IV Holder**

Separate IV-stand makes life complicated. Not anymore.
Part No 1531407

**Stand assist incl. sling**

For prescribed standing periods
Fits hydraulic and manual version.

Part No 1531449 (standing plate + knee support + sling)

Standing plate only

Part No 1535330

**Ball grips**

Suitable for those who have difficulty closing their hand.

Others

Small castor 100 mm diameter.

Models**Manual**

Pull out the pin, move to required height, release the pin.

**Hydraulic**

Pull up the lever to adjust to required height with hydraulic assistance.

**Electric**

Using the controller, set the height with the electric motor.

Technical data

Step Up Manual	58/67/76 cm	70/79/88 cm	87 - 132 cm	34 - 47 cm	62 - 75 cm
Step Up Hydraulic	58/67/76 cm	70/79/88 cm	87 - 132 cm	34 - 47 cm	62 - 75 cm
Step Up Electric	58/67/76 cm	70/79/88 cm	87 - 132 cm	34 - 47 cm	62 - 75 cm



Step Up Manual	81 cm	16 kg	150 kg	125 mm	215 mm
Step Up Hydraulic	81 cm	18 kg	150 kg	125 mm	215 mm
Step Up Electric	81 cm	23 kg	150 kg	125 mm	215 mm

Colour

RAL 1035

Article code

Step Up Manual 1521677
Step Up Hydraulic 1521678
Step Up Electric 1521679

INVACARE® Limited
 Pencoed Technology Park - Pencoed -
 Bridgend CF35 5HZ
 Tel: +44 1 656 776 222
 Fax: +44 1656 776 220
 E-mail: uk@invacare.com
www.invacare.co.uk
 Sales Order E-mail:
ordersuk@invacare.com

INVACARE® Ireland Ltd
 Unit 5 - Seatown Business Campus
 Seatown Road - SWORDS - County
 Dublin
 Tel. +353 1 8107084
 Fax +353 1 8107085
 Email: Ireland@invacare.com
www.invacare.ie

© 2010 Invacare® International Sarl
 All rights reserved.

All information quoted is believed to be
 correct at time of print.
 Invacare® reserves the right to alter
 product specifications without prior
 consultation.

Step Up - UK - 10/2010 - 1538648



Yes, you can.®