

Invacare® | Storm⁴ Max

NEW



Much more than just a larger power chair!

The new **Storm⁴ Max** is the latest addition to the **Storm** range, designed for heavier clients. It features a longer chassis to improve weight distribution yet still provide excellent manoeuvrability outdoors or in. Extra long, sturdy armrests and legrests will support the client, especially when re-positioning themselves within the seat.

For the ultimate in driving performance and traction control, Invacare's **G-Trac™** technology is also available as an option.

Storm⁴ Max offers the same trouble free servicing that is attributable to the **Storm** range.



Yes, you can.®



Features and Options

Tilt and recline

For individual comfort and positioning the **Storm⁴ Max** offers powered tilt and recline options.

The backrest recliner has a range of 0° - 25° with a high pivot point in order to reduce shearing forces.

Ease of use

Just as the entire **Invacare Storm⁴** family is characterised by its ease of use, so is the **Max**. Batteries and electronics are easily accessible, making servicing simple. The depth of the seat can be simply adjusted in 3 steps which are indicated on the seat frame. An adapter to further reduce the seat depth is offered as an option.



Armrest

The strong and stylish armrest can be mounted in several positions and are width adjustable.



Footboard

Flip-up footboard supports various positions of the feet. It's strong enough to allow clients to reposition themselves within the seat.




Calf pads

The extra wide Max calf pads are individually adjustable in depth and angle to give optimum support.

Technical data

Storm⁴ Max	490 - 590 mm 550 - 650 mm 610 - 710 mm	560, 595, 630 mm ¹	480, 510 mm	570 mm ^{2,3} , 630 mm ³	360 - 480 mm	0° - 16° ⁴ , 0° - 9° ⁵	0 - 25°	1390 mm	2080 mm	
Storm⁴ Max	60 mm	6°(10,5 %)	350W	31 km ⁶	6 km/h	220 kg				

Standard colours

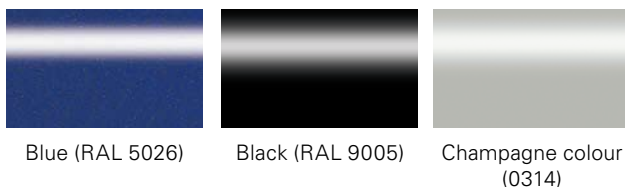

- 1/ Adapter to reduce seat depth by 130 mm available
- 2/ For small seat frame (490 mm) only
- 3/ Additional 50 mm of height adjustment
- 4/ Powered seat tilt
- 5/ Man. setttable seat tilt
- 6/ according to ISO 7176-4:2008

Colour options



Invacare Limited
Pencoed Technology Park - Pencoed -
Bridgend CF35 5AQ - United Kingdom
Tel: +44 1 656 776 222
Fax: +44 1656 776 220
E-mail: uk@invacare.com
Sales Order E-mail:
ordersuk@invacare.com
www.invacare.co.uk

Invacare Ireland Ltd
Unit 5 - Seatown Business Campus
Seatown Road - SWORDS - County
Dublin - Ireland
Tel. +353 1 8107084
Fax +353 1 8107085
Email: Ireland@invacare.com
www.invacare.ie

© 2011 Invacare International Sàrl
All rights reserved.

All information quoted is believed to be correct at time of print.
Invacare reserves the right to alter product specifications without prior consultation.

Storm4 Max - UK - 09/2011 - 1551547



Yes, you can.®



More detailed information, technical data, references for transport or description of certain options can be found in the wheelchair's user manual.