

**Micro, light, portable: The Invacare<sup>®</sup> Lynx**



**Yes, you can.<sup>®</sup>**

## Invacare® **Lynx**: Probably the most portable micro scooter in the world



### Simple to Transport

Designed for easy transport, the **Lynx** simply fits in any car boot. You will be amazed by how quick and simple it is to dismantle and reassemble. No need for tools, no confusing cables or fiddly parts: simply made of 4 main pieces, it offers a trouble free dismantling. You will be ready to go before you know it!





### **At Ease Anywhere**

As the ideal companion for shopping and daily activities, the *Invacare® Lynx* takes you anywhere you want. Very easy to manoeuvre, the portable micro scooter can also be discreet, especially when space is tight.

### **Compact but Powerful**

Compact and elegant, the *Lynx* has enough power and capacity to convince you. In utmost comfort, the *Lynx* will drive you quickly and safely at a speed of 6km/h. Designed to withstand up to 135 kg user weight, it is a strong and robust scooter.

### **Easy to Manage, Easy to Use**

With certainty, you will love its many advantages. A comfortable padded seat that swivels 360° and is adjustable in height, an adjustable tiller and armrests which are adjustable in width: the *Lynx* is user-friendly. Thanks to its straightforward control panel, it is easy to manoeuvre: you will enjoy a smooth ride !



**For all these reasons,  
the *Invacare® Lynx*  
is possibly the best  
portable micro scooter  
of its generation.**

## Safety Features

### Control panel

- Easily visible intuitive keypad
- Battery level indicator



### Reverse audible alert

When reversing the scooter will automatically alert people behind you by giving an audible warning signal (beep).

### User weight

The Invacare® **Lynx** has been tested to withstand a maximum user weight of 135 kg. It is a strong and robust scooter ensuring greater durability.

### Dependable electrics

Only high quality wiring and electric connections are used in order to optimize reliability.

### Anti-corrosion chassis

Anti-corrosion treated chassis for extended life of the scooter.



## Accessories



### Crutch & Cane holder

Keeps your crutch or cane in place during trips



### Rear basket

Holds your extra shopping



### Storage cover

Protects from dirt and water during storage



# Safe, solid, secure: *The Invacare® Leo*



## **ROBUST, RELIABLE AND EASY-TO-USE**

With unrivalled levels of safety and solidity, “**Leo**” will quickly become your reliable companion for your essential daily outings, leisure excursions and relaxing breaks away from home. With a robust chassis, powerful motor, intelligent electronics and distinctive looks.



**Yes, you can.®**



**Leo** is the best choice of mini scooter in its class currently on the market. **Leo** will help you rediscover your independence, giving you the freedom to travel wherever and whenever it suits you – in comfort, safety and in style.

Planning a short break or weekend away? **Leo** conveniently dismantles into five separate parts. **Leo** neatly fits in the boot of most cars and can be easily dismantled and re-assembled without tools.



**PERFECT CHOICE**

A new range of user-friendly features makes **Leo** the perfect choice for anyone seeking a reassuringly safe and stable drive. With such high standards of design and technical engineering, we are confident that no other mini scooter can currently match **Leo**'s wide array of safety features.

# THE INVACARE LEO OFFERS MANY SOPHISTICATED SAFETY FEATURES:



## Rear brake light

This new feature warns pedestrians that you are slowing down, thus preventing collisions.

## Intelligent lighting system

Should a light bulb burn out, the effect remains localised. This ensures maximum visibility of the scooter as all other lights still function normally.

## Dependable electrics

For maximum reliability, only the best-quality wiring and electrical connections are used onboard the **Leo** scooter:

## 2-step disengaging lever

The disengaging lever protects your scooter from being accidentally knocked into freewheel mode. This key safety feature prevents any unintentional disengaging of the drive mode.

## Anti-corrosion chassis

The chassis is finished with an anti-corrosion treatment, protecting against rust and prolonging the life of the scooter:

## Reliability guaranteed

Maintenance and easy serviceability were the key criteria throughout the entire **Leo** design process. To ensure high levels of reliability, several features such as anti-splash guards, quality cabling, fuses and connectors are included as standard to protect the electronics and motor from water damage and corrosion; the main causes of failures and break downs. Simply removing the battery cover offers direct access to the scooter batteries and electronics in case of service or maintenance.

## Intuitive control panel

**Leo** features an easily visible, intuitive keypad with domed buttons for accurate operation. Visual LED indication and audio feedback indicates which buttons have been activated.

## Low battery warning

You are alerted to a low battery level in two ways; a visual warning (flashing light) and an audible alert (beep) informs you when the battery needs recharging.

## Strength and durability

The **Leo** scooter has been tested to a maximum user weight of 136 kg, offering optimum durability and strength.

## Full lighting system

Modern 'Xenon'-style headlights provide increased visibility and reduced power consumption compared to traditional scooters, allowing for longer periods between battery charging. Integrated indicator lights inform road users and pedestrians that you are about to turn. Indicators automatically switch off after 30 seconds should you forget to turn them off yourself.





### **Mirror**

Moveable mirror can be mounted to the left or right.



### **Crutch & Cane holder**

Keeps your crutch or cane stable during the drive.



### **Oxygen bottle holder**

Safe storage of oxygen during transit.



### **Seat bag**



### **Lockable front box**

Store personal valuables safely while driving.



### **Walking frame holder**

Easy storage for your walking frame.



### **Off-board charger kit**

Once removed, the batteries can be charged away from the scooter.



### **Seat belt**

Ensures you are kept securely and comfortably in place.



### **Rear basket**

Perfect for bulky or extra shopping.



### **Storage cover**

Protects from dirt and water during storage.

**ACCESSORIES**



# Stable, Trusty, Durable: The *Invacare*® Orion



Imagine being able to undertake long journeys in complete safety and comfort. Whether you enjoy a daily trip around the shops or love to travel further afield, the *Invacare*® **Orion** gives you freedom of independence. Large batteries ensure plenty of capacity to enjoy excursions away from home. Whatever the terrain, you'll enjoy a controlled and sturdy drive: let the powerful motor, fully suspended chassis and big wheels take you smoothly over obstacles or uneven ground.



**Yes, you can.®**





## RELIABLE AND EASY TO MAINTAIN

The **Orion** is equipped with multiple features to ensure a safe, smooth ride. Intelligent electronics provide warnings when the batteries need charging. Lights operate in such a way that all bulbs function even if one light blows. An electronic speed reduction feature automatically slows you down when going round bends, ensuring safe manoeuvres at all times. Top quality components mean the **Orion** is reliable and easy to service. Batteries and electronics can be accessed easily and all parts are protected from corrosion, protecting the scooter against dirt and water and preventing any unnecessary malfunction.

## COMFORT IN MIND

Ergonomic steering on the **Orion** is specifically designed to offer relief to hands and arms during longer trips. User-friendly adjustment of the tiller is positioned opposite the handbrake and within easy reach. Together with the soft, adjustable swivel seat, this makes for an effortless, enjoyable drive. Combining ergonomics with style, the **Orion** is the natural choice for a relaxing, reliable drive.





# Rugged, Robust, Reliable: The Invacare® Comet

ENJOY A SAFE DRIVE



If you're looking for a safe, fast and enjoyable drive, the Invacare® **Comet** will give you a star performance! Whether you're stocking up at the supermarket or putting its outdoor performance to the test, you can rely on the **Comet** to get you where you want to go, quickly and in comfort.



Yes, you can.®





## FOR A SAFE, CONTROLLED DRIVE

With large wheels, a powerful motor and adjustable suspension, the **Comet** simply invites adventure. Easily tackling uneven terrain and obstacles, the **Comet** gives you the peace of mind and ability to go anywhere and everywhere, whilst high speeds and a large battery means you can go further, faster. Automatic speed reduction round bends ensures safe manoeuvres at all times. Intelligent electronics inform you when the battery needs recharging and indicators switch on and off automatically. Dual bulb headlights and a brake light all add up to extra safety.

## RELAXING AND RELIABLE

The **Comet** has been built with comfort in mind. Ergonomic handles help you execute precise manoeuvres in tight spaces and offer relief to hands and arms during longer trips. A soft, swivelling seat offers an optimal driving position and allows you to move around while seated. Constructed from top-quality components, the **Comet** is not only reliable, but easy to service. Simply remove the top cover to access the battery, electronics, motor and cables. All parts are protected from corrosion and water to safeguard the vital workings of the scooter, preventing unnecessary malfunctions.







## THE INVACARE ORION AND COMET OFFER A WIDE RANGE OF SAFETY FEATURES




### **Hand brake**

Ensures immediate braking if required.



### **Automatic speed reduction**

This feature automatically reduces your speed around bends.

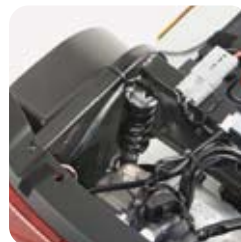

### **Full Suspension**

A fully-suspended chassis provides a comfortable drive at all times.



**Yes, you can.®**



### **Full lighting system with dual bulb headlights**

Dual headlights provide increased visibility. Indicators automatically go off after 20 seconds. Indicators can easily be programmed with or without an audible signal.



### **Low battery warning**

A visual warning on the battery level display and an audio warning (beep), tells you when the battery needs recharging.

### **Large wheels**

Provides greater feeling of safety when negotiating kerbs and obstacles.



### **Integrated brake light**

Heavy duty brake lighting warns people around you that the scooter is slowing down.

### **Audible reverse alert**

An audible signal (beep) lets people know that the scooter is reversing.

### **Anti-tip wheels**

Prevents the scooter from tipping backwards when climbing steep slopes or negotiating high obstacles.

### **Anti-corrosion chassis**

Anti-corrosion-treated chassis prolongs the life of the scooter.



### **Built-in anti-splash guards**

Protects electronics and transaxle from dirt and water.

### **Intelligent lighting system protection**

Only the bulb that has blown goes out: all others remain lit (scooter is still visible).

### **Dependable electrics**

Best-in-class quality wiring and electrical connections optimise reliability.

### **Guaranteed reliability**

Maintenance and easy serviceability are key design features. Batteries and electronics are easily accessible.



### **Two-step disengaging lever**

Prevents the scooter from accidentally free-wheeling.

### **Finite tiller movement**

The tiller can be easily adjusted with a lever to suit your needs.

### **Control panel**

A clear, intuitive keypad with domed buttons ensures accurate operation. LED and audio signals indicate when each function is activated.



### **Ergonomic steering**

Prevents fatigue in arms and hands during long trips.



# THE INVACARE ORION & COMET ACCESSORIES

Only available on Comet



## Seat raiser

An electronic mechanism raises the seat to your referred height.



## Crutch & Cane holder

Keeps your crutch or cane stable when on the move.



## Lockable front box

Secure box for personal valuables.



## Seat suspension module

Adds extra comfort during driving.



## Cup holder

Stores water or a soft drink.



## Heavy-duty storage cover

Keeps the scooter clean and dry.




## Off-board charging kit

Simply remove batteries to charge with the off-board charger kit.



### Mirror

Moveable mirror can be mounted to the left or right.

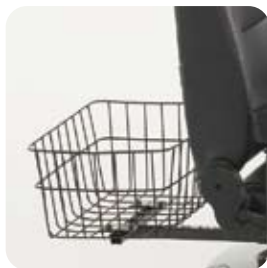

### Rear basket

Perfect for bulky or extra shopping.



### Oxygen bottle holder

Safe storage of oxygen during transit.



### Seat belt (Orion)

Ensures you are kept securely and comfortably in place.



### Safety vest



### Seat bag



### Walking frame holder

Easy storage for your walking frame.



### Hour counter



# Technical datas, Safety features and Accessories

This document entails all the comparative technical elements that you will need to select the scooter that will best suit your needs.

## ACCESSORIES

	LYNX	LEO	ORION	COMET
<b>Cane and crutch holder</b>	•	•	•	•
<b>Walking frame holder</b>		•	•	•
<b>Cup holder</b>			•	•
<b>Oxygen bottle holder</b>		•	•	•
<b>Lockable front box</b>		•	•	•
<b>Rear Basket</b>	•	•	•	•
<b>Seat bag</b>			•	•
<b>Charging off board kit</b>		•	•	•
<b>Mirror</b>		•	•	•
<b>Seat belt</b>		•	•	•
<b>Safety vest</b>			•	•
<b>Seat cover</b>			•	•
<b>Storage cover</b>	•	•	•	•
<b>Hour counter</b>			•	•
<b>Seat raiser</b>				•
<b>Seat suspension</b>			•	•
<b>Colour of shrouds</b>	 Moonstone Silver  Ruby Red  Sapphire Blue	 Moonstone Silver  Sapphire Blue	 Moonstone Silver  Onyx Blue	 Moonstone Silver  Onyx Blue
<b>Upholstery colour</b>	 Black vinyl	 Grey/Black vinyl	 Black vinyl	 Black vinyl

ASPHITE 334 860 265 R.C.S. TOURS - Photos : Paul Prunier - © Invacare® International Srl 2008. All rights reserved. Lym, Leo, Orion & Comet are registered trade marks of Invacare.

Invacare Ltd, South Road, Bridgend Industrial Estate,  
Bridgend, CF31 3PY - United Kingdom - Tel : (44) 1656 664 321  
[www.invacare.co.uk](http://www.invacare.co.uk) - [uk@invacare.com](mailto:uk@invacare.com)



All information quoted is believed to be correct at time of print. Invacare® reserves the right to alter product specifications without prior consultation.  
Please note colours may vary slightly from those displayed above.

Yes, you can.®

## TECHNICAL DATAS

• LYNX • LEO • ORION 3 wheels/\* 4 wheels • COMET 3 wheels/\* 4 wheels



- 43 cm
- 47 cm
- 51 cm
- 51 cm




- 38 cm
- 41 cm
- 47 cm
- 49 cm



- 37 - 42 cm
- 46 - 56 cm
- 42 - 49,5 cm
- 44 - 51,5 cm



- 40 cm
- 47,5 cm
- 51 cm
- 55 cm

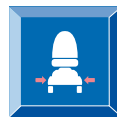

- 50 cm
- 59 cm
- 62,5 cm
- 66 cm



- 101 cm
- 122 cm
- 127/\*129 cm
- 145/\*147,5 cm




- 39 kg (excl. batteries)
- 49 kg (incl. batteries)
- 59 kg (excl. batteries)
- 83 kg (incl. batteries)
- 71/\*78 kg (excl. batteries)
- 103/\*110 kg (incl. batteries)
- 81/\*84 kg (excl. batteries)
- 133/\*136 kg (incl. batteries)



- 135 kg
- 136 kg
- 136 kg
- 160 kg



- 2 x 12 Ah
- 2 x 36 Ah
- 2 x 50 Ah
- 2 x 75 Ah

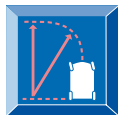

- 127 cm
- 131 cm
- 112,5/\*140 cm
- 156,5/\*163 cm



- 5 cm
- 6 cm
- 8 cm
- 10 cm



- 4 cm
- 7 cm
- 9 cm
- 13,5 cm



- 8°
- 10°
- 8°/\*12°
- 10°



- 150 W
- 240 W
- 6 km/h : 240 W - 500 W<sup>(1)</sup>
- 10 km/h : 240 W - 600 W<sup>(1)</sup>
- 12 km/h : 250 W - 1200 W<sup>(1)</sup>
- 10 km/h : 250 W - 1200 W<sup>(1)</sup>
- 12,8 km/h : 250 W - 1200 W<sup>(1)</sup>
- 15 km/h : 250 W - 1200 W<sup>(1)</sup>



- 16 km
- 36 km
- 42 km
- 60<sup>(2)</sup> - 55<sup>(3)</sup> km



- 8 km/h
- 8 km/h
- 6-10-12 km/h
- 10-12,8-15 km/h

<sup>(1)</sup> Peak performance - <sup>(2)</sup> 10 km/h - <sup>(3)</sup> 12,8 km/h / 15 km/h

## SAFETY FEATURES

	LYNX	LEO	ORION	COMET
<b>Control panel:</b> Clear and intuitive keypad	•	•	•	•
<b>Control panel:</b> Domed buttons, LED and audio signals		•	•	•
<b>Low battery visual and audio warning</b>		•	•	•
<b>Handbrake</b>			•	•
<b>Automatic speed reduction around bends</b>			•	•
<b>Full suspension: front and rear</b>			•	•
<b>Large wheels</b>		•	•	•
<b>2-step disengaging lever</b>		•	•	•
<b>Full lighting system – DUAL bulb headlight</b>		•*	•	•
<b>Full lighting system – ‘Xenon’- style lights</b>		•		
<b>Integrated brake light</b>		•	•	•
<b>Reverse audible warning alert</b>	•	•	•	•
<b>Anti-tip wheels</b>	•	•	•	•
<b>Anti-corrosion chassis</b>	•	•	•	•
<b>Intelligent lighting system protection</b>		•	•	•
<b>Splash protection</b>		•	•	•
<b>Dependable electrics</b>	•	•	•	•
<b>Reliability guaranteed</b>	•	•	•	•

\* only available in Germany