



INVACARE  
**Softform**

**MATTRESSES AND ACCESSORIES**

# **Softform**

## **THE INVACARE SOFTFORM RANGE**

The Invacare Softform range offers a complete solution to effective pressure reduction management. Designed in consultation with clinical professionals, these products are constructed from high specification materials to meet the demands of modern healthcare environments.

**SECTION A PRESSURE-REDUCING MATTRESSES**

**SECTION B BIRTHING CARE AND SPECIAL CARE BABY RANGE**

**SECTION C ACCESSORIES**

**SECTION D TECHNICAL SPECIFICATIONS AND WARRANTIES**

# PRESSURE-REDUCING MATTRESSES

comfortable





# Softform Premier Glide

RISK	Very High Risk
WATERLOW SCALE	25*
FOAM WARRANTY	8 Years
COVER WARRANTY	4 Years
REPLACEMENTS	All Parts
WEIGHT LIMIT	39 stone/247.6 kg
DIMENSIONS	197cm x 88cm x 15.2cm
PRODUCT CODE	SPGM197

**Introducing The Invacare Softform Premier Glide mattress, designed to conform with electric profiling beds.**

#### THE HIGH RISK MATTRESS FOR THE MODERN WARD

The Invacare Softform Premier Glide is an exciting new development from Invacare - a mattress that provides exceptional pressure reduction and patient comfort, and which uniquely conforms to the profile of an electric ward bed. The patented design of this mattress allows the top surface supporting the patient to remain in position when profiling. This significantly reduces inappropriate and unwanted patient movement. Such inappropriate movement may contribute to tissue damage resulting from shear and friction forces, particularly on vulnerable areas. Furthermore, a reduction in inappropriate patient movement delivers a benefit to Moving and Handling which would positively contribute to both pressure ulcer prevention and nursing resource management.





## FEATURES

## BENEFITS

Patented Glide system allows mattress to fully conform to the bedframe

Optimises the properties of a profiling bed

Top surface of the Glide remains under the patient throughout the profiling of an electric bed

Significantly reduces shear and friction forces, most particularly on the heel areas, which are especially vulnerable to pressure ulcer development

Cover fabric conforms to deep contoured foam castellations on upper surface of the insert pad when mattress is profiled

Minimises “rucking” of cover fabric on top surface, reducing potential for tissue damage

Suitable for Very High Risk patients (when used with a profiling bed)

Can significantly contribute to pressure ulcer prevention

Does not require turning or rotating

Significantly reduces manual handling issues, and helps to facilitate staff management resources

Two-way stretch PU waterproof cover

Provides excellent recovery and durability

High frequency welded seams

Reduces risk of fluid ingress, and positively contributes to infection control policy

Full length flap conceals the zip

Reduces possible fluid ingress through the zip

Toughened PU coated base

Prolongs longevity of the mattress, and reduces the risk of strikethrough, hydrolysis or delamination

Deep contoured insert pad for enhanced pressure area care

Maximises body contact area and optimises pressure reduction

Single piece foam “U” core

Facilitates patient transfer by providing effective side-wall support

Does not impede CPR procedures

Clinically assessed to ensure this does not occur

Double-coated PU sides

Reduces potential for physical damage

Top cover is fully detachable


Aids inspection and allows ease of removal for laundering in cases of heavy contamination. Also facilitates cost benefit of allowing the replacement of individual parts, rather than the whole mattress


## NO NEED TO TURN OVER

\*Essential nursing care is pivotal in pressure ulcer prevention. This mattress will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. The Waterlow score is a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.

**INVACARE**  
**Softform**  
*Superior*

<b>RISK</b>	Very High Risk
<b>WATERLOW SCALE</b>	25*
<b>FOAM WARRANTY COVER WARRANTY</b>	8 Years 4 Years
<b>REPLACEMENTS</b>	All Parts
<b>WEIGHT LIMIT</b>	39 stone/247.6kg
<b>DIMENSIONS</b>	197cm x 88cm x 15.2cm
<b>PRODUCT CODE</b>	SSM197/GB

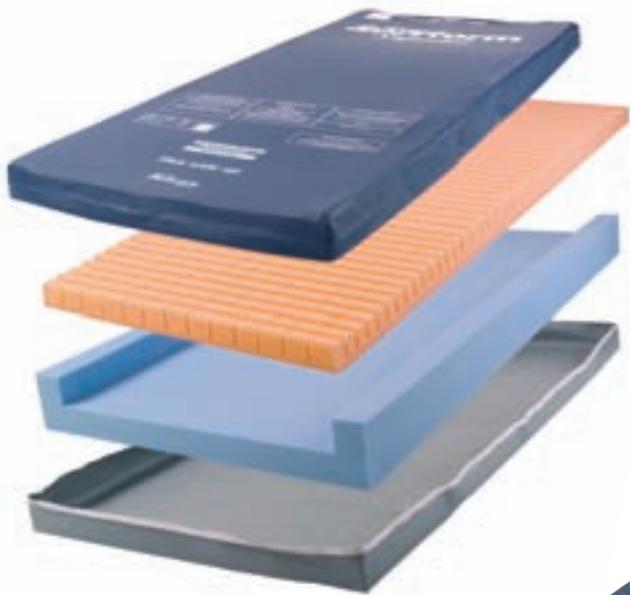

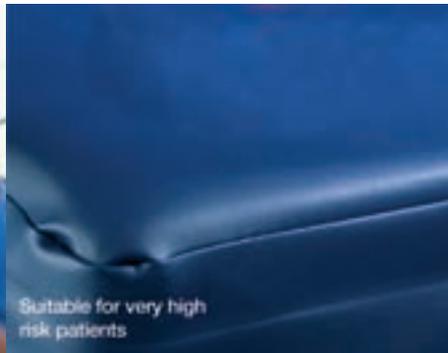

Handles offer improved transportation

Suitable for very high risk patients

## FEATURES

## BENEFITS

Two-way stretch cover

Aids pressure reduction

Optional handles

Offers improved transportation of the mattress (handles are fitted to one side only)

High frequency welded zip, centrally located, with concealing flap

Reduces risk of fluid ingress

Highly visible, easy to follow instructions

Offers best practice advice on moving and cleaning

Tapering in heel zones

To improve pressure reduction in the high-risk heel zone

Tapered top surfaces

Enables easy use of bedframe backrests

Deep contoured insert pad for enhanced pressure area care

Maximises body contact area and optimises pressure reduction

Single piece foam "U" core

Facilitates patient transfer by providing effective side-wall support

Does not impede CPR procedures

Clinically assessed to ensure this does not occur

Modular construction

Enables easy parts replacement, increasing product lifetime

Double-coated PU sides

Reduces risk of strikethrough, hydrolysis or delamination

Does not require rotation

Reduces moving and handling risks for nursing staff

Toughened PU-coated base


Promotes longevity of the mattress


## NO NEED TO TURN OVER

\*Essential nursing care is pivotal in pressure ulcer prevention. This mattress will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. The Waterlow score is a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.


**INVACARE**  
**Softform**  
*Premier*

<b>RISK</b>	High / Very High Risk
<b>WATERLOW SCALE</b>	20*
<b>FOAM WARRANTY</b>	8 Years
<b>COVER WARRANTY</b>	4 Years
<b>REPLACEMENTS</b>	All Parts
<b>WEIGHT LIMIT</b>	39 stone / 247.6 kg
<b>DIMENSIONS</b>	197cm x 88 cm x 15.2 cm
<b>PRODUCT CODE</b>	SPM197

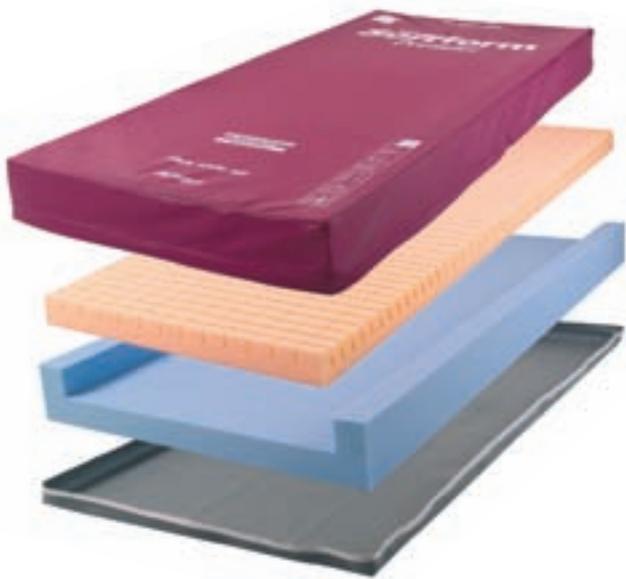

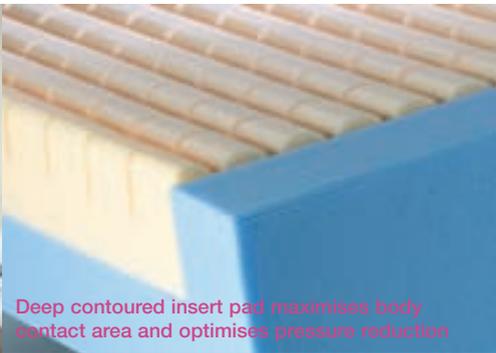

Deep contoured insert pad maximises body contact area and optimises pressure reduction

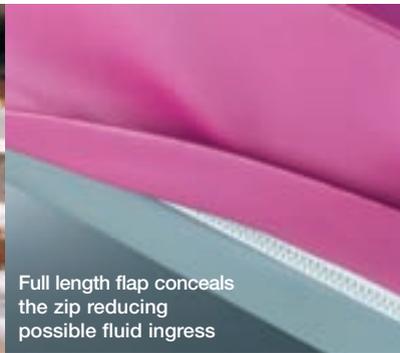

Full length flap conceals the zip reducing possible fluid ingress

## FEATURES

## BENEFITS

Suitable for high risk patients (very high risk when used with a profiling bed)

Offers high levels of protection for patients

Does not require turning

Reduces manual handling problems

Two-way stretch PU waterproof cover

Provides excellent recovery and durability, and helps reduce shear and friction forces

High frequency welded seams and zip

Reduces risk of fluid ingress

Full length flap conceals the zip

Reduces possible fluid ingress through the zip

Toughened PU coated base

Prolongs longevity of the mattress

Deep contoured insert pad for enhanced pressure area care

Maximises body contact area and optimises pressure reduction

Single piece foam "U" core

Facilitates patient transfer by providing effective side-wall support

Does not impede CPR procedures

Clinically assessed to ensure this does not occur

Double-coated PU sides

Reduces risk of strikethrough, hydrolysis or delamination

Top cover is fully detachable

Aids inspection and allows ease of removal for laundering in cases of heavy contamination

Modular construction

All parts are replaceable

## NO NEED TO TURN OVER

\* Essential nursing care is pivotal in pressure ulcer prevention. This mattress will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. The Waterlow score is a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.

**INVACARE**  
**Softform**  
*Premier Visco*

<b>RISK</b>	High Risk
<b>WATERLOW SCALE</b>	20*
<b>FOAM WARRANTY</b>	4 Years
<b>COVER WARRANTY</b>	4 Years
<b>REPLACEMENTS</b>	All Parts
<b>WEIGHT LIMIT</b>	39 stone / 247 kg
<b>DIMENSIONS</b>	197cm x 88cm x 15.2cm
<b>PRODUCT CODE</b>	SPM197 / VIS

**THE HIGH RISK VISCO MATTRESS**

The Softform Premier Visco mattress features a deep contoured visco insert pad, which is high-density and temperature-sensitive to the body, allowing it to conform to the user's shape. Many nursing staff consider the properties of visco elastic foam to be highly beneficial for patients at "High Risk" of developing pressure ulcers. The Softform Premier Visco mattress effectively combines the successful design of the Premier Mattress with the highest quality visco elastic foam to deliver a support surface with exceptional levels of comfort.





## FEATURES

## BENEFITS

Visco elastic foam insert

This high-density and temperature-sensitive foam warms to conform to the user's shape, and maximises the body contact area. This affords exceptional comfort and facilitates pressure reduction by minimising peak pressures

Suitable for High Risk patients

Can significantly contribute to pressure ulcer prevention

Does not require turning

Requires horizontal rotate method only, which reduces manual handling problems

Two-way stretch PU waterproof cover

Provides excellent recovery and durability

High frequency welded seams

Reduces risk of fluid ingress, and positively contributes to infection control policy

Full length flap conceals the zip

Reduces possible fluid ingress through the zip

Toughened PU coated base.

Prolongs longevity of the mattress, and reduces the risk of strikethrough, hydrolysis or delamination.

Insert is interchangeable with Softform Premier's deep contoured insert pad

Insert can be changed to provide visco elastic or regular foam, depending on the individual patient's needs

Single piece foam "U" core

Facilitates patient transfer by providing effective side-wall support

Does not impede CPR procedures

Clinically assessed to ensure this does not occur

Double-coated PU sides

Reduces potential for physical damage

Top cover is fully detachable

Aids inspection and allows ease of removal for laundering in cases of heavy contamination. Also facilitates cost benefit of allowing the replacement of individual parts, rather than the whole mattress

## NO NEED TO TURN OVER

\* Essential nursing care is pivotal in pressure ulcer prevention. This mattress will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. The Waterlow score is a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.


**INVACARE**  
**Softform**  
*Original*

<b>RISK</b>	High Risk
<b>WATERLOW SCALE</b>	20*
<b>FOAM WARRANTY</b>	4 Years
<b>COVER WARRANTY</b>	4 Years
<b>REPLACEMENTS</b>	All Parts
<b>WEIGHT LIMIT</b>	26 Stone/165kg
<b>DIMENSIONS</b>	195cm x 88cm x 17.5cm
<b>PRODUCT CODE</b>	SF100

**THE PROVEN MATTRESS REPLACEMENT**

The highly successful Invacare Softform Original mattress has established itself as a leading mattress of choice within the UK healthcare marketplace. The Invacare Softform Original combines a unique design concept with high quality materials to offer protection to those at “High Risk” of pressure damage.



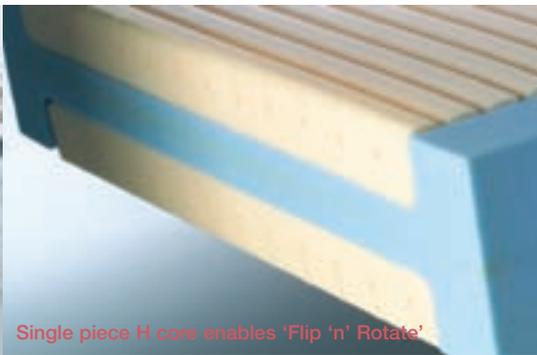

Single piece H core enables 'Flip 'n' Rotate'

## FEATURES

Two way stretch PU waterproof cover

Single piece H core

Two contoured insert pads for enhanced pressure area care

Modular construction

## BENEFITS

Provides excellent recovery and durability and helps reduce shear and friction forces

Enables "Flip'n'Rotate", effectively prolonging mattress life. Also facilitates patient transfer by providing side-wall support

Maximises body contact area and optimises pressure reduction

All parts are replaceable

## FOUR-WAY TURNING

\*Essential nursing care is pivotal in pressure ulcer prevention. This mattress will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. The Waterlow score is a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value

**INVACARE**  
**Softform**  
*Excel*

<b>RISK</b>	At Risk
<b>WATERLOW SCALE</b>	15*
<b>FOAM WARRANTY</b>	2 Years
<b>COVER WARRANTY</b>	2 Years
<b>REPLACEMENTS</b>	All Parts
<b>WEIGHT LIMIT</b>	22 Stone/140kg
<b>DIMENSIONS</b>	195cm x 88cm x 15.2cm
<b>PRODUCT CODE</b>	SXL195


**PRESSURE REDUCING MATTRESS**

The Invacare Softform Excel mattress provides a high level of patient comfort. The mattress core comprises a contour-cut high specification foam which provides even weight distribution, greater patient comfort and minimises the risk of bottoming out. The cover is made of a high quality two-way stretch vapour permeable PU material. Designed as a non-turning mattress, it has a tough PU coated base to promote longevity.





Toughened PU coated base able to withstand the challenges of modern wards

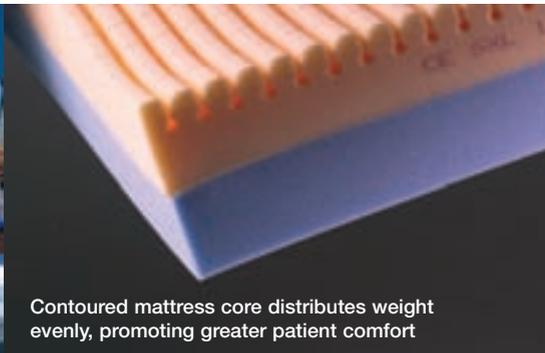

Contoured mattress core distributes weight evenly, promoting greater patient comfort



## FEATURES

## BENEFITS

Two way stretch PU cover

Provides excellent recovery and durability and helps reduce shear and friction forces

Contoured mattress core

Distributes patient weight evenly; promotes greater patient comfort

Base foam layer

Provides substantial support

Partially concealed zip around three sides

Allows easy inspection of core

Toughened PU coated base

Able to withstand challenges of modern wards

## NO NEED TO TURN OVER

\*Essential nursing care is pivotal in pressure ulcer prevention. This mattress will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. The Waterlow score is a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value



# Softform

## A&E Trolley Mattress

<b>RISK</b>	At Risk
<b>WATERLOW SCALE</b>	15*
<b>WARRANTY</b>	1 Year
<b>REPLACEMENTS</b>	All Parts
<b>WEIGHT LIMIT</b>	TM20 & TM40: 17 Stone/108 kg TM60: 22 Stone /140 kg
<b>DIMENSIONS</b>	182.9 x 55.9 x 5/10/11.5cm
<b>PRODUCT CODES</b>	TM20/TM40/TM60/TM60VEL



### PRESSURE REDUCING TROLLEY MATTRESS

The Invacare Softform A&E Trolley mattress has been designed specifically to meet the demands of the modern A&E Department whilst providing comfort for the patient. The ergonomically designed, profiled surface of the mattress offers significant reduction of pressure and shear and friction forces. These features combine to effectively reduce the potential incidence of pressure ulcers in the A&E department.

## FEATURES

Anti-static 2 way stretch vapour permeable cover

Profiled top surface interlocked into a higher density supportive foam base core

High frequency welded flap-concealing zip

Two inspection zips, with 5cm anti-contamination cover flap

Options of either securing straps at head and foot ends, or velcro fastening attachments

Visco-elastic temperature sensitive version also available

## BENEFITS

Provides excellent recovery and durability

Reduces the possibility of bottoming out when a patient is in the semi-recumbent position

Reduces risk of fluid ingress

Allows easy inspection at audit

Prevents mattresses from slipping whilst on the trolley

Visco responds to variations in body heat, ensuring a more even redistribution of pressure over the total body surface

\*Essential nursing care is pivotal in pressure ulcer prevention. This mattress will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. The Waterlow score is a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.

**BIRTHING CARE  
AND SPECIAL CARE  
BABY RANGE**

promotive

**INVACARE**  
**Softform**  
*Birthing Mattress*

<b>RISK</b>	At Risk (with the Cap-on Cover)
<b>WATERLOW SCALE</b>	15* (with the Cap-on Cover)
<b>FOAM WARRANTY</b>	4 Years
<b>COVER WARRANTY</b>	1 Year
<b>REPLACEMENTS</b>	All Parts
<b>WEIGHT LIMIT</b>	39 Stone/247.6kg
<b>DIMENSIONS</b>	Customisation to size available
<b>PRODUCT CODE</b>	SBM

**THE CUSTOMISED BIRTHING MATTRESS SOLUTION**

Designed in consultation with some of the leading maternity units in the UK, the Invacare Softform Birthing Mattress has been designed to aid pressure reduction and to be highly durable. For added assurance in infection control, all seams are manufactured with an extremely strong, 10mm wide, high-frequency weld and audit inspection zips.





## FEATURES

## BENEFITS

Fully customised sizing

Dimensions of this product can be fully customised to fit most birthing beds

High frequency welded seams and zip

Designed to greatly reduce risk of fluid ingress

Zip concealed within a 3cm deep high-frequency welded flap

Reduces risk of fluid ingress and adds extra protection from zip leakage

Two-way stretch high quality waterproof vapour permeable PU cover

Provides excellent recovery and durability

Interior and exterior PU with toughened grey base

Delivers enhanced durability and prevents mechanical damage

Optional heavy duty base board

Delivers added structural support

High quality RX39/200 foam core

Firm foam interior suitable for purpose, while delivering improved comfort for expectant mothers

## NO NEED TO TURN OVER

\*Essential nursing care is pivotal in pressure ulcer prevention. This mattress will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. The Waterlow score is a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.




# Softform

## New-born

<b>RISK</b>	Low Risk / Comfort
<b>CHEMICAL FREE</b>	No chemicals or minerals in cover or foam*
<b>WARRANTY</b>	12 Months
<b>FIRE REGULATION</b>	Cigarette and Match
<b>FOAM INSERT</b>	Flat foam
<b>DIMENSIONS</b>	Can be customised up to maximum dimensions of 50.8cm x 86.4cm
<b>CODE</b>	SNM

### CRIB / BASSINET MATTRESS

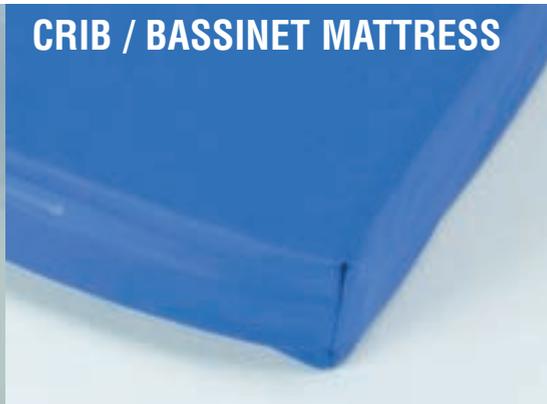

### FEATURES

High frequency welded seams and zip

Zip concealed within a 3cm deep high-frequency welded flap

Two-way stretch high quality waterproof vapour permeable PU cover

### BENEFITS

Designed to greatly reduce risk of fluid ingress

Reduces risk of fluid ingress and adds extra protection from zip damage

Provides excellent recovery and durability and helps reduce shear and friction forces

\*No Chemicals denotes that no additional chemicals or minerals have been added during the manufacturing of these products. Invacare has recently introduced the manufacture of mattresses without additional chemicals. This action has been taken in light of the studies showing a possible correlation between toxic gas positioning and sudden infant death (cot death). A clinical information pack is available on request, citing arguments for and against toxic gas poisoning.

**INVACARE**  
**Softform**  
*Incubator Pad*

<b>RISK</b>	High Risk
<b>CHEMICAL FREE</b>	No Chemicals or minerals in cover or foam*
<b>WARRANTY</b>	12 Months
<b>FIRE REGULATION</b>	Cigarette and Match
<b>FOAM FORMAT</b>	Profiled foam
<b>DIMENSIONS</b>	Can be customised up to maximum dimensions of 122cm x 71cm
<b>CODE</b>	SIP



**INVACARE**  
**Softform**  
*Cot Mattress*

<b>RISK</b>	Low Risk / Comfort (High Risk option also available)
<b>CHEMICAL FREE</b>	No Chemicals or minerals in cover or foam*
<b>WARRANTY</b>	12 Months
<b>FIRE REGULATION</b>	Cigarette and Match
<b>DIMENSIONS</b>	Can be customised up to maximum dimensions of 89cm x 140cm
<b>CODE</b>	SCM



## FEATURES

High frequency welded seams and zip

Zip concealed within a 3cm deep  
high-frequency welded flap

Two-way stretch high quality waterproof  
vapour permeable PU cover

## BENEFITS

Designed to greatly reduce risk of fluid ingress


Reduces risk of fluid ingress and adds extra protection from  
zip damage

Provides excellent recovery and durability and helps reduce shear  
and friction forces

\*No Chemicals denotes that no additional chemicals or minerals have been added during the manufacturing of these products. Invacare has recently introduced the manufacture of mattresses without additional chemicals. This action has been taken in light of the studies showing a possible correlation between toxic gas positioning and sudden infant death (cot death). A clinical information pack is available on request, citing arguments for and against toxic gas poisoning.

# ACCESSORIES

DELIVERING IMPROVED LEVELS OF PROTECTION FOR VULNERABLE AREAS

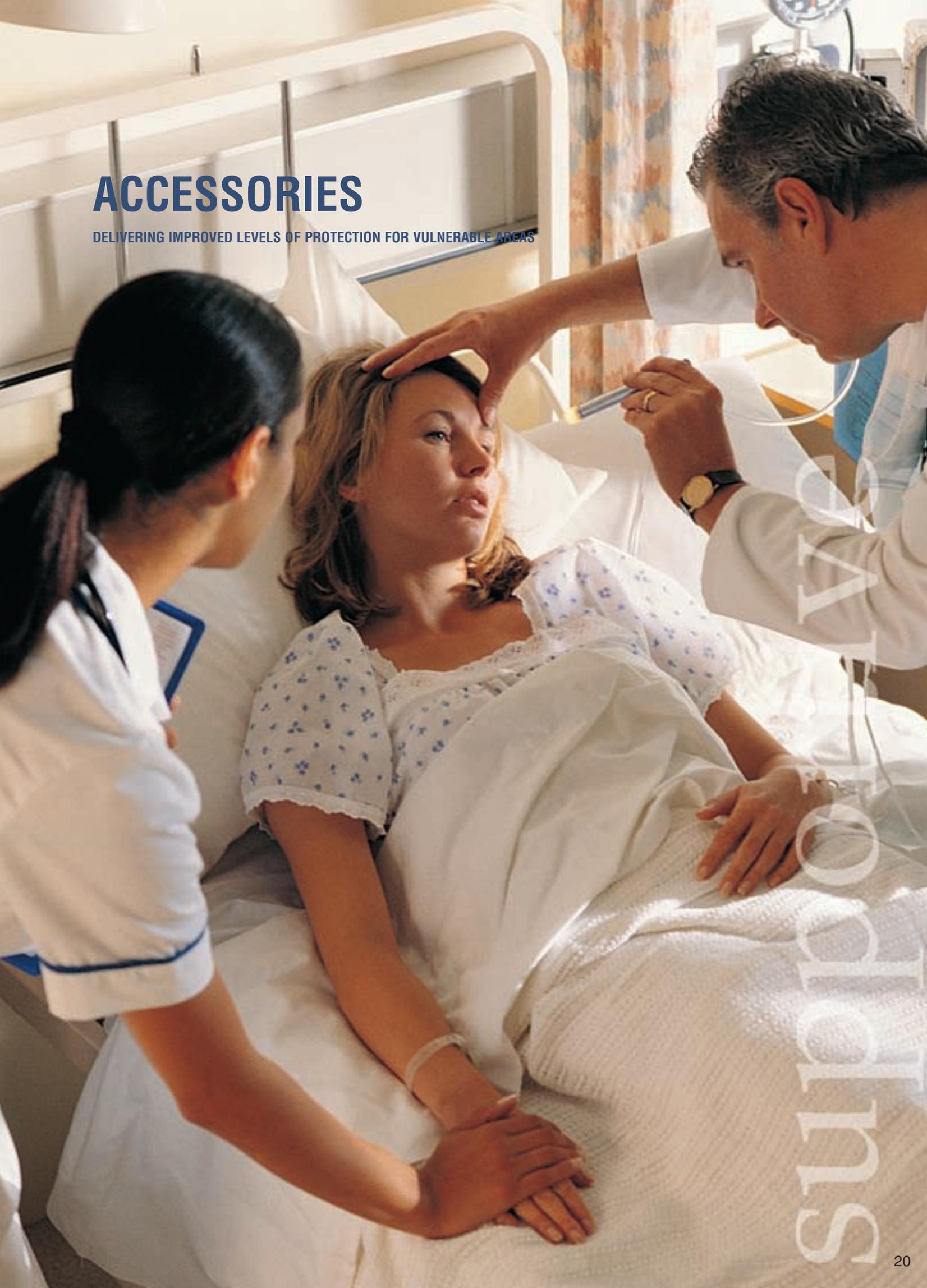

W  
V  
S  
U  
P  
P  
L  
I  
E  
S

**INVACARE**  
**Softform**  
*Heelpad*

The Invacare Softform Heelpad is designed to reduce pressure on the vulnerable heel area. The controlled volume of gel in the sacs allows pressure displacement and even weight distribution. Pressure on a patient's heels are kept to a minimum when the heels are in contact with the gel sacs.

<b>WARRANTY</b>	12 Months
<b>DIMENSIONS</b>	62cm x 48cm x 3.5cm (HP100) 32cm x 48cm x 3.5cm (HP50)
<b>PRODUCT CODES</b> Double Heel Pad Single Heel Pad	HP100 HP50
<b>REPLACEMENT PLATILON COVERS</b> (Set of five) Double Heel Pad Single Heel Pad	HP100 / C5 HP50 / C5



# **Softform** *Odstock Wedge*

<b>WARRANTY</b>	12 Months
<b>DIMENSIONS</b>	86.5cm x 57cm x 11.5cm
<b>PRODUCT CODE</b>	OSW100

The Odstock Wedge reduces pressure under the heels without increasing pressure on the thighs and calves. It improves and maintains posture in the semi recumbent position; the support under the thigh reduces the tendency for patients to slip down the bed, thus reducing shearing forces on the sacrum. By tilting the pelvis, it also reduces lumbar lordosis – desirable following back surgery and ideal for improving the sleeping posture of people with back problems. Pressure on the heels can be more than halved in comparison with a standard Kings Fund mattress.



FEATURES	BENEFITS
Twin gel sacs on Invacare Softform Heelpad (Single heel version also available)	Provide excellent immersion for high-risk heel areas
High frequency welded seams protects products from contamination	Ensures increased product integrity
Vapour permeable stretch cover with PU base	High quality materials promote product longevity
Invacare Softform Heelpad weighs 3.58kg Invacare Odstock Wedge weighs 1.5kg	Comparatively lightweight - ensures safe handling; easily transferred between Acute and Community use
Products designed to be used in conjunction to provide improved support to lower limbs	Reduces requirement for patient handling
Heelpads supplied with three platilon covers	Platilon covers are disposable, facilitating improved infection control



# **Softform** *Flexipad*

## THE INVACARE FLEXIPAD RANGE


Invacare' Flexipads offer a range of solutions for facilitating a comfortable support surface. Flexipads are gel sacs, with PU covers. They can be used to support limbs, joints, the torso or head, and are available in a range of formats:

- **SINGLE FLEXIPAD** – fitted with a strap for fastening to chairs, the Flexipad increases the support surface area under elbows, arm or back.
- **DOUBLE FLEXIPAD** – gel sacs with a foam core – highly versatile in a number of applications including supporting large bony prominences.
- **FLEXIPAD ULTRA** – a lighter, more flexible format, suitable for a wide range of applications, available with or without straps.
- **FLEXIPAD MULTIPACK A** – Four Flexipad Ultras in one pack (No straps)
- **FLEXIPAD MULTIPACK B** – Two Flexipad Ultras (No straps), Two Flexipad Singles (with straps) and one Double in one pack. Additional adjustable extension straps are included.
- **FLEXIPAD MULTIPACK C** – Four Flexipad Doubles in a pack

The Flexipad range is available in either PU wipedown or Black Smooth formats.

<b>WARRANTY</b>	12 Months	
<b>REPLACEMENTS</b>	Outer cover and gel sac	
<b>DIMENSIONS</b>		
Single Flexipad	22 x 33 x 5cm	
Double Flexipad	22 x 33 x 5cm	
Flexipad Ultra	22 x 33 x 3.5cm	
<b>PRODUCT CODES</b>		
Softform Flexipad Single	FGP50	
Softform Flexipad Double	FGP100	
Softform Flexipad Ultra	FGP75	
Softform Flexipad Ultra (with Straps)	FGP75/S	
Softform Flexipad Multipack A	FGPA	
Softform Flexipad Multipack B	FGPB	
Softform Flexipad Multipack C	FGPC	



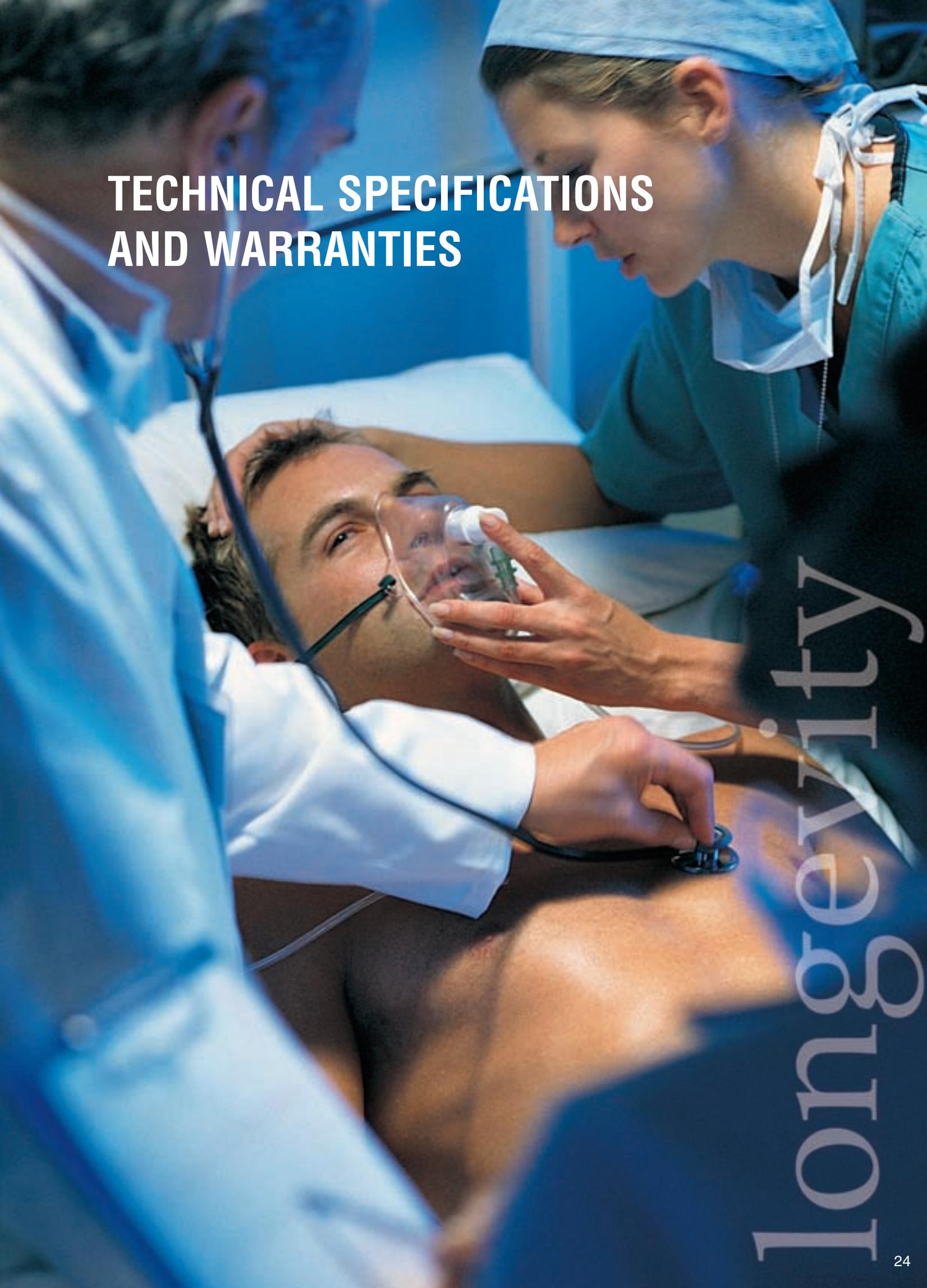

**TECHNICAL SPECIFICATIONS  
AND WARRANTIES**

Longevity

## PRODUCT

## WARRANTY

## FIRE TESTING

## GRADE REF & COLOUR

<b>GLIDE</b>	4 Years cover 8 Years foam	BS 7177: Crib 5	Gr. RX 39/120 Pink Gr. RX 39/200 Blue
<b>SUPERIOR</b>	4 Years cover 8 Years foam	BS 7177: Crib 5	Gr. RX 39/120 Pink Gr. RX 39/200 Blue
<b>PREMIER</b>	4 Years cover 8 Years foam	BS 7177: Crib 5	Gr. RX 39/120 Pink Gr. RX 39/200 Blue
<b>PREMIER VISCO</b>	4 Years*	BS 5852: Crib 5	Gr. RX 39/120 Pink Gr. RX 39/200 Blue Gr. VC5580 White
<b>ORIGINAL</b>	4 Years cover 8 Years foam	BS 7177: Crib 5	Gr. RX 39/120 Pink Gr. RX 39/200 Blue
<b>EXCEL</b>	2 Years	BS 7177: Crib 5	Gr. RX 39/120 Pink Gr. RX 39/200 Blue
<b>A&amp;E TROLLEY MATTRESS</b>	1 Year	BS 5852: Crib 5	Gr. RX 36/125 Beige Gr. RX 39/200 Blue
<b>BIRTHING MATTRESS</b>	1 Years cover 4 Years foam	BS 7177: Crib 5	Gr. RX 39/200 Blue
<b>NEWBORN</b>	1 Year	BS EN 597-1,2	Gr. RN 26/100 White
<b>INCUBATOR PAD</b>	1 Year	BS EN 597-1,2	Gr. RN 26/100 White
<b>COT MATTRESS</b>	1 Year	BS EN 597-1,2	Gr. RN 26/100 White
<b>HEEL PAD</b>	1 Year	BS 5852: Crib 5	Gr. RX 33/125 White
<b>ODSTOCK WEDGE</b>	1 Year	BS 5852: Crib 5	Gr. RX 33/125 White

## CLEANING OF INVACARE PRODUCTS

**GENERAL CLEANING** – to prevent delamination, use a solution of neutral detergent and hot water, or alternatively a 0.1% chloride solution. Apply with single use wipe. Following application, the cover should be rinsed with clean water (damp cloth) and **thoroughly** dried.

Removable contaminated covers should be laundered at 80°C with a dilute detergent solution. Please ensure covers are thoroughly rinsed and dried before placing back on the Invacare product.

(Invacare recommends that nursing staff should always actively liaise with a qualified Infection Control Nurse regarding decontamination protocols.)

\*Standard Configuration

NOMINAL DENSITY RANGE (Kg/M)	NOMINAL HARDNESS RANGE (N)	FATIGUE CLASS	WEIGHT OF PRODUCT (Kg)
38-40	105-135	Very Severe	13
38-40	180-200	Very Severe	
38-40	105-135	Very Severe	13
38-40	180-200	Very Severe	
38-40	105-135	Very Severe	13
38-40	180-200	Very Severe	
38-40	105-135	Very Severe	13
38-40	180-200	Very Severe	
55	80	N/A	
38-40	105-135	Very Severe	12.5
38-40	180-200	Very Severe	
35-37	105-135	Very Severe	10
38-40	180-200	Very Severe	
35-37	110-140	Very Severe	5.5
38-40	180-200	Very Severe	
38-40	180-200	Very Severe	15
25-27	85-110	Average	0.5
25-27	86-110	Average	0.5
25-27	87-110	Average	2.7
32-34	110-140	Severe	3
32-34	110-140	Severe	2.5



## CLINICAL BIBLIOGRAPHY RELATED TO SOFTFORM PRODUCTS

*A Static-Led Approach to Pressure Ulcers On Evaluation After Three Years* Pt No. (1504043)  
Jane James and Janet Thomas, Tissue Viability Nurses, West Wales  
General Hospital, Carmarthenshire

*Evaluating Pressure Reducing Mattresses and Electric Bed Frames* Pt No. (1492885)  
David Gray et al Tissue Viability Nurse, Aberdeen Royal Hospital

*A Study Of The Performance Of A Pressure Reducing Foam Mattress  
After Three Years Of Use* Pt No. (1492978)  
David Gray et al Tissue Viability Nurse, Aberdeen Royal Hospital

*A Randomised Clinical Trial of Two Types of Foam Mattresses* Pt No. (1492979)  
David Gray et al Tissue Viability Nurse, Aberdeen Royal Hospital

*Softform Birthing Mattress Product Focus* Pt No. (1492983)  
Verna Pritchard, Tissue Viability Nurse, Warwick Hospital

*Evaluation Of Static Support Surfaces* Pt No. (1492981)

*Understanding Pressure Ulcers* Pt No. (1492223)

Copies of the above clinical papers are available on-line. Please log on to the Invacare Limited website [www.invacare.co.uk](http://www.invacare.co.uk) and go to the Antidecubitus & Positioning section for mattresses. Alternatively you can order from Customer Services using Email: [ordersuk@invacare.com](mailto:ordersuk@invacare.com).





**Yes, you can.<sup>®</sup>**

Invacare Limited.

South Road

Bridgend Industrial Estate

Bridgend

CF31 3PY, UK

United Kingdom

**Tel:** +44 (0)1656 647327

**Fax:** +44 (0)1656 649016

**Email:** uk@invacare.com

**Orders:** ordersuk@invacare.com

**Web:** www.invacare.co.uk

Invacare Ireland Limited.

Unit 5 Seatown Campus

Seatown Road

SWORDS

County Dublin

Ireland


**Tel:** +(353) 1 810 7084

**Fax:** +(353) 1 810 7085

**Email:** ireland@invacare.com

**Web:** www.invacare.ie

Distributed by:



Sofform is a registered Trade Mark of Invacare Limited.

Invacare Limited reserve the right to change products without prior notification.

V7.2.08



1488715