

Invacare® | *Softform Premier Active*




Warranty
Foam 8 Years
Cover 4 Years
Pump 2 years
Insert 2 years

Static Dynamic pressure reducing mattress

The unique patented design of the *Invacare Softform Premier Active* features an alternating air insert beneath the castellated foam of the tried and trusted *Softform Premier* mattress. The mattress retains its properties as a 'Very High Risk'* static surface, but should a patient require stepping up to a dynamic surface, a discreet pump can be fitted to the air insert, transforming the mattress into an alternating surface delivering additional levels of pressure relief. There is no requirement to transfer the client, thereby minimising any discomfort.



Yes, you can.®

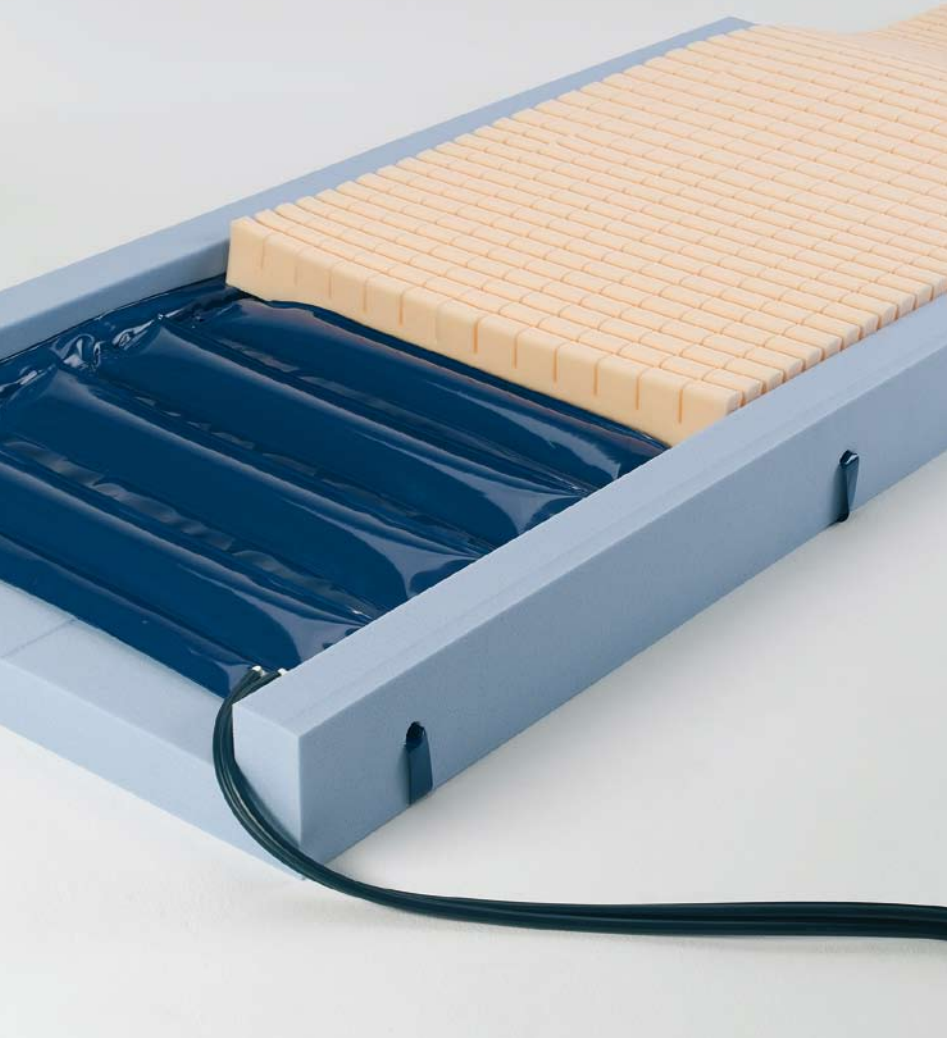

Static Dynamic mattress

The air insert contained within the mattress foam construction allows stepping up to and down between static and dynamic alternating systems with minimum patient disruption. The air cells alternate on a ten minute cycle beneath the foam insert, delivering a very low impact at the point of change, thereby delivering high levels of comfort in addition to improved pressure reduction.

Pump automatically assesses patient weight

Sensory software within the pump can gauge the patients weight, enabling pressure to be maintained accordingly at an appropriate level. The pump is also fully alarmed in case of power failure or fault.

Multi stretch cover

The mattress cover is multi stretch and conforms to the deep contoured castellations on the upper surface of the mattress, minimising "rucking" and further reducing the potential for tissue damage.

Features



Air insert within Softform Premier mattress foam construction

Allows stepping up and down between static and dynamic alternating systems with minimum patient disruption.



Inflation of air cells alternates on a 10 minute cycle

Delivers an extremely comfortable mattress surface, while reducing interface pressures on any one area of the body.



'Very High Risk' static mattress


In the event of a power failure, the patient is still resting on a 'Very High Risk' pressure reducing surface, minimising the potential for pressure damage occurring.



Features



The *Invacare **Softform Premier*** is a leading pressure-reducing static mattress.

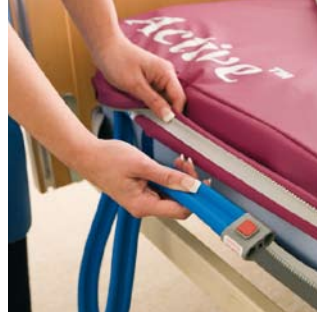

To step up a patient, retract the pipes from the mattress.



Fit pipes to the Active pump.



Simply switch on, no additional adjustments required.



The client is now on an alternating surface.



Technical data

Softform Premier Active						
	Ranging from 80 - 90 cm	Ranging from 181 x 210 cm	15.2 cm	13 kg	247.6 kg (39 stone)	Cleaning: 80°

Manufactured to European Standards: BS EN 597-1 Cig, BS EN 597-2 Match, BS7177:1996 Crib Source 5. Tested to pass ISO 1419:1995.

Colour



Fuchsia (cover)

INVACARE Limited
Pencoed Technology Park - Pencoed -
Bridgend CF35 5AQ
Tel: +44 1 656 776 222
Fax: +44 1656 776 220
E-mail: uk@invacare.com
www.invacare.co.uk
Sales Order E-mail:
ordersuk@invacare.com

INVACARE Ireland Ltd
Unit 5 - Seatown Business Campus
Seatown Road - SWORDS - County
Dublin
Tel. +353 1 8107084
Fax +353 1 8107085
Email: Ireland@invacare.com
www.invacare.ie

© 2011 Invacare International Sàrl
All rights reserved.

All information quoted is believed to be correct at time of print.
Invacare reserves the right to alter product specifications without prior consultation.

Softform Premier Active - UK
07/2011 - 1505223

* Essential nursing care is pivotal in pressure ulcer prevention. This mattress will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action based planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. A range of assessment scales can be used as a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.



Yes, you can.®