

Invacare® | **Comet™**



 REHACARE
DESIGN AWARD
2008/2009
GOOD DESIGN

Rugged, Robust, Reliable: The *Invacare Comet*

The **Comet** is packed with features that enable a safe, fast and enjoyable drive. Reliability remains the hallmark of this secure and speedy scooter. Whether stocking up at the supermarket or challenging its outdoor performance, you can rely on the **Comet** to get you where you want to go, quickly, in comfort and with peace of mind. A soft swivelling seat offers you an optimal driving position, and allows you to move around and enjoy activities while seated. By combining rugged outdoor capabilities with ergonomics and good looks, the *Invacare Comet* is the ultimate solution for your transportation needs.



Yes, you can.®



For a safe, controlled drive

With its large wheels, a powerful motor and adjustable suspension system in both the front and rear of the chassis, the **Comet** simply invites adventure. The **Comet** is equipped with an electronic speed reduction feature, which automatically decreases speed when going round bends: this ensures safe, comfortable manoeuvres at all times. Intelligent electronics offer a visual and audio signal when the batteries need recharging. The unique lighting system includes an easily visible DUAL-bulb headlight, indicators that automatically switch on and off and an integrated heavy-duty brake light for safety.

Compact for convenience

Ergonomic handles help you execute precise manoeuvres in tight spaces and offer relief to hands and arms during longer trips. The adjustable tiller lever is positioned opposite the handbrake and is within easy reach. A two-step disengaging lever protects the scooter from free-wheeling if accidentally knocked out of drive mode.

Simplicity and ease of use

Using only quality components not only makes the **Comet** reliable, but easy to service. Simply remove the top cover to access batteries, electronics, motor and cables. The motor and electronics are protected from corrosion and water (splash protection) which safeguard the vital workings of the scooter preventing

any unnecessary malfunction. The lights operate in such a way that lights still function if one bulb blows.

The **Comet** has undergone extensive testing in the German independent TÜV test laboratory and complies with the Medical Device Standards for Europe.

Features and Accessories



Easily adjust the angle of the tiller



Ergonomic steering
Prevents fatigue during long trips.



Hand brake
Ensures immediate braking if required.



Brake light
Heavy duty brake lighting warns people that the scooter is slowing down.

Features and Accessories



Two-step disengaging lever
Prevents the scooter from accidentally free-wheeling.



Built-in anti-splash guards
Protects electronics and transaxle from dirt and water.



Automatic speed reduction
Speed lowers during turns to ensure a safe manoeuvre.



Suspension
Can be adjusted at front and rear for optimal driving comfort.



Seat Bag
Fixes over the back of the seat to provide additional storage.



Lockable front box
Secure box for personal valuables.



Seat suspension module
For increased comfort on the move.



Cup holder
Stores water or soft drinks.



Heavy-duty storage cover
Keeps the scooter clean and dry.




Off-board charging
Simply remove batteries to charge with the off-board charger kit.

Technical data


					
Invacare Comet 510 mm	490 mm	440 - 515 mm	550 mm	660 mm	1475 mm
					
Invacare Comet 136 kg	84 kg ¹	160 kg	2 x 75 Ah	3260 mm	100 mm
					
Invacare Comet 135 mm	10°	55 km	12,8 km/h	250 W / 1200 W ²	

1/ Weight of heaviest part
2/ Peak performance

Colour of shrouds

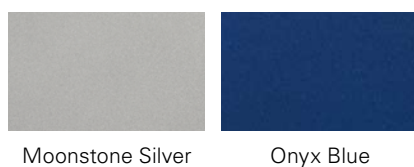

Upholstery

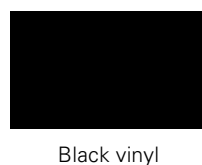

Invacare Limited
Pencoed Technology Park - Pencoed -
Bridgend CF35 5AQ - United Kingdom
Tel: +44 1 656 776 222
Fax: +44 1656 776 220
E-mail: uk@invacare.com
Sales Order E-mail:
ordersuk@invacare.com
www.invacare.co.uk

Invacare Ireland Ltd
Unit 5 - Seatown Business Campus
Seatown Road - SWORDS - County
Dublin - Ireland
Tel: +353 1 8107084
Fax +353 1 8107085
Email: Ireland@invacare.com
www.invacare.ie

© 2012 Invacare International Sàrl
All rights reserved.

All information quoted is believed to be
correct at time of print.
Invacare reserves the right to alter
product specifications without prior
consultation.

Comet - UK - 02/2012 - 1506451