



Freedom[®] NXT[™]

Generation NEXT Folding Tilt-in-Space Mobility



Yes, you can.[®]



Creative thinking
beyond tilt-in-
space...

A hybrid "by design,"
the NXT is strong and
durable— yet lightweight,
compact, and easy to use.



For kids as they
grow...

Three fasteners adjust seat
frame width 50mm (2") in
most sizes.

Two fasteners adjust seat
depth another 100mm (4") in
each size.

Two fasteners adjust back
angle from 80° to 120°.



For kids on the go,
we put our best
foot forward...

Foot-activated tilt
with flip-up pedal.

Foot-activated wheel lock.

Tipping lever / anti-tipper.



Transforms, transfers, and transports...

With the NXT, the "ins and outs" of travel have never been easier.



Transforms as easy as 1, 2, 3...

- 1: Quick release wheels and detachable legrests.
- 2: Flip back armrests and fold down back.
- 3: Pull up on strap to fold and stow for transport or storage.




Transfers as easy as A, B, C...

- A: Travels light, weighing just 10 kilos.
- B: Folds into a small package.
- C: Fits into small spaces.



Transports safely...

The NXT and Freedom Seating system is fully crash tested. Tie-downs for transportation are on the strongest part of frame. Tie-down points are highly visible and easily accessible.

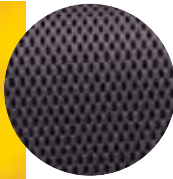

NEW Air Mesh upholstery

Allows for air to move freely between the seating surface and the body. Air Mesh is designed to help protect sensitive skin by cooling the body to reduce heat and moisture. Also available with an incontinence resistance interface.



Freedom Select positioning accessories

Anterior trunk supports, chest straps and pelvic positioning belts are available in virtually any style you prefer.



Freedom® NXT... Generation NEXT Folding Tilt-in-Space Mobility

Specifications

The NXT is versatile, quickly transforming for ever-changing kids



Seat Frame Widths

The NXT is available in 11 sizes from 255mm (10") wide to 510mm (20") wide.

Each size from 255mm (10") wide to 405mm (16") wide adjusts an additional 50mm (2") in width.

The 430mm (17") wide size adjusts an additional 25mm (1") in width to 460mm (18") wide.

The 460mm, 485mm & 510mm (18", 19" & 20") wide sizes are set.

Seat Frame Depths

The NXT is available in 3 sizes: short, medium or long.

Each size adjusts an additional 100mm (4") in depth.

Seat frame depth is measured without seating, from the front of the back post to the front end of the seat frame.

The short range is 305mm to 405mm (12" - 16").

The medium range is 355mm to 460mm (14" - 18").

The long range is 405mm to 510mm (16" - 20").

Seat Depth Extensions

For additional seat frame depths, extensions can be added to the back of the seat frame. There is a 50mm, 75mm or a 100mm (2", 3" or 4") extension to choose from. These are available through our i-Custom department.

Seat to Floor Height

The NXT can be configured in a variety of seat to floor heights.

Selections in front castors and rear wheel sizes determine a range of choices from 318mm low to 470mm high (12.5" to 18.5"), and most everything in between.

Rider Capacity

91kg (14.3 stone)
See NXT prescription form for further details.

Tilt Range

0-45 degrees.



Invacare Ltd
Unit 4 Pencoed Technology Park
Pencoed CF35 5AQ
T: 01656 776222 F: 01656 776220
e: uk@invacare.com
www.invacare.co.uk

Invacare Ireland
Unit 5 Seatown Business Campus
Seatown Road
Swords
County Dublin
T: +353 1 8107084 F: +353 1 8107085
e: ireland@invacare.com
www.invacare.ie



Crash Tested to ISO 7176 Part 19