




"The ideal match for  
my active lifestyle"

Sarah Williams

## The perfect fit for your lifestyle

The *Invacare Action<sup>5</sup>* has been designed for the more active individual, looking for an affordable chair with great drive performance and the ability to fold easily for handling, storage and transportation purposes. With its two new advanced innovations (patent pending) and lightweight construction, the new *Action<sup>5</sup>* gives individuals a chair they can enjoy and rely on. Key features include:

- **Unique folding mechanism**
- **Stability of a rigid chair**
- **Quick step-less adjustability**
- **Highly configurable**

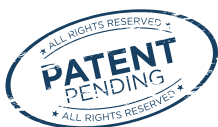

reddot award 2015  
winner



**Yes, you can.®**





### Unique folding mechanism

- The innovative design of the chair allows the frame to greatly reduce in size and the backrest to fold on top of the seat, making the chair compact.
- No additional handling is required to fold the chair, making it effortless for storage and transportation.




### All the Stability of a rigid chair

- The 'H-Shaped' horizontal folding mechanism (patent pending) provides an improved level of stability compared to other folding chairs on the market.
- The self-locking function ensures a secure seat position and makes the chair far easier to propel effectively improving drive performance.

### Adjust the chair whilst the individual remains seated

- Our unique, clever design allows backrest angle adjustments to be made whilst the person remains in the chair (patent pending).
- The step-less adjustability makes assessments and fittings both quick and easy to ensure the perfect fit for individuals.
- The built in scale allows for a more accurate backrest angle to be set for postural comfort.

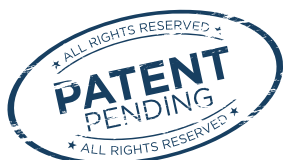

### Make it personal

- There are a variety of vibrant colours available to choose from, allowing individuals to design the chair around their personality.
- The **Action<sup>5</sup>** can also be configured depending on individual needs, starting from selecting a more active look with the fixed frame, or the swing away version for even easier folding and transportation.
- The **Action<sup>5</sup>** is cross compatible with other Invacare products including the Alber power add on range.



### Award winning

- Already an award winning chair, the Invacare **Action<sup>5</sup>** has won a Red Dot Award for Product Design.
- The Red Dot Award is not only the largest product competition in the world, but also one of the most prestigious design awards.



**reddot** award 2015  
winner





## CONFIGURABILITY



## LEGREST OPTIONS



## ACCESSORIES



### CONFIGURABILITY

#### A) Tension and depth adjustable seat

Adjustments can be made via easy to use straps.

#### B) Centre of gravity

100mm centre of gravity adjustments can be made in 25mm increments which can even be fine-tuned to 12mm.

#### C) Camber

Standard adjustments of 0-2° available with the option of 3 & 4° adjustments if required.

### LEGREST OPTIONS

#### D) Swing away legrests

Allows the legrests to swing either inwards or outwards.

#### E) Fixed front

An 80° fixed front frame is available for individuals who are more independent.

#### F) Selection of legrest options

Choose from 70, 80 or 90° legrests, elevating legrests, or an amputee support.

### ACCESSORIES

#### G) Backrest cover

Includes a backrest pocket and reflective strip for better visibility.

#### H) Active anti-tipper

Can be operated by both the individual and the carer for additional support.

#### I) Front wheels

A selection of wheels are available from 4" to 7", including the stylish aluminium soft roll option.

## ARMRESTS



A



B



C

## SWING AWAY AND REMOVABLE ARMREST




D



E



F

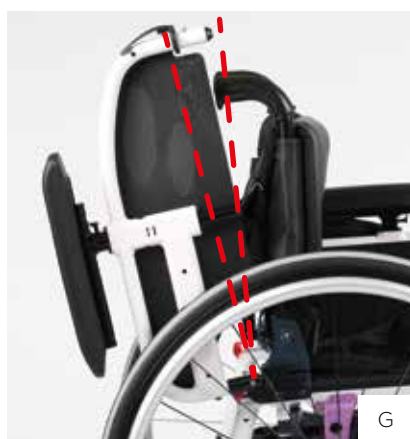

G



H

### ARMRESTS

#### A) Mudguard

The stylish mudguards are lightweight keeping the overall weight of the wheelchair to a minimum.

#### B) Conventional height adjustable armrest

Depth adjustable by sliding on the seat rail.

#### C) Hemi-armrest

This optional armrest provides additional comfort for the individual with its height, depth and angle adjustability.

### SWING AWAY & REMOVABLE ARMREST

#### D) Swing away and removable armrest

Easily accessible and can be flipped backwards to assist with transfers.

#### E) Height adjustment

Simple push button providing easy adjustment.

#### F) Release lever

Even with limited hand function, the armrest can easily be released and flipped backwards.






















#### G) Flip back position

A simple screw setting allows the position of the armrest to be fixed for easy access when flipped backwards.

#### H) Reversible armpad

The position of the armrest can be easily reversed providing a short and long armpad in one.

## TECHNICAL DATA

							
<b>Invacare Action<sup>5</sup></b>	305 / 330 / 355 / 380 / 405 / 430 / 455 / 480 / 505 / 530 mm	350 / 375 / 400 / 425 / 450 / 475 / 500 mm	0 - 9 °	Front 360-535 mm in steps of 25 mm Rear: 360-510 mm in steps of 25 mm	Fixed 90° Or stepless 90° -15°/+15°	Fixed 300 / 325 / 350 / 400 / 425 / 450 / 475 / 500 mm Adjustable 360-400 and 410-490 mm	Fixed front 80° Swing away legrest Elevating legrest 0-80°
							
<b>Invacare Action<sup>5</sup></b>	Swing away frame: 200-510 mm Fixed frame: 300-530 mm	186-227 mm 186-266 mm 205-280 mm	SW +160 mm To SW+ 250 mm	320 -380 mm	845 -1205 mm	670 -860 mm	701 -1108 mm
							
<b>Invacare Action<sup>5</sup></b>	7°	1050 -1380 mm	100 & 130 kg	From 12,3 kg	From 8,3 kg		

This product was tested and measured according to EN 12183:2014. Dimensions may vary by 15 mm / 2°. In addition, our wheelchairs comply with ISO 7176:2008 standards and with ISO 10542, and have additional features providing increased levels of occupant security and safety whilst travelling in a motor vehicle. Wheelchair users should however transfer to the vehicle seat and use the vehicle installed restraint system whenever it is feasible. For more pre-sales information about this wheelchair, please consult the [www.invacare.co.uk](http://www.invacare.co.uk) website. For information about typical user description, intended use, adjustments and transportation, please refer to the user guide.

## FRAME COLOURS

					
Dark Red	Pop Blue	Anthracite Black	Pearl White	Electric Green	Hypnotic Purple
					
Pearl Grey	Matt Black	Intense Yellow	Vibrant Orange		

## UPHOLSTERY



Black Nylon /  
Grey mesh

## COLOUR ACCENTS



Silver

Blue

Purple

Please note colours may vary slightly from those displayed above.

Invacare Limited  
Pencoed Technology Park -  
Pencoed - Bridgend CF35 5AQ -  
United Kingdom  
Tel: +44 1 656 776 222  
Fax: +44 1656 776 220  
E-mail: [uk@invacare.com](mailto:uk@invacare.com)  
Sales Order E-mail:  
[ordersuk@invacare.com](mailto:ordersuk@invacare.com)  
[www.invacare.co.uk](http://www.invacare.co.uk)

Invacare Ireland Ltd  
Unit 5 - Seatown Business Campus  
Seatown Road - SWORDS - County  
Dublin - Ireland  
Tel. +353 1 8107084  
Fax +353 1 8107085  
Email: [Ireland@invacare.com](mailto:Ireland@invacare.com)  
[www.invacare.ie](http://www.invacare.ie)

© 2015 Invacare International Sàrl  
All rights reserved.

All information quoted is believed to be correct at time of print. Invacare reserves the right to alter product specifications without prior consultation.

Action5 - UK - 09/2015 - 1588788



Follow us on Twitter:  
[@Think\\_Mobility](https://twitter.com/Think_Mobility)



Follow us on YouTube:  
Invacare Ltd - Think Mobility



**Yes, you can.®**