

Invacare® **Softform Premier Visco**

Pressure reducing Visco elastic mattress



The *Invacare **Softform Premier Visco*** is an anti-decubitus mattress featuring a deep visco insert pad, which is high-density and temperature-sensitive to the body, allowing it to conform to the user's shape. This affords exceptional comfort and facilitates pressure reduction by minimising peak pressures.



Contact: Ciarán McPhillips

Product Specialist Mobile:

087 943 8811

**HSE no.1 Mattress
Contract Award 7433-16
Lot3**

Trusted Solutions,
Passionate People



Pressure ulcer prevention

The temperature-sensitive visco elastic insert pad within the *Invacare Softform Premier Visco* warms to conform to the user's shape, maximising the body contact area and thereby facilitating good pressure reduction. The mattress cover is two-way stretch and waterproof, and positively helps to reduce shear and friction forces.

Affords exceptional comfort

The high quality visco elastic foam featured in this mattress delivers an exceptionally comfortable support surface.

Features and Options




Cover is resistant to fluid ingress

Features include two-way stretch waterproof polyurethane cover with high frequency welded seams and a full length flap concealing the zip. Also has a toughened PU coated base to prolong the longevity of the mattress.



Modular Construction

Enables easy parts replacement, increasing product lifetime. Insert pad can be interchanged with regular castellated foam, if required.

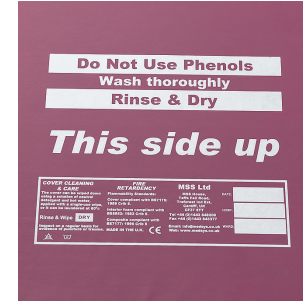

Does not require turning

Mattress does not require turning, significantly reducing the manual handling issues.

Technical Data

	Width	Length	Height	Total product weight	Max. user weight	Max. cleaning temperature	
							Manufactured to European Standards: BS EN 597-1 Cig, BS EN 597-2 Match, BS7177: 1996 Crib Source 5
Softform Premier Visco	Single size: 660 - 900 mm Double size: 910 - 1370 mm King size: 1380 - 1525 mm	1410 - 2080 mm	152 mm	13 kg	247.6 kg	80°	Available in Crib 5 & Crib 7.

Cover Options



*available at an additional charge

Warranty

Foam 4 years
Cover 4 years



'Strikethrough Resistant Technology'

The cover is made from a highly chemically resistant polyurethane polymer making it much more resistant to abrasion and damage and swells less with water contact. When cleaning a cover with 'Strikethrough Resistant Technology', the material changes less and thus reverts to its natural state very quickly. The time in which the cover is more prone to physical damage is significantly reduced.

Please contact your regional Product Specialist Ciarán McPhillips for further information Mobile: 087 943 8811

Phone: 01 810 7084
Fax: 01 810 7085
www.invacare.ie
Email: ordersireland@invacare.com

Invacare Ireland Ltd,
Unit 5, Seatown Business campus,
Seatown Road, Swords, Co. Dublin,
Ireland

