



Invacare®

# Matrix® Flo-tech™ Contour

Contoured moulded foam for improved posture

**matrix**  
SEATING SERIES

Trusted Solutions,  
Passionate People



The surface of the *Invacare **Matrx Flo-tech*** Contour cushion is contoured into sections. Its design supports the client and increases pressure redistribution. The combination of improved posture and increased support surface area will contribute significantly to pressure ulcer prevention. This deep cushion is suitable for heavier people.

### Moulded foam cushion

Moulding allows *Matrx Flo-tech* cushions to provide a three-dimensional form, which cannot be emulated in cut foam. Moulded foam is highly durable, increasing the product's longevity.

### Contoured seating surface

Offers improved stability and postural support; castellations reduce shear and friction forces on tissue.

## Features and Options



### Ingress Resistant Cover

Cover is two-way stretch, water-resistant and vapour-permeable. Cover features an anti-slip base with carry handle for easy transportation.

### Available in 83 sizes

The *Invacare **Matrx Flo-tech Contour*** is available in a range of 83 sizes, from 330 x 360 to 600 x 550 mm. (See size chart for details).

### Integral Sag Rigidizer option

Option of sag compensator to accommodate the sag in wheelchair seats, and ensure the cushion is correctly presented.

### Ultra-fresh™ added to all cushions for antimicrobial and odour protection

Antimicrobial treatment that has been added to all *Flo-tech* foam to help make it resistant to the growth of odour causing bacteria and helps provide longer lasting freshness.

## Technical data

For more comprehensive pre-sales information about this product, including the product's user manual, please see the [Invacare website](#).

	Various sizes	Height	Total product weight	Max. user weight <sup>1</sup>	Max. cleaning temperature	Warranty	
<i>Matrx Flo-tech Contour</i>							Manufactured to European Standards: BS EN 1021-1 Cig, BS EN 1021-1 Match, and BS 7176 (Crib 5) <b>Disclaimer</b> 1. Guide maximum user weight limits are set assuming that the cushion selected is of appropriate width and depth for the intended user/s.
	A range of 83 sizes	103 mm	1.38 kg	152 kg (24 stone) (191 kg on sizes 560 x 560 and 560 x 610 mm) (30 stone)	80°	3 years	



## Colours

Black Smooth

Black wipedown

Black spacer fabric

Baltic Blue

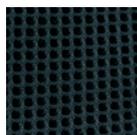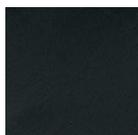

Tel: +44 1 656 776 222  
 Fax: +44 1656 776 220  
 E-mail: [uk@invacare.com](mailto:uk@invacare.com)  
 Sales Order E-mail:  
[ordersuk@invacare.com](mailto:ordersuk@invacare.com)  
[www.invacare.co.uk](http://www.invacare.co.uk)

Invacare Limited  
 Pencoed Technology Park  
 Pencoed - Bridgend CF35  
 5AQ - United Kingdom

Tel. +353 1 8107084  
 Fax +353 1 8107085  
 Email: [ireland@invacare.com](mailto:ireland@invacare.com)  
[www.invacare.ie](http://www.invacare.ie)

Invacare Ireland Ltd  
 Unit 5 - Seatown  
 Business Campus  
 Seatown Road  
 SWORDS - County  
 Dublin - Ireland



Yes, you can.®

