



Invacare[®] | **Matrix[®] Flo-tech[™] Lite**

Slimline contoured cushion affording exceptional comfort and improved posture

matrix[®]
SEATING SERIES

Trusted Solutions,
Passionate People



The *Invacare Matrix Flo-tech Lite* cushion is a slimline contoured cushion ideally suited for people at low risk of developing pressure ulcers and for comfort purposes only. The contoured seating surface offers improved stability and postural support. With sizes starting from 25 x 25 cm (250 x 250 mm) the *Matrix Flo-tech Lite* cushion is particularly suited to paediatric clients.

Features and Options



Ingress Resistant Cover

Cover is two-way stretch, water-resistant and vapour-permeable. Cover features an anti-slip base with carry handle for easy transportation.

Available in 61 sizes

The *Invacare Matrix Flo-tech Lite* is available in a range of 61 sizes, from 250 x 250 to 480 x 530 mm. (See prescription chart for details)

Integral sag Rigidizer option

Option to compensate the sag in wheelchair seats, and ensure the cushion is correctly presented.

Ultra-fresh™ added to all cushions for antimicrobial and odour protection

Antimicrobial treatment that has been added to all *Flo-tech* foam to help make it resistant to the growth of odour causing bacteria helps provide longer lasting freshness.

Slimline foam base made from high quality resilient, moulded foam

The high quality resilient foam is manufactured using a unique formula which offers high levels of pressure redistribution and improves product longevity. Moulding allows the *Matrix Flo-tech Lite* cushion to provide a 3D form, which cannot be emulated in cut foam. The pommel at the front of the cushion ensures femurs are in a midline position for good posture and allows air to circulate.

Contoured seating surface


Offers improved stability and postural support. Also delivers a greater degree of weight distribution compared to that of a flat surface. The castellations of the cushion allow slight movement which in turn reduces shear and friction forces.

Technical data

For more comprehensive pre-sales information about this product, including the product's user manual, please see the [Invacare website](#).

Various sizes	Height	Total product weight	Max. user weight ¹	Max. cleaning temperature	Warranty	
						Manufactured to European Standards: BS EN 1021-1 Cig, BS EN 1021-2 Match and BS 7176 (Crib 5)
A range of 61 sizes	78 mm	1.5 kg	127 kg (20 stone)	80°	3 years	Disclaimer 1. Guide maximum user weight limits are set assuming that the cushion selected is of appropriate width and depth for the intended user/s.

Colours

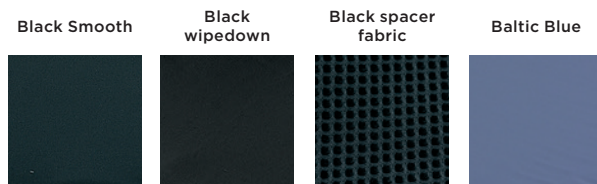

Tel: +44 1 656 776 222
Fax: +44 1656 776 220
E-mail: uk@invacare.com
Sales Order E-mail: ordersuk@invacare.com
www.invacare.co.uk

Invacare Limited
Pencoed Technology Park
Pencoed - Bridgend CF35
5AQ - United Kingdom

Tel. +353 1 8107084
Fax +353 1 8107085
Email: ireland@invacare.com
www.invacare.ie

Invacare Ireland Ltd
Unit 5 - Seatown
Business Campus
Seatown Road
SWORDS - County
Dublin - Ireland



Yes, you can.®

