

Invacare® | **Esprit Action⁴_{NG}**

NEW



Easy to transport, lightweight powerchair

Invacare **Esprit Action⁴_{NG}** is a compact and lightweight power wheelchair for those that like to get about easily. It's tool free removable wheels and compact battery pack make it ideal for transporting (heaviest item weighing only 15 Kg). It has lightweight wheels with an option for self propel and a reduced overall width (seat plus 170 mm) making it narrow enough to manoeuvre around the tightest of spaces.

The **Esprit Action⁴_{NG}** provides all the advantages of driving a power wheelchair indoor and outdoor, without the size and bulk of an average power wheelchair.



Yes, you can.®



Compact and lightweight

Thanks to the new chassis structure (available in 5 colours) and overall width, everything is smaller and narrower for getting around indoors. And with a specially designed frame the total width of the chair does not add more than 170 mm to the width of the seating. The **Esprit Action^{®4} NG** is easy to dismantle making it highly transportable.

Highly Configurable

The adult version has eight seat widths options from 380 to 605 mm and three seat depths from 400 to 500 mm.

The junior version offers four different seat widths of 305 to 380 mm and seat depths from 325 - 375 mm.

More options available

The armrest width allows for configuration adjustments of up to 50 mm (maximum) allowing for more customised seating and frame positions.

It can be equipped with a variety of configured options including backrest for Matrix positioning. **Esprit Action^{®4} NG** also provides the option of an attendant

control should this be required. There is a range of colours to choose from including Ocean Blue, Slate Grey, Pearl White, and Anthracite Black for the adult version.

For the junior version there is Ocean Blue, Pearl White, Bubble Gum Pink and Anthracite Black.

Features and Options



Dual intuitive pressure sensors

This allows an attendant to drive the wheelchair from behind in an easy yet intuitive manner.



Alber anti-tippers as standard

With Jack up function to help the quick release of wheels.



Welded castor mounts provide greater rigidity on the front of the chair.



Wheel protection bumper

Features and Options



Adaptable arm widths gives up to 50mm adjustment.



Simply push and turn to remove wheel.



Disassembled and folded chair

The heaviest part weighs only 15 kg (chassis).



Overall width reduced

The frames allow for the motors to sit inside and reduce the width.



Self Propel Wheels (Option)

Convert the chair to a manual chair.



Sleek yet Powerful Design

The **Esprit Action^{®4} NG** has been designed in such a way to allow all the configurable options you would expect in a manual chair, but incorporating the powerful drive wheels discreetly into the hubs.



Junior Version

Esprit Action^{®4} NG is available in two different sizes : adult or junior version.

Technical data

Esprit Action⁴ NG 380 - 605 mm	400 - 500 mm	460 or 485 mm	430 - 510 mm	350 - 480 mm	3°	0° - 30°
Esprit Action⁴ NG Junior 305-380 mm	325-375 mm	460 mm	430 - 510 mm	220 - 415 mm	3°	0° - 30°
Esprit Action⁴ NG Seat width + 170 mm	897 - 1002 mm	900 - 1100 mm	640 - 840 mm	45 kg*	125 kg	12 Ah (17 Ah)
Esprit Action⁴ NG Junior + 170 mm	897 - 977 mm	900 - 1100 mm	640 - 840 mm	45 kg*	80kg	12 Ah (17 Ah)
Esprit Action⁴ NG 1860 mm	30 mm	5° / 10%	150 W	14 km (12 Ah) 30 km (17 Ah)	6 km/h	Adjustable armrests up to 50 mm
Esprit Action⁴ NG Junior 1860 mm	30 mm	5° / 10%	150 W	14 km (12 Ah) 30 km (17 Ah)	6 km/h	up to 50 mm

* Weight of heaviest part is the frame at 15kg

Colours



Ocean Blue
(Black legrests on Junior Version)



Slate Grey
(Adult Version Only)



Pearl White
(Black legrests on Junior Version)



Anthracite Black



Bubble Gum Pink
(Junior Version Only)

Invacare Limited
Pencoed Technology Park - Pencoed - Bridgend CF35 5AQ - United Kingdom
Tel: +44 1 656 776 222
Fax: +44 1656 776 220
E-mail: uk@invacare.com
Sales Order E-mail: ordersuk@invacare.com
www.invacare.co.uk


Invacare Ireland Ltd
Unit 5 - Seatown Business Campus
Seatown Road - SWORDS - County Dublin - Ireland
Tel. +353 1 8107084
Fax +353 1 8107085
Email: Ireland@invacare.com
www.invacare.ie

© 2013 Invacare International Sàrl
All rights reserved.

All information quoted is believed to be correct at time of print.
Invacare reserves the right to alter product specifications without prior consultation.

Esprit Action4NG - UK - 04/2013 - 1561754

Upholstery

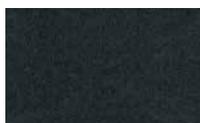

Black Nylon

Manufacture special options : various specials on the chair are available (size, accessories, colors with RAL). Please contact Customer Service. For more information, technical details, references and other options, please refer to the user manual wheelchair.



Yes, you can.®