

# Action<sup>3</sup><sub>NG</sub> Mid Wheel Propulsion



Optimum stability - Clever 6 wheel design

LIGHTWEIGHT | MANOEUVRABLE | FOLDABLE



The **Action<sup>3</sup><sub>NG</sub> Mid Wheel Propulsion (MWP)** is a 6 wheeled manual wheelchair, specially designed for those individuals who need increased stability, without compromising on manoeuvrability and comfort. A specific addition to the **Action family** with an extensive range of cross-compatible options and accessories, this highly configurable and lightweight foldable chair can be set-up to meet many individual needs. The forward position of the rear wheels provides easier access to handrims and enhances manoeuvrability. This outstanding design offers a unique resting position, with two additional rear castors and can easily be tilted back to achieve a 25 degree seat tilt.

- ▶ Enhanced stability and manoeuvrability compared to a standard 4 wheeled manual chair
- ▶ A smooth motion, mechanical tilt for positioning, that can be achieved individually or by a caregiver
- ▶ The two extra rear wheels ensure a very stable tilt position for safety
- ▶ Easy to store and transport with a lightweight foldable frame. Total wheelchair weight of 17.6kg. Transport weight of 11.2kg.

## Enhanced manoeuvrability

With its clever design and central wheel position the turning radius is shorter compared to a standard manual wheelchair. Therefore creating a responsive and easy to handle driving performance for effortless use indoors.

### Ergonomic efficiency:

- ▶ Easier access to the handrims
- ▶ Position of upper limbs are in a more natural and ergonomic position
- ▶ Less energy required for propelling
- ▶ Less stress on the shoulders with the inline position of the persons arms



### Integrated resting position

The **Action<sup>3</sup><sub>NG</sub> MWP** benefits from two 5" rear castors, allowing this foldable manual wheelchair to tilt to 25 degrees. This smart innovation enables a caregiver to tilt the lightweight chair effortlessly and also encourages independence by allowing the individual to 'self tilt' using their wheelchair skills if required.

Featuring this tilt function means any areas within the region of the seating position that require pressure relief,\* can be addressed quickly and easily, helping to prevent any potential pressure points arising.



\*It is widely accepted that there is a correlation between the use of tilt for pressure relief which is supported by a body of peer reviewed clinical evidence. Details of which can be found on [www.ncbi.nlm.nih.gov](http://www.ncbi.nlm.nih.gov).

# 6

## EXCLUSIVE FEATURES

### 1 Mechanical tilt position

- ▶ 25° tilt allows a comfortable resting position

### 3 Swing away detachable armrest

- ▶ Ergonomic foam handle armrests, facilitate lateral transfers.

### 2 Tilt security lever

- ▶ When activated, enables individuals to overcome small obstacles up to 65mm in height safely.
- ▶ When deactivated it allows the chair to be put into a tilting position.



### 4 Enhanced manoeuvrability

- ▶ Compared to a standard 4 wheel manual chair

### 5 Cable actuated brakes

- ▶ Follows the central wheel position
- ▶ **Action<sup>3</sup>NG MWP** can be securely parked even in a tilted position



### A compact design and a controlled weight

## 6

- ▶ The driving unit width and the total width folded are the same as the **Action<sup>3</sup>NG**, when folded:
  - Driving unit width: seat width + 190 mm
  - Total width folded: 295 mm

Total weight is 17.6 kg, therefore only + 3.4 kg heavier than **Action<sup>3</sup>NG**. The total length has been optimised to equal the **Action<sup>3</sup>NG** with anti-tippers.




## Technical data ▼

For more comprehensive pre-sales information about this product, including the product's user manual, please visit [www.invacare.co.uk](http://www.invacare.co.uk)

Seat width ▼	Seat depth ▼	Seat angle ▼	Seat height ▼	Backrest angle ▼	Backrest height ▼	Front frame angle ▼	
<b>Action<sup>3</sup>NG MWP</b>	380 / 405 / 430 / 455 / 480 / 505 mm	400 / 450 mm	3°	460 - 485 mm	7° / 0° - 30° / - 15° - +15°	350 - 510 mm	80° / 0 - 75°
Length of footrests ▼	Armrest height ▼	Driving unit width ▼	Total width, folded ▼	Total length inc. footrests ▼	Total length without footrests ▼	Total height ▼	
<b>Action<sup>3</sup>NG MWP</b>	200 - 400 / 330 - 460 / 440 - 510 mm	230 / 170 - 295 mm	Seat width + 190 mm	295 mm	1036 - 1124 mm	816 - 890 mm	810 - 1040 mm
Max. safe slope ▼	Radius ▼	Max. user weight ▼	Total weight ▼	Mass of the heaviest part ▼			
<b>Action<sup>3</sup>NG MWP</b>	4°	800 mm	125 kg	17,6 kg	11,2 kg		

## Frame colours ▼

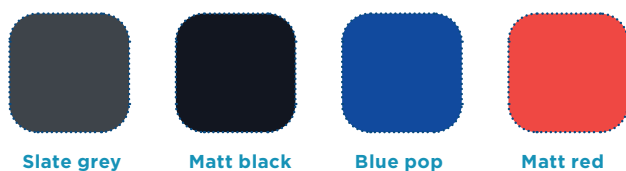

## Upholstery ▼



Please note: colours may vary slightly from those displayed here.

Invacare Limited  
Pencoed Technology Park - Pencoed  
Bridgend CF35 5AQ - United Kingdom  
Tel: +44 1 656 776 222  
Fax: +44 1656 776 220  
E-mail: [uk@invacare.com](mailto:uk@invacare.com)  
Sales Order E-mail: [ordersuk@invacare.com](mailto:ordersuk@invacare.com)  
[www.invacare.co.uk](http://www.invacare.co.uk)

Invacare Ireland Ltd  
Unit 5 - Seatown Business Campus  
Seatown Road - SWORDS  
County Dublin - Ireland  
Tel. +353 1 8107084  
Fax +353 1 8107085  
Email: [Ireland@invacare.com](mailto:Ireland@invacare.com)  
[www.invacare.ie](http://www.invacare.ie)



**Yes, you can.®**



© 2019 Invacare International GmbH. All rights reserved. All information quoted is believed to be correct at time of print. Invacare reserves the right to alter product specifications without prior consultation.  
Action 3NG MWP - 1631793 - UK - 10/2019