

AQUATEC®

Ocean Ergo



Modular personal care solutions

SAFE | INDEPENDENT | ERGONOMIC



Design-focused on user experience

to enhance independence, safety and flexibility.

Personal care is highly private activity and is an important part of a person's life and wellbeing. The **Aquatec Ocean Ergo** shower chair commodes have been designed with user comfort, independence and dignity in mind. Based on the well-established Aquatec Ocean range, the Ocean Ergo comes packed with a range of new features, designed to make personal care activities easier for both the user and carer.

Safety, independence and dignity

A newly contoured seat surface has been developed to encourage a more upright and stable seating position, promoting independence and overall

comfort. Enhanced stability at the front of the base allows for a smaller footprint and offers users added safety and reassurance when seated in the chair. Furthermore, the optimised wheelbase makes it even easier for a carer to manoeuvre the chair in small spaces as well as allowing for a better fit over standard toilets for independent toileting.

Modular and flexible design

The modular design of the Ocean Ergo means that the chairs can be adapted to suit a wide range of user needs by making some simple adjustments, or by adding some of our easy to fit, off-the-shelf accessories.

NEW Contoured, ergonomic seat plate

- ▶ Based on the physical imprint left behind when seated
- ▶ A 5° seat tilt lifts the knees into a natural squat position for toileting
- ▶ Transfer supports at the front to allow the user to position themselves into a comfortable position and aid side transfers for active users



OTHER features include:

- ▶ Stable, stainless steel frame
- ▶ Machine washable, tension adjustable back with quick-release buckles
- ▶ Easy disassembly for decontamination, transport and storage

NEW FEATURES

1 Innovative ergonomic seat

- ▶ Insight-led design encourages a more upright seating position for enhanced independence and comfort

2 Greater adaptability



- ▶ With an increased height adjustment range on the footrests and armrests

4 Compact and easy to handle

- ▶ Optimised wheel base makes it easy to manoeuvre in small spaces

3 Extremely stable

- ▶ Enhanced stability at the front for additional safety and reassurance



Highly Configurable

Available with a extensive range of positioning, safety, toileting and seating accessories that can be used to support a wide range of conditions and postural needs.

See full range of accessories online at www.invacare.co.uk



Positioning



Safety



Toileting



Seating and Backrests

Technical data ▼

For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.



Seat width ▼



Seat depth ▼



Total width ▼



Total depth ▼



Width between armrests ▼

Ocean Ergo	480 mm	480 mm	565 mm	905 mm	455 mm
Ocean Ergo XL	480 mm	480 mm	645 mm	905 mm	535 mm
Ocean 24" Ergo	480 mm	480 mm	670 mm	1050 mm	455 mm
Ocean 24" Ergo XL	480 mm	480 mm	670 mm	1050 mm	535 mm




Weight ▼



Load capacity ▼

Ocean Ergo	17 kg	150 kg
Ocean Ergo XL	18 kg	180 kg
Ocean 24" Ergo	23 kg	150 kg
Ocean 24" Ergo XL	24 kg	180 kg

Frame colour ▼

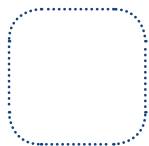

White

Regulatory data ▼



The Aquatec Ocean Ergo is part of the Safe Patient Handling range.



Safe Patient Handling

Invacare Limited
Pencoed Technology Park - Pencoed
Bridgend CF35 5AQ - United Kingdom
Tel: +44 1 656 776 222
Fax: +44 1656 776 220
E-mail: uk@invacare.com
Sales Order E-mail: ordersuk@invacare.com
www.invacare.co.uk

Invacare Ireland Ltd
Unit 5 - Seatown Business Campus
Seatown Road - SWORDS
County Dublin - Ireland
Tel. +353 1 8107084
Fax +353 1 8107085
Email: Ireland@invacare.com
www.invacare.ie



YouTube

