


A dual surface specialist mattress designed for spinal injuries

The *Invacare Softform Premier Spinal* is a high specification replacement mattress designed for those with spinal conditions or injury. It is suitable for those considered at **'Very High Risk'**¹ of developing pressure ulcers. The **Spinal** mattress provides optional support to cover a wide range of comfort and spinal alignment requirements. Following a spinal injury, a firm support surface may be necessary to assist in immobilising the spine. The unique dual layer design provides the option of a softer, **comfort foam** offering greater pressure redistribution, which may be more appropriate at a later stage in the rehabilitation process.

Trusted Solutions, Passionate People



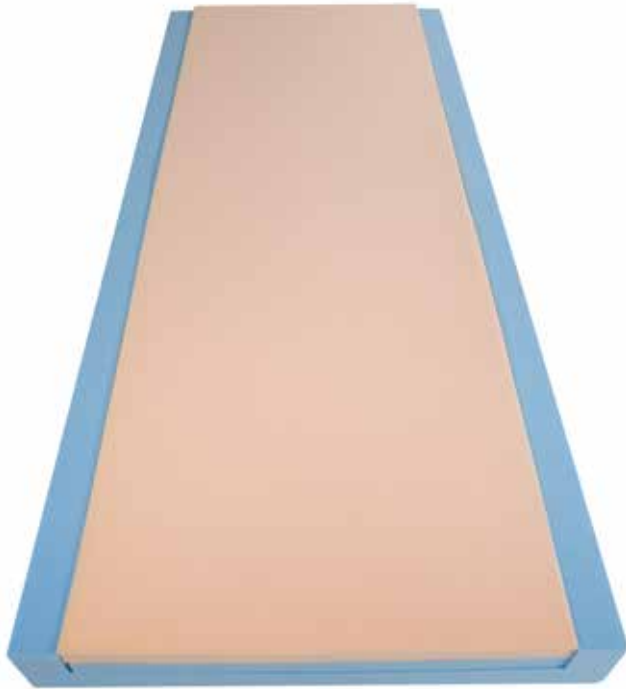

Comfort surface



Firm surface

Specialist stable support surface

Back injury and pain can vary greatly and is specific to each individual. The **Softform Premier Spinal** mattress offers a dual interface option which provides a bespoke solution without the need to change to an alternative mattress.

The interface foam can be quickly and easily changed at the bedside to reduce the risk of waiting time for the patient and limit the need for lengthy transfer.


The specifically manufactured foam sits level with firm side supports and promotes effective patient care, handling and transfer, without compromise to spinal alignment or pressure redistribution.

The *Invacare Softform Premier Spinal*'s unique dual foam insert provides a choice of two stable support surfaces with excellent comfort capabilities. The choice of which is individual to each user, shear reduction and pressure redistribution is addressed with either surface.

The design of the mattress allows for excellent conformity if profiled on electric bed frames. Practical and durable, the **Softform Premier Spinal** mattress delivers an exceptional level of comfort and pressure redistribution, and can be a valuable contributor to the rehabilitation process.

Mattress Support


surface too soft

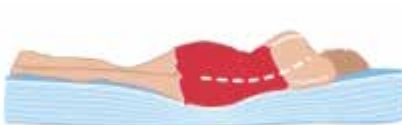

surface too hard



Correct support & spinal alignment

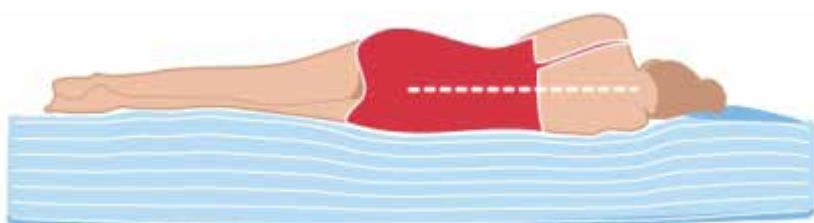

Features



Infection Control

An integral plaiton liner forms part of the interior of the **Softform Premier Spinal** mattress. As the risks of contact with bodily fluids may be increased, the liner which encases the internal foam attributes to infection control by reducing the risk of fluid ingress. In combination with the **SRT** (Strike through Resistant Technology)² cover, this secondary protection offers a double barrier to mattress contamination.

SRT covers are unique to Invacare, and have been designed to be “strike through” resistant. Strike through occurs when the mattress cover fibres break down or separate, and small cracks develop resulting in contaminants seeping through to the inner foam usually as a result of damage by rigorous cleaning. These specially coated **SRT** covers are designed to withstand the stringent but essential cleaning regimes in place in healthcare settings and as a result, reduce mattress failure rates.




Unique dual foam insert provides two surface options

The blue insert provides a high density surface for additional stability, when the orange insert is placed in the upper position, the patient immerses into the mattress providing comfort and pressure reduction.




Does not require turning or rotating

This mattress does not require turning or rotating, significantly minimising manual handling issues.

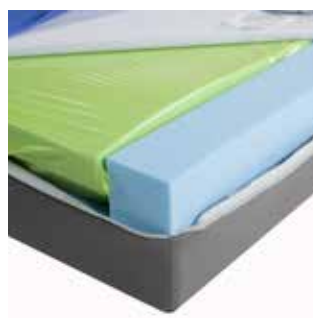

Single piece foam “U” core

Facilitates patient transfer by providing effective side-wall support. A robust foam base prevents the mattress from bottoming out.

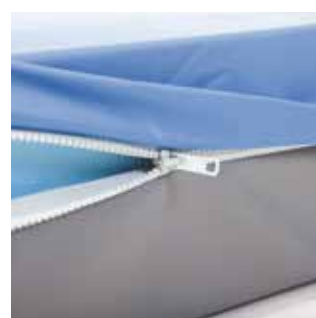

Full length flap concealing zip

A full length flap protects the zip. This reduces the possible risk of fluid ingress, preventing contamination of the foam core.

Multi stretch, breathable, waterproof cover with corner welds, reducing the opportunity for accumulation of dust and liquid residues.

Technical data

						
Single	660-900mm Double and king size sizes available on request	1250-2100mm	152 mm	14 kg	Max. User Weight: 247.6 kg (39 Stone)	Cleaning: 80°

Manufactured to European Standards: BS EN 597-1 Cig, BS EN 597-2 Match.
Foam Interior: BS 5852: 1990 Crib Source 5; Cover BS7175: 1989 Crib Source 5.
Composite tested to BS 7177 Crib 5.

Warranty

Foam 4 years
Cover 2 years

Colour



Violet



Invacare Limited
Pencoed Technology Park -
Pencoed - Bridgend CF35 5AQ
- United Kingdom
Tel: +44 1 656 776 222
Fax: +44 1656 776 220
E-mail: uk@invacare.com
Sales Order E-mail:
ordersuk@invacare.com
www.invacare.co.uk

Invacare Ireland Ltd
Unit 5 - Seatown Business
Campus
Seatown Road - SWORDS -
County Dublin - Ireland
Tel. +353 1 8107084
Fax +353 1 8107085
Email: Ireland@invacare.com
www.invacare.ie

© 2016 Invacare International Sàrl
All rights reserved.

All information quoted is believed to be correct at time of print.
Invacare reserves the right to alter product specifications without prior consultation.

Softform Premier Spinal - UK -
08/2016 - 1576072



Yes, you can.®

- 1/ Essential nursing care is pivotal in pressure ulcer prevention. This mattress will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action based planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. A range of assessment scales can be used as a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value. The risk category in this instance is based on the Waterlow Scale.
- 2/ Please note: The cover is not Strikethrough proof. Misuse or incorrect care and handling can cause damage to the mattress cover.