



Invacare[®]

Aero-tech Uno

Free movement of air

Trusted Solutions,
Passionate People



Clinical Advantage

The *Invacare Aero-tech Uno* is a single valve air cushion which is specifically designed for those considered to be at 'Very High Risk' of developing pressure ulcers. Individual interconnected air cells allow free movement of air within the cushion, which allows the level of air to be individually adjusted to body shape thus maximising immersion. This air fluidised movement minimises changes in peak pressures, therefore providing excellent protection against pressure damage.

The natural rubber material and vulcanised bonded base, offers excellent product longevity and may be suitable for those who experience allergic reactions to latex¹.

The *Invacare Aero-tech Uno* addresses critical pressure points providing the client with a high level of skin protection. In addition the continual movement of air provides protection for the ischial tuberosities and coccyx area. This cushion provides an infinitely adjustable, low friction and low shear environment which can protect vulnerable clients against pressure damage.

Features and Options



Air Valve Adjustment

Incremental air valve adjustment allows the level of air to be adjusted to suit individuals' needs for pressure reduction, stability and comfort.



Custom Sizes Available

A range of cushion sizes are available, ensuring a perfect fit for all makes and models of wheelchair. (For size guide, please refer to order form).



Rounded Cell Shape

The rounded shape of the base of the cells ensures good immersion for individuals.



Bonded Base Layer

The vulcanized bonded base layer of the cushion ensures a robust seal, prolonging longevity.

Technical data


For more comprehensive pre-sales information about this product, including the product's user manual, please see the [Invacare website](#).

	Width	Length	Height	Min. / Max. user weight	Max. cleaning temperature	Warranty	
							Manufactured to European Standards: BS EN 1021-1 Cig, BS EN 1021-1 Match, and BS 7176 (Crib 5)
Aero-tech Uno	380 - 480 mm	340 - 530 mm	50 & 100 mm	20 kg / 150 kg (3-23 stone)	60°	2 years	

Colours

Green rubber

Black stretch fabric

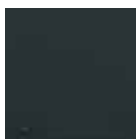

1. The **Aero-tech Uno** is free of latex allergen (ref: BGFA Bochum). Risk assessment must be carried out prior to use.
2. Risk Category - Moderate/High Risk

Tel: +44 1 656 776 222
Fax: +44 1656 776 220
E-mail: uk@invacare.com
Sales Order E-mail:
ordersuk@invacare.com
www.invacare.co.uk

Invacare Limited
Pencoed Technology Park
Pencoed - Bridgend CF35
5AQ - United Kingdom

Tel. +353 1 8107084
Fax +353 1 8107085
Email: ireland@invacare.com
www.invacare.ie

Invacare Ireland Ltd
Unit 5 - Seatown
Business Campus
Seatown Road
SWORDS - County
Dublin - Ireland



Yes, you can.®

