



*Invacare*<sup>®</sup> | **Aero-tech 4**

For customised pressure care and positioning

Trusted Solutions,  
Passionate People



The *Invacare Aero-tech 4* is a four valve air cushion which is specifically designed for those considered to be at 'Very High Risk' of developing pressure ulcers but require improved stability. Individual interconnected air cells allow free movement of air within the cushion, which allows the level of air to be individually adjusted to body shape/weight thus maximising immersion. This air fluidised movement minimises changes in peak pressures, therefore providing excellent protection against pressure damage.

The natural rubber material and vulcanised bonded base, offers excellent product longevity and may be suitable for those who experience allergic reactions to latex<sup>1</sup>.

## Clinical Advantage

The *Invacare Aero-tech 4* addresses critical pressure points providing the client with a high level of skin protection. In addition the continual movement of air provides protection for the ischial tuberosities and coccyx area. This cushion provides an infinitely adjustable, low friction and low shear environment which can protect vulnerable clients against pressure damage. The four separate air chambers improve stability, and can be individually adjusted to accommodate mild postural requirements.

## Features and Options



### Air Valve Adjustment

Incremental air valve adjustment allows the level of air to be adjusted to suit individuals' needs for pressure reduction, stability and comfort.



### Custom Sizes Available

A range of cushion sizes are available, ensuring a perfect fit for all makes and models of wheelchair. (For size guide, please refer to order form).



### Rounded Cell Shape

The rounded shape of the base of the cells ensures good immersion for individuals.



### Optional Cover Plates

Rigid plastic cover plates are available in a range of sizes and can be used to suppress cells.









### Bonded Base Layer

The vulcanized bonded base layer of the cushion ensures a robust seal, prolonging longevity.

## Technical data


For more comprehensive pre-sales information about this product, including the product's user manual, please see the [Invacare website](#).

|                    | Width   | Length  | Height  | Min. / Max. user weight   | Max. cleaning temperature   | Warranty   |   |
|--------------------|---|---|---|---|---|--|---|
| <b>Aero-tech 4</b> | <br>380 - 470 mm | <br>330 - 530 mm | <br>100 mm | <br>20 kg / 150 kg<br>(3-23 stone) | <br>60° | <br>2 years | Manufactured to European Standards:<br>BS EN 1021-1 Cig, BS EN 1021-1 Match, and BS 7176 (Crib 5) |

## Colours

Green rubber

Black stretch fabric

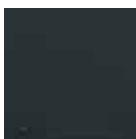

### Disclaimer

- The **Aero-tech 4** is free of latex allergen (ref: BGFA Bochum). Risk assessment must be carried out prior to use.

Tel: +44 1 656 776 222  
Fax: +44 1656 776 220  
E-mail: [uk@invacare.com](mailto:uk@invacare.com)  
Sales Order E-mail:  
[ordersuk@invacare.com](mailto:ordersuk@invacare.com)  
[www.invacare.co.uk](http://www.invacare.co.uk)

Invacare Limited  
Pencoed Technology Park  
Pencoed - Bridgend CF35  
5AQ - United Kingdom

Tel. +353 1 8107084  
Fax +353 1 8107085  
Email: [ireland@invacare.com](mailto:ireland@invacare.com)  
[www.invacare.ie](http://www.invacare.ie)

Invacare Ireland Ltd  
Unit 5 - Seatown  
Business Campus  
Seatown Road  
SWORDS - County  
Dublin - Ireland



**Yes, you can.®**

